

Acumed Implants in the MR Environment



Table of Contents

Introduction	4
Terminology	5
Warnings & Precautions	6
IFU and Surgical Technique	6
Anatomic Radial Head Portfolio	7
Anatomic Radial Head System	7
Anatomic Radial Head Solutions	8
Anatomic Radial Head Solutions 2 (stems and heads only)	9
Intramedullary Nails Portfolio	14
Polarus® Humeral Rod System	14
Polarus® Plus Humeral Rod System	15
Polarus® 3 Solution (Locking and Proximal Locking Nails)	16
Fibula Rod System	18
Forearm Fracture Solutions (Radius & Ulna Rods)	19
Fibula Nail 2 System	20
Ulna Nail 2 System	23
Hand and Wrist Portfolio	27
Acu-Loc® Wrist Plating System	27
Acu-Loc® 2 Wrist Plating System	29
Acu-Loc® Wrist Spanning Plate	35
Forearm Fracture Solutions (plates)	38
Hand Fracture System	40
Modular Hand System	43
Osteotomy System	44
Total Wrist Fusion Plating System	45
Screw and Pin Portfolio	48
Acutrak® Headless Compression Screw System	48
Acutrak 2® Headless Compression Screw System	50
Acutrak® 3 Headless Compression Screw System	52
AcuTwist® Acutrak Compression Screw System	56
Biotrak® Headless Resorbable Compression Screw System	57
Cannulated Screw System	57
Extremity Screw System	64
Tension Band Pin System / Tension Band Pin System 2	66

Shoulder and Elbow Portfolio	67
Elbow Plating System	67
Congruent Elbow Plating System	75
Radial Head Plating System	79
Clavicle Plating System	80
Clavicle Hook Plating System	87
Locking Clavicle Plating System	89
Acu-Sinch® Repair System	93
Dual-Trak Clavicle Screw System	93
Polarus® 3 Solution (plates)	94
Polarus® Proximal Humeral Plating System	96
Scapula Plating System	98
Foot and Ankle Portfolio	104
Ankle Plating System	104
Ankle Plating System 3	110
Calcaneal Plating System	114
Forefoot-Midfoot Plating System	121
Small Fragment Base Set	129
Hip and Pelvic Portfolio:	137
Pelvic Plating System	137
External Fixation Systems	141
Stableloc External Fixation System	141
Small Bone Distractor	141
Small Bone Fixator	141
Scapholunate Repair System	142
Scapholunate Repair System	142
Tendon and Ligament Portfolio	143
Acu-Sinch Knotless Ankle Syndesmosis Repair System	143
PEEK Suture Anchor System	144
Chest Wall Portfolio	145
RibLoc U Plus Chest Wall Plating System	145

Introduction

Magnetic Resonance Imaging (MRI) is one of many imaging modalities used by radiologists to diagnose and treat various conditions within the body. MRI is often favored for its detail and clarity when imaging soft tissues, and doing so without delivering the ionizing radiation inherent with X-ray and CT. Because of the strong, static magnetic field, and varying magnetic and electric fields during scanning, the safety of patients in the MR environment is of great importance to the healthcare professionals who perform these procedures.

Safety concerns arise when a patient has orthopaedic implants manufactured from metal or carbon fiber, or when external fixation devices are in place, because of the potential for movement or heating of the devices in the MR environment. In some cases, safety concerns may be greater due to the type, location, or configuration of the implant or the presence of multiple implants in close proximity to one another. Consequently, the use of MRI technology may not always be possible.

Purpose

This document will communicate MRI safety information regarding Acumed implants in the MR environment. While we cannot make specific MR claims for all of our implants at this time, Acumed will present the results of our testing in this document as additional information becomes available. When allowed by regulatory authorities, Acumed will use one of the internationally recognized ASTM symbols in our implant labeling.

Implant Materials

Acumed implants are commonly manufactured from commercial pure (CP) titanium (Ti) and titanium alloys, such as Ti-6Al-4V ELI. A recent FDA article¹ states:

“...[Test] results indicate that for alloys Ti-6Al-7Nb, CP titanium Grades 1 through 4, Ti-6Al-4V ELI, Ti-15Mo, Biodur 108, and Co-28Cr-6Mo ..., magnetically induced force is less than the device weight in any current clinical 1.5, 3.0, and 7.0-T MR system.”

Accordingly, many Acumed devices do not pose significant risk for magnetically induced force or torque. Refer to the device label and Instructions for Use to determine the implant material.

Terminology

ASTM International is an international standards organization that has published the following voluntary consensus standards for the testing and labeling of devices that may be used in MRI:

- ▶ F2052 Test Method for Measurement of Magnetically Induced Displacement Force on Medical Devices in the Magnetic Resonance Environment
- ▶ F2119 Test Method for Evaluation of MR Image Artifacts from Passive Implants
- ▶ F2182 Test Method for Measurement of Radio Frequency Induced Heating On or Near Passive Implants During Magnetic Resonance Imaging
- ▶ F2213 Test Method for Measurement of Magnetically Induced Torque on Medical Devices in the Magnetic Resonance Environment
- ▶ F2503 Practice for Marking Medical Devices and Other Items for Safety in the Magnetic Resonance Environment

ASTM F2503 defines the three classifications for products used in MRI as MR Safe, MR Conditional, and MR Unsafe. Products in these classifications may be labeled with the following symbols:



MR Safe

MR Safe is an item that poses no known hazards resulting from exposure to any MR environment. MR Safe items are composed of materials that are electrically nonconductive, nonmetallic, and nonmagnetic. The MR Safe icon consists of the letters “MR” surrounded by a green square.



MR Conditional

MR Conditional is an item with demonstrated safety in the MR environment within defined conditions. At minimum, the conditions of the static magnetic field, the switched gradient magnetic field, and the radiofrequency fields must be disclosed. Additional conditions may be required. The MR Conditional icon consists of the letters “MR” within a yellow equilateral triangle with a thick black band around the perimeter.



MR Unsafe

MR Unsafe is an item that poses unacceptable risks to the patient, medical staff, or other persons within the MR environment. The MR Unsafe marking consists of the letters “MR” surrounded by a red circle with a diagonal red bar across the letters extending from the upper left quadrant to the lower right quadrant of the circle and oriented at 45° from the horizontal.

The symbols may appear in black and white when not practical to print in the standard colors.

Warnings & Precautions

Unless explicitly stated in this document, Acumed devices have not been tested for heating, migration, or image artifact in the magnetic resonance (MR) environment and the safety of these devices in the MR environment is unknown. Scanning a patient with untested implants may result in patient injury.

Patients with implants should seek a medical opinion before entering the MR environment.

For Devices That Have Been Tested

Warning

- ▶ Failure to follow the specific MR conditions may result in injury to the patient.
- ▶ Unconscious patients or those insensitive to heat or pain may not inform the technologist that the implant is being affected by the MR environment.
- ▶ Long scanning time or higher specific absorption rate (SAR) increases the risk of heating. Pause the scan periodically and solicit feedback from the patient.
- ▶ Unknown and untested implant configurations or the addition of implants from other manufacturers may increase the risks of heating and image artifact.

Caution

- ▶ Local transmit coils are not recommended for use on the area with implants as their effects have not been tested.
- ▶ Maximize the distance from the implant to the image area to minimize any potential artifact, when possible.
- ▶ Maximize the distance between the implant and sides of the bore, coils, and conductors to minimize the risk of heating, when possible.
- ▶ The implants have not been assessed for compatibility with MRI systems above 3 Tesla.

Important

- ▶ An increase in ventilation through the bore may help to lessen the increase in body temperature during scanning and reduce the effects of implant heating.
- ▶ Nonclinical testing does not account for the variations in MRI systems, pulse sequences, multi-channel RF coils, and local transmit coils that might be used. Some MRI systems have image artifact reduction sequences if imaging near a device is necessary.
- ▶ The use of low field MRI systems (less than 1 T) can mitigate risks for patients with implants by lowering the frequency and SAR.

IFU and Surgical Technique

Acumed offers Instructions for Use (IFU) and Surgical Techniques to promote the safe and effective use of our systems. Visit www.acumed.net for the most current versions.

Important: IFUs and Surgical Techniques may contain important safety information. They are intended for the Operating Surgeon and supporting Healthcare Professionals.

Anatomic Radial Head Portfolio

- ▶ Anatomic Radial Head System
- ▶ Anatomic Radial Head Solutions
- ▶ Anatomic Radial Head Solutions 2

Anatomic Radial Head System

Sterile Part Number	Non-Sterile Part Number	Description
TR-H200L-S	N/A	20.0 mm Head, Left
TR-H200R-S	N/A	20.0 mm Head, Right
TR-H220L-S	N/A	22.0 mm Head, Left
TR-H220R-S	N/A	22.0 mm Head, Right
TR-H240L-S	N/A	24.0 mm Head, Left
TR-H240R-S	N/A	24.0 mm Head, Right
TR-H260L-S	N/A	26.0 mm Head, Left
TR-H260R-S	N/A	26.0 mm Head, Right
TR-H280L-S	N/A	28.0 mm Head, Left
TR-H280R-S	N/A	28.0 mm Head, Right
TR-S0600-S	N/A	6.0 mm x 0.0 mm Stem
TR-S0602-S	N/A	6.0 mm x 2.0 mm Stem
TR-S0604-S	N/A	6.0 mm x 4.0 mm Stem
TR-S0606-S	N/A	6.0 mm x 6.0 mm Stem
TR-S0608-S	N/A	6.0 mm x 8.0 mm Stem
TR-S0700-S	N/A	7.0 mm x 0.0 mm Stem
TR-S0702-S	N/A	7.0 mm x 2.0 mm Stem
TR-S0704-S	N/A	7.0 mm x 4.0 mm Stem
TR-S0706-S	N/A	7.0 mm x 6.0 mm Stem
TR-S0708-S	N/A	7.0 mm x 8.0 mm Stem
TR-S0800-S	N/A	8.0 mm x 0.0 mm Stem
TR-S0802-S	N/A	8.0 mm x 2.0 mm Stem
TR-S0804-S	N/A	8.0 mm x 4.0 mm Stem
TR-S0806-S	N/A	8.0 mm x 6.0 mm Stem
TR-S0808-S	N/A	8.0 mm x 8.0 mm Stem
TR-S0900-S	N/A	9.0 mm x 0.0 mm Stem
TR-S0902-S	N/A	9.0 mm x 2.0 mm Stem

Anatomic Radial Head Portfolio [continued]

Anatomic Radial Head System

Sterile Part Number	Non-Sterile Part Number	Description
TR-S0904-S	N/A	9.0 mm x 4.0 mm Stem
TR-S0906-S	N/A	9.0 mm x 6.0 mm Stem
TR-S0908-S	N/A	9.0 mm x 8.0 mm Stem
TR-S1000-S	N/A	10.0 mm x 0.0 mm Stem
TR-S1002-S	N/A	10.0 mm x 2.0 mm Stem
TR-S1004-S	N/A	10.0 mm x 4.0 mm Stem
TR-S1006-S	N/A	10.0 mm x 6.0 mm Stem
TR-S1008-S	N/A	10.0 mm x 8.0 mm Stem

Anatomic Radial Head Solutions

Sterile Part Number	Non-Sterile Part Number	Description
TR-H200L-S	N/A	20.0 mm Head, Left
TR-H200R-S	N/A	20.0 mm Head, Right
TR-H220L-S	N/A	22.0 mm Head, Left
TR-H220R-S	N/A	22.0 mm Head, Right
TR-H240L-S	N/A	24.0 mm Head, Left
TR-H240R-S	N/A	24.0 mm Head, Right
TR-H260L-S	N/A	26.0 mm Head, Left
TR-H260R-S	N/A	26.0 mm Head, Right
TR-H280L-S	N/A	28.0 mm Head, Left
TR-H280R-S	N/A	28.0 mm Head, Right
TR-S0600-S	N/A	6.0 mm x 0.0 mm Stem
TR-S0602-S	N/A	6.0 mm x 2.0 mm Stem
TR-S0604-S	N/A	6.0 mm x 4.0 mm Stem
TR-S0606-S	N/A	6.0 mm x 6.0 mm Stem
TR-S0608-S	N/A	6.0 mm x 8.0 mm Stem
TR-S0700-S	N/A	7.0 mm x 0.0 mm Stem
TR-S0702-S	N/A	7.0 mm x 2.0 mm Stem
TR-S0704-S	N/A	7.0 mm x 4.0 mm Stem

Anatomic Radial Head Portfolio [continued]

Anatomic Radial Head Solutions

Sterile Part Number	Non-Sterile Part Number	Description
TR-S0706-S	N/A	7.0 mm x 6.0 mm Stem
TR-S0708-S	N/A	7.0 mm x 8.0 mm Stem
TR-S0800-S	N/A	8.0 mm x 0.0 mm Stem
TR-S0802-S	N/A	8.0 mm x 2.0 mm Stem
TR-S0804-S	N/A	8.0 mm x 4.0 mm Stem
TR-S0806-S	N/A	8.0 mm x 6.0 mm Stem
TR-S0808-S	N/A	8.0 mm x 8.0 mm Stem
TR-S0900-S	N/A	9.0 mm x 0.0 mm Stem
TR-S0902-S	N/A	9.0 mm x 2.0 mm Stem
TR-S0904-S	N/A	9.0 mm x 4.0 mm Stem
TR-S0906-S	N/A	9.0 mm x 6.0 mm Stem
TR-S0908-S	N/A	9.0 mm x 8.0 mm Stem
TR-S1000-S	N/A	10.0 mm x 0.0 mm Stem
TR-S1002-S	N/A	10.0 mm x 2.0 mm Stem
TR-S1004-S	N/A	10.0 mm x 4.0 mm Stem
TR-S1006-S	N/A	10.0 mm x 6.0 mm Stem
TR-S1008-S	N/A	10.0 mm x 8.0 mm Stem
TR-SL06-S	N/A	6 mm Morse Taper Long Stem
TR-SL08-S	N/A	8 mm Morse Taper Long Stem
TR-SL10-S	N/A	10 mm Morse Taper Long Stem
TR-SL12-S	N/A	12 mm Morse Taper Long Stem

Anatomic Radial Head Solutions 2 (stems and heads only)

Sterile Part Number	Non-Sterile Part Number	Description
50-0056-S	N/A	6.0 mm x 0.0 mm Anatomic Radial Head Partial Blast Stem
50-0057-S	N/A	6.0 mm x 2.0 mm Anatomic Radial Head Partial Blast Stem
50-0058-S	N/A	6.0 mm x 4.0 mm Anatomic Radial Head Partial Blast Stem
50-0059-S	N/A	6.0 mm x 6.0 mm Anatomic Radial Head Partial Blast Stem
50-0060-S	N/A	6.0 mm x 8.0 mm Anatomic Radial Head Partial Blast Stem

Anatomic Radial Head Portfolio [continued]

Anatomic Radial Head Solutions 2 (stems and heads only)

Sterile Part Number	Non-Sterile Part Number	Description
50-0061-S	N/A	7.0 mm x 0.0 mm Anatomic Radial Head Partial Blast Stem
50-0062-S	N/A	7.0 mm x 2.0 mm Anatomic Radial Head Partial Blast Stem
50-0063-S	N/A	7.0 mm x 4.0 mm Anatomic Radial Head Partial Blast Stem
50-0064-S	N/A	7.0 mm x 6.0 mm Anatomic Radial Head Partial Blast Stem
50-0065-S	N/A	7.0 mm x 8.0 mm Anatomic Radial Head Partial Blast Stem
50-0066-S	N/A	8.0 mm x 0.0 mm Anatomic Radial Head Partial Blast Stem
50-0067-S	N/A	8.0 mm x 2.0 mm Anatomic Radial Head Partial Blast Stem
50-0068-S	N/A	8.0 mm x 4.0 mm Anatomic Radial Head Partial Blast Stem
50-0069-S	N/A	8.0 mm x 6.0 mm Anatomic Radial Head Partial Blast Stem
50-0070-S	N/A	8.0 mm x 8.0 mm Anatomic Radial Head Partial Blast Stem
50-0071-S	N/A	9.0 mm x 0.0 mm Anatomic Radial Head Partial Blast Stem
50-0072-S	N/A	9.0 mm x 2.0 mm Anatomic Radial Head Partial Blast Stem
50-0073-S	N/A	9.0 mm x 4.0 mm Anatomic Radial Head Partial Blast Stem
50-0074-S	N/A	9.0 mm x 6.0 mm Anatomic Radial Head Partial Blast Stem
50-0075-S	N/A	9.0 mm x 8.0 mm Anatomic Radial Head Partial Blast Stem
50-0076-S	N/A	10.0 mm x 0.0 mm Anatomic Radial Head Partial Blast Stem
50-0077-S	N/A	10.0 mm x 2.0 mm Anatomic Radial Head Partial Blast Stem
50-0078-S	N/A	10.0 mm x 4.0 mm Anatomic Radial Head Partial Blast Stem
50-0079-S	N/A	10.0 mm x 6.0 mm Anatomic Radial Head Partial Blast Stem
50-0080-S	N/A	10.0 mm x 8.0 mm Anatomic Radial Head Partial Blast Stem
5001-0518L-S	N/A	Anatomic Radial Head Solutions 2 Head 18 mm, Left
5001-0518R-S	N/A	Anatomic Radial Head Solutions 2 Head 18 mm, Right
5001-0520L-S	N/A	Anatomic Radial Head Solutions 2 Head 20 mm, Left
5001-0520R-S	N/A	Anatomic Radial Head Solutions 2 Head 20 mm, Right
5001-0522L-S	N/A	Anatomic Radial Head Solutions 2 Head 22 mm, Left
5001-0522R-S	N/A	Anatomic Radial Head Solutions 2 Head 22 mm, Right
5001-0524L-S	N/A	Anatomic Radial Head Solutions 2 Head 24 mm, Left
5001-0524R-S	N/A	Anatomic Radial Head Solutions 2 Head 24 mm, Right
5001-0526L-S	N/A	Anatomic Radial Head Solutions 2 Head 26 mm, Left

Anatomic Radial Head Portfolio [continued]

Anatomic Radial Head Solutions 2 (stems and heads only)

Sterile Part Number	Non-Sterile Part Number	Description
5001-0526R-S	N/A	Anatomic Radial Head Solutions 2 Head 26 mm, Right
5001-0528L-S	N/A	Anatomic Radial Head Solutions 2 Head 28 mm, Left
5001-0528R-S	N/A	Anatomic Radial Head Solutions 2 Head 28 mm, Right
50-0109-S	N/A	11.0 mm x 0.0 mm Anatomic Radial Head Partial Blast Stem
50-0110-S	N/A	11.0 mm x 2.0 mm Anatomic Radial Head Partial Blast Stem
50-0111-S	N/A	11.0 mm x 4.0 mm Anatomic Radial Head Partial Blast Stem
50-0112-S	N/A	11.0 mm x 6.0 mm Anatomic Radial Head Partial Blast Stem
50-0113-S	N/A	11.0 mm x 8.0 mm Anatomic Radial Head Partial Blast Stem
50-0114-S	N/A	12.0 mm x 0.0 mm Anatomic Radial Head Partial Blast Stem
50-0115-S	N/A	12.0 mm x 2.0 mm Anatomic Radial Head Partial Blast Stem
50-0116-S	N/A	12.0 mm x 4.0 mm Anatomic Radial Head Partial Blast Stem
50-0117-S	N/A	12.0 mm x 6.0 mm Anatomic Radial Head Partial Blast Stem
50-0118-S	N/A	12.0 mm x 8.0 mm Anatomic Radial Head Partial Blast Stem
TR-S0600-S	N/A	6.0 mm x 0.0 mm Stem
TR-S0602-S	N/A	6.0 mm x 2.0 mm Stem
TR-S0604-S	N/A	6.0 mm x 4.0 mm Stem
TR-S0606-S	N/A	6.0 mm x 6.0 mm Stem
TR-S0608-S	N/A	6.0 mm x 8.0 mm Stem
TR-S0700-S	N/A	7.0 mm x 0.0 mm Stem
TR-S0702-S	N/A	7.0 mm x 2.0 mm Stem
TR-S0704-S	N/A	7.0 mm x 4.0 mm Stem
TR-S0706-S	N/A	7.0 mm x 6.0 mm Stem
TR-S0708-S	N/A	7.0 mm x 8.0 mm Stem
TR-S0800-S	N/A	8.0 mm x 0.0 mm Stem
TR-S0802-S	N/A	8.0 mm x 2.0 mm Stem
TR-S0804-S	N/A	8.0 mm x 4.0 mm Stem
TR-S0806-S	N/A	8.0 mm x 6.0 mm Stem
TR-S0808-S	N/A	8.0 mm x 8.0 mm Stem
TR-S0900-S	N/A	9.0 mm x 0.0 mm Stem

Anatomic Radial Head Portfolio [continued]

Anatomic Radial Head Solutions 2 (stems and heads only)

Sterile Part Number	Non-Sterile Part Number	Description
TR-S0902-S	N/A	9.0 mm x 2.0 mm Stem
TR-S0904-S	N/A	9.0 mm x 4.0 mm Stem
TR-S0906-S	N/A	9.0 mm x 6.0 mm Stem
TR-S0908-S	N/A	9.0 mm x 8.0 mm Stem
TR-S1000-S	N/A	10.0 mm x 0.0 mm Stem
TR-S1002-S	N/A	10.0 mm x 2.0 mm Stem
TR-S1004-S	N/A	10.0 mm x 4.0 mm Stem
TR-S1006-S	N/A	10.0 mm x 6.0 mm Stem
TR-S1008-S	N/A	10.0 mm x 8.0 mm Stem
TR-S1100-S	N/A	11.0 mm x 0.0 mm Stem
TR-S1102-S	N/A	11.0 mm x 2.0 mm Stem
TR-S1104-S	N/A	11.0 mm x 4.0 mm Stem
TR-S1106-S	N/A	11.0 mm x 6.0 mm Stem
TR-S1108-S	N/A	11.0 mm x 8.0 mm Stem
TR-S1200-S	N/A	12.0 mm x 0.0 mm Stem
TR-S1202-S	N/A	12.0 mm x 2.0 mm Stem
TR-S1204-S	N/A	12.0 mm x 4.0 mm Stem
TR-S1206-S	N/A	12.0 mm x 6.0 mm Stem
TR-S1208-S	N/A	12.0 mm x 8.0 mm Stem
TR-SL06-S	N/A	6 mm Morse Taper Long Stem
TR-SL07-S	N/A	7 mm Morse Taper Long Stem
TR-SL08-S	N/A	8 mm Morse Taper Long Stem
TR-SL09-S	N/A	9 mm Morse Taper Long Stem
TR-SL10-S	N/A	10 mm Morse Taper Long Stem
TR-SL11-S	N/A	11 mm Morse Taper Long Stem
TR-SL12-S	N/A	12 mm Morse Taper Long Stem

Anatomic Radial Head Portfolio [continued]

Conditions for safe scanning in the MRI

A patient with the head or stem implants from the Anatomic Radial Head, Anatomic Radial Head Solutions, and Anatomic Radial Head Solutions 2 systems may be safely scanned under the following conditions. Failure to follow these conditions may result in injury to the patient.

Nominal values(s) of Static Magnetic Field	1.5 T or 3 T
Maximum Spatial Field Gradient	19 T/m (1900 gauss/cm)
RF Excitation	Circularly Polarized (CP)
RF Transmit Coil Type	Whole body transmit coil, Head RF transmit-receive coil
Maximum Whole Body SAR at 1.5 T and 3.0 T	2.0 W/kg
Limits on Scan Duration 1.5 T	2.0 W/kg whole body average SAR for 15 minutes of continuous RF (a sequence or back to back series/scan without breaks)
Limits on Scan Duration 3 T	2.0 W/kg whole body average SAR for 12 minutes of continuous RF (a sequence or back to back series/scan without breaks)
MR Image Artifact	The presence of this implant may produce an image artifact of 71 mm.
If information about a specific parameter is not included, there are no conditions associated with that parameter.	

Intramedullary Nails Portfolio

- ▶ Polarus® Humeral Rod System
- ▶ Polarus Plus Humeral Rod
- ▶ Polarus 3 Solution (nails)
- ▶ Fibula Rod System
- ▶ Forearm Fracture Solutions (nails)
- ▶ Fibula Nail 2 System
- ▶ Ulna Nail 2 System

Polarus® Humeral Rod System

Sterile Part Number	Non-Sterile Part Number	Description
HCA5200-S	N/A	5.0 mm x 20.0 mm Cancellous Screw
HCA5250-S	N/A	5.0 mm x 25.0 mm Cancellous Screw
HCA5300-S	N/A	5.0 mm x 30.0 mm Cancellous Screw
HCA5325-S	N/A	5.0 mm x 32.5 mm Cancellous Screw
HCA5350-S	N/A	5.0 mm x 35.0 mm Cancellous Screw
HCA5375-S	N/A	5.0 mm x 37.5 mm Cancellous Screw
HCA5400-S	N/A	5.0 mm x 40.0 mm Cancellous Screw
HCA5450-S	N/A	5.0 mm x 45.0 mm Cancellous Screw
HCA5500-S	N/A	5.0 mm x 50.0 mm Cancellous Screw
HCA5550-S	N/A	5.0 mm x 55.0 mm Cancellous Screw
HCA-L5325-S	N/A	5.3 mm x 25 mm Cancellous Locking Screw
HCA-L5330-S	N/A	5.3 mm x 30 mm Cancellous Locking Screw
HCA-L5335-S	N/A	5.3 mm x 35 mm Cancellous Locking Screw
HCA-L5340-S	N/A	5.3 mm x 40 mm Cancellous Locking Screw
HCA-L5345-S	N/A	5.3 mm x 45 mm Cancellous Locking Screw
HCO3175-S	N/A	3.5 mm x 17.5 mm Cortical Screw
HCO3200-S	N/A	3.5 mm x 20.0 mm Cortical Screw
HCO3225-S	N/A	3.5 mm x 22.5 mm Cortical Screw
HCO3250-S	N/A	3.5 mm x 25.0 mm Cortical Screw
HCO3275-S	N/A	3.5 mm x 27.5 mm Cortical Screw
HCO3300-S	N/A	3.5 mm x 30.0 mm Cortical Screw
HCO3325-S	N/A	3.5 mm x 32.5 mm Cortical Screw
HCO3350-S	N/A	3.5 mm x 35.0 mm Cortical Screw
HCO3400-S	N/A	3.5 mm x 40.0 mm Cortical Screw
HCO3450-S	N/A	3.5 mm x 45.0 mm Cortical Screw
HCO3500-S	N/A	3.5 mm x 50.0 mm Cortical Screw
HCO3550-S	N/A	3.5 mm x 55.0 mm Cortical Screw

Intramedullary Nails Portfolio [continued]

Polarus® Humeral Rod System

Sterile Part Number	Non-Sterile Part Number	Description
HCO3600-S	N/A	3.5 mm x 60.0 mm Cortical Screw
HR-0001-S	N/A	Polarus Cap Screw Assembly
HR-0050-S	N/A	Polarus Cap Screw Assembly
HR-L1115-S	N/A	11 mm x 150 mm Polarus Locking Humeral Rod

Polarus® Plus Humeral Rod System

Sterile Part Number	Non-Sterile Part Number	Description
HCA5200-S	N/A	5.0 mm x 20.0 mm Cancellous Screw
HCA5250-S	N/A	5.0 mm x 25.0 mm Cancellous Screw
HCA5300-S	N/A	5.0 mm x 30.0 mm Cancellous Screw
HCA5325-S	N/A	5.0 mm x 32.5 mm Cancellous Screw
HCA5350-S	N/A	5.0 mm x 35.0 mm Cancellous Screw
HCA5375-S	N/A	5.0 mm x 37.5 mm Cancellous Screw
HCA5400-S	N/A	5.0 mm x 40.0 mm Cancellous Screw
HCA5450-S	N/A	5.0 mm x 45.0 mm Cancellous Screw
HCA5500-S	N/A	5.0 mm x 50.0 mm Cancellous Screw
HCA5550-S	N/A	5.0 mm x 55.0 mm Cancellous Screw
HCO3175-S	N/A	3.5 mm x 17.5 mm Cortical Screw
HCO3200-S	N/A	3.5 mm x 20.0 mm Cortical Screw
HCO3225-S	N/A	3.5 mm x 22.5 mm Cortical Screw
HCO3250-S	N/A	3.5 mm x 25.0 mm Cortical Screw
HCO3275-S	N/A	3.5 mm x 27.5 mm Cortical Screw
HCO3300-S	N/A	3.5 mm x 30.0 mm Cortical Screw
HCO3325-S	N/A	3.5 mm x 32.5 mm Cortical Screw
HCO3350-S	N/A	3.5 mm x 35.0 mm Cortical Screw
HCO3400-S	N/A	3.5 mm x 40.0 mm Cortical Screw
HCO3450-S	N/A	3.5 mm x 45.0 mm Cortical Screw
HCO3500-S	N/A	3.5 mm x 50.0 mm Cortical Screw
HCO3550-S	N/A	3.5 mm x 55.0 mm Cortical Screw

Intramedullary Nails Portfolio [continued]

Polarus® Plus Humeral Rod System

Sterile Part Number	Non-Sterile Part Number	Description
HCO35WA	N/A	3.5 mm Cortical Screw Washer
HCO35WA-S	N/A	3.5 mm Cortical Screw Washer
HCO3600-S	N/A	3.5 mm x 60.0 mm Cortical Screw
HR-0001-S	N/A	Polarus Cap Screw Assembly
HR-0050-S	N/A	Polarus Cap Screw Assembly
HR-0820-S	N/A	8 mm x 200 mm Polarus Plus Humeral Rod
HR-0822-S	N/A	8 mm x 220 mm Polarus Plus Humeral Rod
HR-0824-S	N/A	8 mm x 240 mm Polarus Plus Humeral Rod
HR-0826-S	N/A	8 mm x 260 mm Polarus Plus Humeral Rod
HR-0828-S	N/A	8 mm x 280 mm Polarus Plus Humeral Rod

Polarus® 3 Solution (Locking and Proximal Locking Nails)

Sterile Part Number	Non-Sterile Part Number	Description
4001-1015L-S	N/A	Polarus 3 Proximal Locking Nail 150 mm Left
4001-1015R-S	N/A	Polarus 3 Proximal Locking Nail 150 mm Right
4002-10200-S	N/A	Polarus 3 Locking Nail 200 mm
4002-10220-S	N/A	Polarus 3 Locking Nail 220 mm
4002-10240-S	N/A	Polarus 3 Locking Nail 240 mm
4002-10260-S	N/A	Polarus 3 Locking Nail 260 mm
4002-10280-S	N/A	Polarus 3 Locking Nail 280 mm
4004-10000-S	N/A	Polarus 3 Cap Screw 0 mm
4004-10002-S	N/A	Polarus 3 Cap Screw 2 mm
4004-10004-S	N/A	Polarus 3 Cap Screw 4 mm
4004-10006-S	N/A	Polarus 3 Cap Screw 6 mm
3011-43018-S	3011-43018	4.3 mm x 18 mm Low-Profile Hexalobe Screw
3011-43020-S	3011-43020	4.3 mm x 20 mm Low-Profile Hexalobe Screw
3011-43022-S	3011-43022	4.3 mm x 22 mm Low-Profile Hexalobe Screw
3011-43024-S	3011-43024	4.3 mm x 24 mm Low-Profile Hexalobe Screw
3011-43026-S	3011-43026	4.3 mm x 26 mm Low-Profile Hexalobe Screw

Intramedullary Nails Portfolio [continued]

Polarus® 3 Solution (Locking and Proximal Locking Nails)

Sterile Part Number	Non-Sterile Part Number	Description
3011-43028-S	3011-43028	4.3 mm x 28 mm Low-Profile Hexalobe Screw
3011-43030-S	3011-43030	4.3 mm x 30 mm Low-Profile Hexalobe Screw
3011-43032-S	3011-43032	4.3 mm x 32 mm Low-Profile Hexalobe Screw
3011-43034-S	3011-43034	4.3 mm x 34 mm Low-Profile Hexalobe Screw
3011-43036-S	3011-43036	4.3 mm x 36 mm Low-Profile Hexalobe Screw
3011-43038-S	3011-43038	4.3 mm x 38 mm Low-Profile Hexalobe Screw
3011-43040-S	3011-43040	4.3 mm x 40 mm Low-Profile Hexalobe Screw
3011-43042-S	3011-43042	4.3 mm x 42 mm Low-Profile Hexalobe Screw
3011-43044-S	3011-43044	4.3 mm x 44 mm Low-Profile Hexalobe Screw
3011-43046-S	3011-43046	4.3 mm x 46 mm Low-Profile Hexalobe Screw
3011-43048-S	3011-43048	4.3 mm x 48 mm Low-Profile Hexalobe Screw
3011-43050-S	3011-43050	4.3 mm x 50 mm Low-Profile Hexalobe Screw
3011-43052-S	3011-43052	4.3 mm x 52 mm Low-Profile Hexalobe Screw
3011-43054-S	3011-43054	4.3 mm x 54 mm Low-Profile Hexalobe Screw
3011-43056-S	3011-43056	4.3 mm x 56 mm Low-Profile Hexalobe Screw
3011-43058-S	3011-43058	4.3 mm x 58 mm Low-Profile Hexalobe Screw
3011-43060-S	3011-43060	4.3 mm x 60 mm Low-Profile Hexalobe Screw
3011-43062-S	3011-43062	4.3 mm x 62 mm Low-Profile Hexalobe Screw
3011-43064-S	3011-43064	4.3 mm x 64 mm Low-Profile Hexalobe Screw
3025-35018-S	3025-35018	3.5 x 18 mm Nonlocking Low-Profile Hexalobe Screw
3025-35020-S	3025-35020	3.5 x 20 mm Nonlocking Low-Profile Hexalobe Screw
3025-35022-S	3025-35022	3.5 x 22 mm Nonlocking Low-Profile Hexalobe Screw
3025-35024-S	3025-35024	3.5 x 24 mm Nonlocking Low-Profile Hexalobe Screw
3025-35026-S	3025-35026	3.5 x 26 mm Nonlocking Low-Profile Hexalobe Screw
3025-35028-S	3025-35028	3.5 x 28 mm Nonlocking Low-Profile Hexalobe Screw
3025-35030-S	3025-35030	3.5 x 30 mm Nonlocking Low-Profile Hexalobe Screw
3025-35032-S	3025-35032	3.5 x 32 mm Nonlocking Low-Profile Hexalobe Screw
3025-35034-S	3025-35034	3.5 x 34 mm Nonlocking Low-Profile Hexalobe Screw
3025-35036-S	3025-35036	3.5 x 36 mm Nonlocking Low-Profile Hexalobe Screw

Intramedullary Nails Portfolio [continued]

Polarus® 3 Solution (Locking and Proximal Locking Nails)

Sterile Part Number	Non-Sterile Part Number	Description
3025-35038-S	3025-35038	3.5 x 38 mm Nonlocking Low-Profile Hexalobe Screw
3025-35040-S	3025-35040	3.5 x 40 mm Nonlocking Low-Profile Hexalobe Screw
3025-35042-S	3025-35042	3.5 x 42 mm Nonlocking Low-Profile Hexalobe Screw
3025-35044-S	3025-35044	3.5 x 44 mm Nonlocking Low-Profile Hexalobe Screw
3025-35046-S	3025-35046	3.5 x 46 mm Nonlocking Low-Profile Hexalobe Screw
3025-35048-S	3025-35048	3.5 x 48 mm Nonlocking Low-Profile Hexalobe Screw
3025-35050-S	3025-35050	3.5 x 50 mm Nonlocking Low-Profile Hexalobe Screw
3025-35052-S	3025-35052	3.5 x 52 mm Nonlocking Low-Profile Hexalobe Screw
3025-35054-S	3025-35054	3.5 x 54 mm Nonlocking Low-Profile Hexalobe Screw
3025-35056-S	3025-35056	3.5 x 56 mm Nonlocking Low-Profile Hexalobe Screw
3025-35058-S	3025-35058	3.5 x 58 mm Nonlocking Low-Profile Hexalobe Screw
3025-35060-S	3025-35060	3.5 x 60 mm Nonlocking Low-Profile Hexalobe Screw
3025-35062-S	3025-35062	3.5 x 62 mm Nonlocking Low-Profile Hexalobe Screw
3025-35064-S	3025-35064	3.5 x 64 mm Nonlocking Low-Profile Hexalobe Screw

Fibula Rod System

Sterile Part Number	Non-Sterile Part Number	Description
40-0026-S	N/A	3.0 mm x 110 mm Fibula Rod
40-0027-S	N/A	3.0 mm x 145 mm Fibula Rod
40-0028-S	N/A	3.0 mm x 180 mm Fibula Rod
40-0029-S	N/A	3.6 mm x 110 mm Fibula Rod
40-0030-S	N/A	3.6 mm x 145 mm Fibula Rod
40-0031-S	N/A	3.6 mm x 180 mm Fibula Rod
CO-3080-S	CO-3080	3.5 mm x 8.0 mm Cortical Screw
CO-3100-S	CO-3100	3.5 mm x 10.0 mm Cortical Screw
CO-3120-S	CO-3120	3.5 mm x 12.0 mm Cortical Screw
CO-3140-S	CO-3140	3.5 mm x 14.0 mm Cortical Screw
CO-3160-S	CO-3160	3.5 mm x 16.0 mm Cortical Screw
CO-3180-S	CO-3180	3.5 mm x 18.0 mm Cortical Screw

Intramedullary Nails Portfolio [continued]

Fibula Rod System

Sterile Part Number	Non-Sterile Part Number	Description
CO-3200-S	CO-3200	3.5 mm x 20.0 mm Cortical Screw
CO-3220-S	CO-3220	3.5 mm x 22.0 mm Cortical Screw
CO-3400-S	CO-3400	3.5 mm x 40.0 mm Cortical Screw
CO-3450-S	CO-3450	3.5 mm x 45.0 mm Cortical Screw
CO-3500-S	CO-3500	3.5 mm x 50.0 mm Cortical Screw
CO-3550-S	CO-3550	3.5 mm x 55.0 mm Cortical Screw
CO-3600-S	CO-3600	3.5 mm x 60.0 mm Cortical Screw
CO-3650-S	CO-3650	3.5 mm x 65.0 mm Cortical Screw

Forearm Fracture Solutions (Radius & Ulna Rods)

Sterile Part Number	Non-Sterile Part Number	Description
CO-3125-S	N/A	3.5 mm x 12.5 mm Cortical Screw
CO-3150-S	N/A	3.5 mm x 15.0 mm Cortical Screw
CO-3175-S	N/A	3.5 mm x 17.5 mm Cortical Screw
CO-3200-S	CO-3200	3.5 mm x 20.0 mm Cortical Screw
CO-3225-S	CO-3225	3.5 mm x 22.5 mm Cortical Screw
CO-3250-S	CO-3250	3.5 mm x 25.0 mm Cortical Screw
CO-3275-S	CO-3275	3.5 mm x 27.5 mm Cortical Screw
CO-3300-S	CO-3300	3.5 mm x 30.0 mm Cortical Screw
CO-3325-S	CO-3325	3.5 mm x 32.5 mm Cortical Screw
RL-3019-S	N/A	3.0 mm x 190 mm Radius Rod, Left
RL-3021-S	N/A	3.0 mm x 210 mm Radius Rod, Left
RL-3023-S	N/A	3.0 mm x 230 mm Radius Rod, Left
RL-3619-S	N/A	3.6 mm x 190 mm Radius Rod, Left
RL-3621-S	N/A	3.6 mm x 210 mm Radius Rod, Left
RL-3623-S	N/A	3.6 mm x 230 mm Radius Rod, Left
RR-3019-S	N/A	3.0 mm x 190 mm Radius Rod, Right
RR-3021-S	N/A	3.0 mm x 210 mm Radius Rod, Right
RR-3023-S	N/A	3.0 mm x 230 mm Radius Rod, Right

Intramedullary Nails Portfolio [continued]

Forearm Fracture Solutions (Radius & Ulna Rods)

Sterile Part Number	Non-Sterile Part Number	Description
RR-3619-S	N/A	3.6 mm x 190 mm Radius Rod, Right
RR-3621-S	N/A	3.6 mm x 210 mm Radius Rod, Right
RR-3623-S	N/A	3.6 mm x 230 mm Radius Rod, Right
UL-3021-S	N/A	3.0 mm x 210 mm Ulna Rod
UL-3023-S	N/A	3.0 mm x 230 mm Ulna Rod
UL-3025-S	N/A	3.0 mm x 250 mm Ulna Rod
UL-3027-S	N/A	3.0 mm x 270 mm Ulna Rod
UL-3621-S	N/A	3.6 mm x 210 mm Ulna Rod
UL-3623-S	N/A	3.6 mm x 230 mm Ulna Rod
UL-3625-S	N/A	3.6 mm x 250 mm Ulna Rod
UL-3627-S	N/A	3.6 mm x 270 mm Ulna Rod

Fibula Nail 2 System

Sterile Part Number	Non-Sterile Part Number	Description
30-0255-S	30-0255	3.5 mm x 8 mm Nonlocking Hexalobe Screw
30-0256-S	30-0256	3.5 mm x 10 mm Nonlocking Hexalobe Screw
30-0257-S	30-0257	3.5 mm x 12 mm Nonlocking Hexalobe Screw
30-0258-S	30-0258	3.5 mm x 14 mm Nonlocking Hexalobe Screw
30-0259-S	30-0259	3.5 mm x 16 mm Nonlocking Hexalobe Screw
30-0260-S	30-0260	3.5 mm x 18 mm Nonlocking Hexalobe Screw
30-0261-S	30-0261	3.5 mm x 20 mm Nonlocking Hexalobe Screw
30-0262-S	30-0262	3.5 mm x 22 mm Nonlocking Hexalobe Screw
30-0263-S	30-0263	3.5 mm x 24 mm Nonlocking Hexalobe Screw
30-0264-S	30-0264	3.5 mm x 26 mm Nonlocking Hexalobe Screw
30-0265-S	30-0265	3.5 mm x 28 mm Nonlocking Hexalobe Screw
30-0266-S	30-0266	3.5 mm x 30 mm Nonlocking Hexalobe Screw
30-0267-S	30-0267	3.5 mm x 32 mm Nonlocking Hexalobe Screw
30-0268-S	30-0268	3.5 mm x 34 mm Nonlocking Hexalobe Screw
30-0269-S	30-0269	3.5 mm x 36 mm Nonlocking Hexalobe Screw

Intramedullary Nails Portfolio [continued]

Fibula Nail 2 System

Sterile Part Number	Non-Sterile Part Number	Description
30-0270-S	30-0270	3.5 mm x 38 mm Nonlocking Hexalobe Screw
30-0271-S	30-0271	3.5 mm x 40 mm Nonlocking Hexalobe Screw
30-0272-S	30-0272	3.5 mm x 45 mm Nonlocking Hexalobe Screw
30-0273-S	30-0273	3.5 mm x 50 mm Nonlocking Hexalobe Screw
30-0274-S	30-0274	3.5 mm x 55 mm Nonlocking Hexalobe Screw
30-0275-S	30-0275	3.5 mm x 60 mm Nonlocking Hexalobe Screw
30-0276-S	30-0276	3.5 mm x 65 mm Nonlocking Hexalobe Screw
N/A	3017-25006	Far End Locking Set Screw, 6 mm
N/A	3017-25007	Far End Locking Set Screw, 7 mm
N/A	3017-25008	Far End Locking Set Screw, 8 mm
N/A	3017-25009	Far End Locking Set Screw, 9 mm
N/A	3017-25010	Far End Locking Set Screw, 10 mm
N/A	3017-25011	Far End Locking Set Screw, 11 mm
N/A	3017-25012	Far End Locking Set Screw, 12 mm
N/A	3017-25013	Far End Locking Set Screw, 13 mm
N/A	3017-25014	Far End Locking Set Screw, 14 mm
N/A	3017-25015	Far End Locking Set Screw, 15 mm
N/A	3017-25016	Far End Locking Set Screw, 16 mm
N/A	3017-65006	Far End Locking Bushing, 6 mm
N/A	3017-65007	Far End Locking Bushing, 7 mm
N/A	3017-65008	Far End Locking Bushing, 8 mm
N/A	3017-65009	Far End Locking Bushing, 9 mm
N/A	3017-65010	Far End Locking Bushing, 10 mm
N/A	3017-65011	Far End Locking Bushing, 11 mm
N/A	3017-65012	Far End Locking Bushing, 12 mm
N/A	3017-65013	Far End Locking Bushing, 13 mm
N/A	3017-65014	Far End Locking Bushing, 14 mm
N/A	3017-65015	Far End Locking Bushing, 15 mm
N/A	3017-65016	Far End Locking Bushing, 16 mm

Intramedullary Nails Portfolio [continued]

Fibula Nail 2 System

Sterile Part Number	Non-Sterile Part Number	Description
N/A	3018-47012	3.5 mm x 12 mm Headless Hexalobe Screw
N/A	3018-47014	3.5 mm x 14 mm Headless Hexalobe Screw
N/A	3018-47016	3.5 mm x 16 mm Headless Hexalobe Screw
N/A	3018-47018	3.5 mm x 18 mm Headless Hexalobe Screw
N/A	3018-47020	3.5 mm x 20 mm Headless Hexalobe Screw
N/A	3018-47022	3.5 mm x 22 mm Headless Hexalobe Screw
N/A	3018-47024	3.5 mm x 24 mm Headless Hexalobe Screw
N/A	3018-47026	3.5 mm x 26 mm Headless Hexalobe Screw
N/A	3018-47028	3.5 mm x 28 mm Headless Hexalobe Screw
N/A	3018-47030	3.5 mm x 30 mm Headless Hexalobe Screw
N/A	3018-47032	3.5 mm x 32 mm Headless Hexalobe Screw
N/A	3018-47034	3.5 mm x 34 mm Headless Hexalobe Screw
N/A	3018-47036	3.5 mm x 36 mm Headless Hexalobe Screw
N/A	3018-47038	3.5 mm x 38 mm Headless Hexalobe Screw
N/A	3018-47040	3.5 mm x 40 mm Headless Hexalobe Screw
N/A	3018-47045	3.5 mm x 45 mm Headless Hexalobe Screw
N/A	3018-47050	3.5 mm x 50 mm Headless Hexalobe Screw
N/A	3018-47055	3.5 mm x 55 mm Headless Hexalobe Screw
N/A	3018-47060	3.5 mm x 60 mm Headless Hexalobe Screw
N/A	3018-47065	3.5 mm x 65 mm Headless Hexalobe Screw
4010-2611N-S	N/A	2.6 mm x 110 mm Fibula Nail
4010-2614N-S	N/A	2.6 mm x 145 mm Fibula Nail
4010-2619N-S	N/A	2.6 mm x 190 mm Fibula Nail
4010-2627N-S	N/A	2.6 mm x 270 mm Fibula Nail
4010-3011N-S	N/A	3.0 mm x 110 mm Fibula Nail
4010-3014N-S	N/A	3.0 mm x 145 mm Fibula Nail
4010-3019N-S	N/A	3.0 mm x 190 mm Fibula Nail
4010-3027N-S	N/A	3.0 mm x 270 mm Fibula Nail
4010-3611N-S	N/A	3.6 mm x 110 mm Fibula Nail

Intramedullary Nails Portfolio [continued]

Fibula Nail 2 System

Sterile Part Number	Non-Sterile Part Number	Description
4010-3614N-S	N/A	3.6 mm x 145 mm Fibula Nail
4010-3619N-S	N/A	3.6 mm x 190 mm Fibula Nail
4010-3627N-S	N/A	3.6 mm x 270 mm Fibula Nail
N/A	4014-0600	Fibula and Forearm Nail +0.4 mm End Cap
N/A	4014-0705	Fibula and Forearm Nail +5 mm End Cap
N/A	4014-0710	Fibula and Forearm Nail +10 mm End Cap
N/A	4014-0715	Fibula and Forearm Nail +15 mm End Cap

Ulna Nail 2 System

Sterile Part Number	Non-Sterile Part Number	Description
30-0255-S	30-0255	3.5 mm x 8 mm Nonlocking Hexalobe Screw
30-0256-S	30-0256	3.5 mm x 10 mm Nonlocking Hexalobe Screw
30-0257-S	30-0257	3.5 mm x 12 mm Nonlocking Hexalobe Screw
30-0258-S	30-0258	3.5 mm x 14 mm Nonlocking Hexalobe Screw
30-0259-S	30-0259	3.5 mm x 16 mm Nonlocking Hexalobe Screw
30-0260-S	30-0260	3.5 mm x 18 mm Nonlocking Hexalobe Screw
30-0261-S	30-0261	3.5 mm x 20 mm Nonlocking Hexalobe Screw
30-0262-S	30-0262	3.5 mm x 22 mm Nonlocking Hexalobe Screw
30-0263-S	30-0263	3.5 mm x 24 mm Nonlocking Hexalobe Screw
30-0264-S	30-0264	3.5 mm x 26 mm Nonlocking Hexalobe Screw
30-0265-S	30-0265	3.5 mm x 28 mm Nonlocking Hexalobe Screw
30-0266-S	30-0266	3.5 mm x 30 mm Nonlocking Hexalobe Screw
30-0267-S	30-0267	3.5 mm x 32 mm Nonlocking Hexalobe Screw
30-0268-S	30-0268	3.5 mm x 34 mm Nonlocking Hexalobe Screw
30-0269-S	30-0269	3.5 mm x 36 mm Nonlocking Hexalobe Screw
30-0270-S	30-0270	3.5 mm x 38 mm Nonlocking Hexalobe Screw
30-0271-S	30-0271	3.5 mm x 40 mm Nonlocking Hexalobe Screw
30-0272-S	30-0272	3.5 mm x 45 mm Nonlocking Hexalobe Screw
30-0273-S	30-0273	3.5 mm x 50 mm Nonlocking Hexalobe Screw

Intramedullary Nails Portfolio [continued]

Ulna Nail 2 System

Sterile Part Number	Non-Sterile Part Number	Description
30-0274-S	30-0274	3.5 mm x 55 mm Nonlocking Hexalobe Screw
30-0275-S	30-0275	3.5 mm x 60 mm Nonlocking Hexalobe Screw
30-0276-S	30-0276	3.5 mm x 65 mm Nonlocking Hexalobe Screw
N/A	3017-25006	Far End Locking Set Screw, 6 mm
N/A	3017-25007	Far End Locking Set Screw, 7 mm
N/A	3017-25008	Far End Locking Set Screw, 8 mm
N/A	3017-25009	Far End Locking Set Screw, 9 mm
N/A	3017-25010	Far End Locking Set Screw, 10 mm
N/A	3017-25011	Far End Locking Set Screw, 11 mm
N/A	3017-25012	Far End Locking Set Screw, 12 mm
N/A	3017-25013	Far End Locking Set Screw, 13 mm
N/A	3017-25014	Far End Locking Set Screw, 14 mm
N/A	3017-25015	Far End Locking Set Screw, 15 mm
N/A	3017-25016	Far End Locking Set Screw, 16 mm
N/A	3017-65006	Far End Locking Bushing, 6 mm
N/A	3017-65007	Far End Locking Bushing, 7 mm
N/A	3017-65008	Far End Locking Bushing, 8 mm
N/A	3017-65009	Far End Locking Bushing, 9 mm
N/A	3017-65010	Far End Locking Bushing, 10 mm
N/A	3017-65011	Far End Locking Bushing, 11 mm
N/A	3017-65012	Far End Locking Bushing, 12 mm
N/A	3017-65013	Far End Locking Bushing, 13 mm
N/A	3017-65014	Far End Locking Bushing, 14 mm
N/A	3017-65015	Far End Locking Bushing, 15 mm
N/A	3017-65016	Far End Locking Bushing, 16 mm
N/A	3018-47012	3.5 mm x 12 mm Headless Hexalobe Screw
N/A	3018-47014	3.5 mm x 14 mm Headless Hexalobe Screw
N/A	3018-47016	3.5 mm x 16 mm Headless Hexalobe Screw
N/A	3018-47018	3.5 mm x 18 mm Headless Hexalobe Screw

Intramedullary Nails Portfolio [continued]

Ulna Nail 2 System

Sterile Part Number	Non-Sterile Part Number	Description
N/A	3018-47020	3.5 mm x 20 mm Headless Hexalobe Screw
N/A	3018-47022	3.5 mm x 22 mm Headless Hexalobe Screw
N/A	3018-47024	3.5 mm x 24 mm Headless Hexalobe Screw
N/A	3018-47026	3.5 mm x 26 mm Headless Hexalobe Screw
N/A	3018-47028	3.5 mm x 28 mm Headless Hexalobe Screw
N/A	3018-47030	3.5 mm x 30 mm Headless Hexalobe Screw
N/A	3018-47032	3.5 mm x 32 mm Headless Hexalobe Screw
N/A	3018-47034	3.5 mm x 34 mm Headless Hexalobe Screw
N/A	3018-47036	3.5 mm x 36 mm Headless Hexalobe Screw
N/A	3018-47038	3.5 mm x 38 mm Headless Hexalobe Screw
N/A	3018-47040	3.5 mm x 40 mm Headless Hexalobe Screw
N/A	3018-47045	3.5 mm x 45 mm Headless Hexalobe Screw
N/A	3018-47050	3.5 mm x 50 mm Headless Hexalobe Screw
N/A	3018-47055	3.5 mm x 55 mm Headless Hexalobe Screw
N/A	3018-47060	3.5 mm x 60 mm Headless Hexalobe Screw
N/A	3018-47065	3.5 mm x 65 mm Headless Hexalobe Screw
4011-3012N-S	N/A	3.0 mm x 120 mm Short Ulna Nail
4011-3017N-S	N/A	3.0 mm x 170 mm Ulna Nail
4011-3019N-S	N/A	3.0 mm x 190 mm Ulna Nail
4011-3021N-S	N/A	3.0 mm x 210 mm Ulna Nail
4011-3023N-S	N/A	3.0 mm x 230 mm Ulna Nail
4011-3025N-S	N/A	3.0 mm x 250 mm Ulna Nail
4011-3027N-S	N/A	3.0 mm x 270 mm Ulna Nail
4011-3612N-S	N/A	3.6 mm x 120 mm Short Ulna Nail
4011-3617N-S	N/A	3.6 mm x 170 mm Ulna Nail
4011-3619N-S	N/A	3.6 mm x 190 mm Ulna Nail
4011-3621N-S	N/A	3.6 mm x 210 mm Ulna Nail
4011-3623N-S	N/A	3.6 mm x 230 mm Ulna Nail
4011-3625N-S	N/A	3.6 mm x 250 mm Ulna Nail

Intramedullary Nails Portfolio [continued]

Ulna Nail 2 System

Sterile Part Number	Non-Sterile Part Number	Description
4011-3627N-S	N/A	3.6 mm x 270 mm Ulna Nail
4011-4012N-S	N/A	4.0 mm x 120 mm Short Ulna Nail
4011-4017N-S	N/A	4.0 mm x 170 mm Ulna Nail
4011-4019N-S	N/A	4.0 mm x 190 mm Ulna Nail
4011-4021N-S	N/A	4.0 mm x 210 mm Ulna Nail
4011-4023N-S	N/A	4.0 mm x 230 mm Ulna Nail
4011-4025N-S	N/A	4.0 mm x 250 mm Ulna Nail
4011-4027N-S	N/A	4.0 mm x 270 mm Ulna Nail
N/A	4014-0600	Fibula and Forearm Nail +0.4 mm End Cap
N/A	4014-0705	Fibula and Forearm Nail +5 mm End Cap
N/A	4014-0710	Fibula and Forearm Nail +10 mm End Cap
N/A	4014-0715	Fibula and Forearm Nail +15 mm End Cap

Conditions for safe scanning in the MRI

A patient with the nail/rod, screw, washer, end cap, bushing, and set screw implants from the Polarus® Humeral Rod, Polarus Plus Humeral Rod, Polarus 3 Solution, Fibula Rod, Forearm Fracture Solutions, Fibula Nail 2, and Ulna Nail 2 Systems may be safely scanned under the following conditions. Failure to follow these conditions may result in injury to the patient.

Nominal values(s) of Static Magnetic Field	1.5 T or 3 T
Maximum Spatial Field Gradient	19 T/m (1900 gauss/cm)
RF Excitation	Circularly Polarized (CP)
RF Transmit Coil Type	Whole body transmit coil, Head RF transmit-receive coil
Maximum Whole Body SAR at 1.5 T and 3.0 T	2.0 W/kg
Limits on Scan Duration	2.0 W/kg whole body average SAR for 6 minutes of continuous RF (a sequence or back to back series/scan without breaks)
MR Image Artifact	The presence of this implant may produce an image artifact of 24 mm.
If more information about a specific parameter is not included, there are no conditions associated with that parameter.	

Hand and Wrist Portfolio

- ▶ Acu-Loc Wrist Plating System
- ▶ Acu-Loc 2 Wrist Plating System
- ▶ Acu-Loc Wrist Spanning Plate
- ▶ Forearm Fracture Solutions (plates)
- ▶ Hand Fracture System
- ▶ Modular Hand System
- ▶ Osteotomy System
- ▶ Total Wrist Fusion Plating System

Acu-Loc® Wrist Plating System

Sterile Part Number	Non-Sterile Part Number	Description
70-0045-S	70-0045	Acu-Loc Volar Distal Ulna Plate Left Standard
70-0046-S	70-0046	Acu-Loc Volar Distal Ulna Plate Right Standard
70-0047-S	70-0047	Acu-Loc Volar Distal Ulna Plate Left Long
70-0048-S	70-0048	Acu-Loc Volar Distal Ulna Plate Right Long
70-0055-S	70-0055	Acu-Loc Dorsal Plate Standard, Left
70-0056-S	70-0056	Acu-Loc Dorsal Plate Standard, Right
70-0057-S	70-0057	Acu-Loc Dorsal Plate Narrow Left
70-0058-S	70-0058	Acu-Loc Dorsal Plate Narrow Right
70-0063-S	70-0063	Acu-Loc Extra Articular Standard
70-0064-S	70-0064	Acu-Loc Extra Articular Narrow
CO-3100-S	CO-3100	3.5 mm x 10.0 mm Cortical Screw
CO-3120-S	CO-3120	3.5 mm x 12.0 mm Cortical Screw
CO-3140-S	CO-3140	3.5 mm x 14.0 mm Cortical Screw
CO-3160-S	CO-3160	3.5 mm x 16.0 mm Cortical Screw
CO-3180-S	CO-3180	3.5 mm x 18.0 mm Cortical Screw
COL-3080-S	COL-3080	3.5 mm x 8.0 mm Locking Cortical Screw
COL-3100-S	COL-3100	3.5 mm x 10.0 mm Locking Cortical Screw
COL-3120-S	COL-3120	3.5 mm x 12.0 mm Locking Cortical Screw
COL-3140-S	COL-3140	3.5 mm x 14.0 mm Locking Cortical Screw
COL-3160-S	COL-3160	3.5 mm x 16.0 mm Locking Cortical Screw
COL-3180-S	COL-3180	3.5 mm x 18.0 mm Locking Cortical Screw
CO-N2308-S	CO-N2308	2.3 mm x 8 mm Non-Toggling Cortical Screw
CO-N2310-S	CO-N2310	2.3 mm x 10 mm Non-Toggling Cortical Screw
CO-N2312-S	CO-N2312	2.3 mm x 12 mm Non-Toggling Cortical Screw
CO-N2314-S	CO-N2314	2.3 mm x 14 mm Non-Toggling Cortical Screw
CO-N2316-S	CO-N2316	2.3 mm x 16 mm Non-Toggling Cortical Screw
CO-N2318-S	CO-N2318	2.3 mm x 18 mm Non-Toggling Cortical Screw

Hand and Wrist Portfolio [continued]

Acu-Loc® Wrist Plating System

Sterile Part Number	Non-Sterile Part Number	Description
CO-N2320-S	CO-N2320	2.3 mm x 20 mm Non-Toggling Cortical Screw
CO-N2322-S	CO-N2322	2.3 mm x 22 mm Non-Toggling Cortical Screw
CO-N2324-S	CO-N2324	2.3 mm x 24 mm Non-Toggling Cortical Screw
CO-N2326-S	CO-N2326	2.3 mm x 26 mm Non-Toggling Cortical Screw
CO-N2328-S	CO-N2328	2.3 mm x 28 mm Non-Toggling Cortical Screw
CO-N2330-S	CO-N2330	2.3 mm x 30 mm Non-Toggling Cortical Screw
CO-N2332-S	CO-N2332	2.3 mm x 32 mm Non-Toggling Cortical Screw
CO-S2308-S	CO-S2308	2.3 mm x 8 mm Locking Cortical Peg
CO-S2310-S	CO-S2310	2.3 mm x 10 mm Locking Cortical Peg
CO-S2312-S	CO-S2312	2.3 mm x 12 mm Locking Cortical Peg
CO-S2314-S	CO-S2314	2.3 mm x 14 mm Locking Cortical Peg
CO-S2316-S	CO-S2316	2.3 mm x 16 mm Locking Cortical Peg
CO-S2318-S	CO-S2318	2.3 mm x 18 mm Locking Cortical Peg
CO-S2320-S	CO-S2320	2.3 mm x 20 mm Locking Cortical Peg
CO-S2322-S	CO-S2322	2.3 mm x 22 mm Locking Cortical Peg
CO-S2324-S	CO-S2324	2.3 mm x 24 mm Locking Cortical Peg
CO-S2326-S	CO-S2326	2.3 mm x 26 mm Locking Cortical Peg
CO-S2328-S	CO-S2328	2.3 mm x 28 mm Locking Cortical Peg
CO-T2308-S	CO-T2308	2.3 mm x 8 mm Locking Cortical Screw
CO-T2310-S	CO-T2310	2.3 mm x 10 mm Locking Cortical Screw
CO-T2312-S	CO-T2312	2.3 mm x 12 mm Locking Cortical Screw
CO-T2314-S	CO-T2314	2.3 mm x 14 mm Locking Cortical Screw
CO-T2316-S	CO-T2316	2.3 mm x 16 mm Locking Cortical Screw
CO-T2318-S	CO-T2318	2.3 mm x 18 mm Locking Cortical Screw
CO-T2320-S	CO-T2320	2.3 mm x 20 mm Locking Cortical Screw
CO-T2322-S	CO-T2322	2.3 mm x 22 mm Locking Cortical Screw
CO-T2324-S	CO-T2324	2.3 mm x 24 mm Locking Cortical Screw
CO-T2326-S	CO-T2326	2.3 mm x 26 mm Locking Cortical Screw
CO-T2328-S	CO-T2328	2.3 mm x 28 mm Locking Cortical Screw

Hand and Wrist Portfolio [continued]

Acu-Loc® Wrist Plating System

Sterile Part Number	Non-Sterile Part Number	Description
PL-DR30L-S	PL-DR30L	Acu-Loc Volar Distal Radius Plate, Narrow, Left
PL-DR30R-S	PL-DR30R	Acu-Loc Volar Distal Radius Plate, Narrow, Right
PL-DR50L-S	PL-DR50L	Acu-Loc Volar Distal Radius Plate, Standard, Left
PL-DR50R-S	PL-DR50R	Acu-Loc Volar Distal Radius Plate, Standard, Right
PL-DR60L-S	PL-DR60L	Acu-Loc Volar Distal Radius Plate, Long, Left
PL-DR60R-S	PL-DR60R	Acu-Loc Volar Distal Radius Plate, Long, Right
PL-DR65L-S	PL-DR65L	Acu-Loc Volar Distal Radius Plate, X-Long, Left
PL-DR65R-S	PL-DR65R	Acu-Loc Volar Distal Radius Plate, X-Long, Right
PL-DR70L-S	PL-DR70L	Acu-Loc Volar Distal Radius Plate, Wide, Left
PL-DR70R-S	PL-DR70R	Acu-Loc Volar Distal Radius Plate, Wide, Right

Acu-Loc® 2 Wrist Plating System

Sterile Part Number	Non-Sterile Part Number	Description
30-0093-S	30-0093	Acu-Loc 2 Volar Distal Radius Plate Extension Link Screw
30-0100-S	30-0100	Acu-Loc 2 Volar Distal Radius Hexalobe Ext Link Screw
30-0218-S	N/A	3.5 mm x 9 mm Locking Hexalobe Screw
30-0219-S	N/A	3.5 mm x 11 mm Locking Hexalobe Screw
30-0220-S	N/A	3.5 mm x 13 mm Locking Hexalobe Screw
30-0221-S	N/A	3.5 mm x 15 mm Locking Hexalobe Screw
30-0222-S	N/A	3.5 mm x 17 mm Locking Hexalobe Screw
30-0223-S	N/A	3.5 mm x 19 mm Locking Hexalobe Screw
30-0224-S	N/A	3.5 mm x 9 mm Nonlocking Hexalobe Screw
30-0225-S	N/A	3.5 mm x 11 mm Nonlocking Hexalobe Screw
30-0226-S	N/A	3.5 mm x 13 mm Nonlocking Hexalobe Screw
30-0227-S	N/A	3.5 mm x 15 mm Nonlocking Hexalobe Screw
30-0228-S	N/A	3.5 mm x 17 mm Nonlocking Hexalobe Screw
30-0229-S	N/A	3.5 mm x 19 mm Nonlocking Hexalobe Screw
30-0232-S	30-0232	3.5 mm x 8 mm Locking Hexalobe Screw

Hand and Wrist Portfolio [continued]

Acu-Loc® 2 Wrist Plating System

Sterile Part Number	Non-Sterile Part Number	Description
30-0233-S	30-0233	3.5 mm x 10 mm Locking Hexalobe Screw
30-0234-S	30-0234	3.5 mm x 12 mm Locking Hexalobe Screw
30-0235-S	30-0235	3.5 mm x 14 mm Locking Hexalobe Screw
30-0236-S	30-0236	3.5 mm x 16 mm Locking Hexalobe Screw
30-0237-S	30-0237	3.5 mm x 18 mm Locking Hexalobe Screw
30-0256-S	30-0256	3.5 mm x 10 mm Nonlocking Hexalobe Screw
30-0257-S	30-0257	3.5 mm x 12 mm Nonlocking Hexalobe Screw
30-0258-S	30-0258	3.5 mm x 14 mm Nonlocking Hexalobe Screw
30-0259-S	30-0259	3.5 mm x 16 mm Nonlocking Hexalobe Screw
30-0260-S	30-0260	3.5 mm x 18 mm Nonlocking Hexalobe Screw
30-0370-S	30-0370	Frag-Loc® Compression Sleeve
30-0371-S	30-0371	Frag-Loc Compression Screw
30-0372-S	30-0372	Frag-Loc Compression Screw, Long
30-2314-S	30-2314	2.3 mm x 14 mm Locking Variable Angle Screw
30-2316-S	30-2316	2.3 mm x 16 mm Locking Variable Angle Screw
30-2318-S	30-2318	2.3 mm x 18 mm Locking Variable Angle Screw
30-2320-S	30-2320	2.3 mm x 20 mm Locking Variable Angle Screw
30-2322-S	30-2322	2.3 mm x 22 mm Locking Variable Angle Screw
30-2324-S	30-2324	2.3 mm x 24 mm Locking Variable Angle Screw
30-2326-S	30-2326	2.3 mm x 26 mm Locking Variable Angle Screw
30-2328-S	30-2328	2.3 mm x 28 mm Locking Variable Angle Screw
N/A	3040-23008	2.7 x 8 mm Locking Low-Profile Hexalobe Screw
N/A	3040-23009	2.7 x 9 mm Locking Low-Profile Hexalobe Screw
N/A	3040-23010	2.7 x 10 mm Locking Low-Profile Hexalobe Screw
N/A	3040-23011	2.7 x 11 mm Locking Low-Profile Hexalobe Screw
N/A	3040-23012	2.7 x 12 mm Locking Low-Profile Hexalobe Screw
N/A	3040-23013	2.7 x 13 mm Locking Low-Profile Hexalobe Screw
N/A	3040-23014	2.7 x 14 mm Locking Low-Profile Hexalobe Screw
N/A	3040-23015	2.7 x 15 mm Locking Low-Profile Hexalobe Screw

Hand and Wrist Portfolio [continued]

Acu-Loc® 2 Wrist Plating System

Sterile Part Number	Non-Sterile Part Number	Description
N/A	3040-23016	2.7 x 16 mm Locking Low-Profile Hexalobe Screw
N/A	3040-23017	2.7 x 17 mm Locking Low-Profile Hexalobe Screw
N/A	3040-23018	2.7 x 18 mm Locking Low-Profile Hexalobe Screw
N/A	3040-23019	2.7 x 19 mm Locking Low-Profile Hexalobe Screw
N/A	3041-23009	2.7 x 9 mm Nonlocking Low-Profile Hexalobe Screw
N/A	3041-23010	2.7 x 10 mm Nonlocking Low-Profile Hexalobe Screw
N/A	3041-23011	2.7 x 11 mm Nonlocking Low-Profile Hexalobe Screw
N/A	3041-23012	2.7 x 12 mm Nonlocking Low-Profile Hexalobe Screw
N/A	3041-23013	2.7 x 13 mm Nonlocking Low-Profile Hexalobe Screw
N/A	3041-23014	2.7 x 14 mm Nonlocking Low-Profile Hexalobe Screw
N/A	3041-23015	2.7 x 15 mm Nonlocking Low-Profile Hexalobe Screw
N/A	3041-23016	2.7 x 16 mm Nonlocking Low-Profile Hexalobe Screw
N/A	3041-23017	2.7 x 17 mm Nonlocking Low-Profile Hexalobe Screw
N/A	3041-23018	2.7 x 18 mm Nonlocking Low-Profile Hexalobe Screw
N/A	3041-23019	2.7 x 19 mm Nonlocking Low-Profile Hexalobe Screw
70-0045-S	70-0045	Acu-Loc Volar Distal Ulna Plate Left Standard
70-0046-S	70-0046	Acu-Loc Volar Distal Ulna Plate Right Standard
70-0047-S	70-0047	Acu-Loc Volar Distal Ulna Plate Left Long
70-0048-S	70-0048	Acu-Loc Volar Distal Ulna Plate Right Long
70-0055-S	70-0055	Acu-Loc Dorsal Plate Standard, Left
70-0056-S	70-0056	Acu-Loc Dorsal Plate Standard, Right
70-0057-S	70-0057	Acu-Loc Dorsal Plate Narrow Left
70-0058-S	70-0058	Acu-Loc Dorsal Plate Narrow Right
70-0063-S	70-0063	Acu-Loc Extra Articular Standard
70-0064-S	70-0064	Acu-Loc Extra Articular Narrow
70-0331-S	70-0331	Divergent Radial Styloid Plate
70-0334-S	70-0334	Volar Lunate Suture Plate
70-0335-S	70-0335	Dorsal Rim Buttress Plate, Right
70-0336-S	70-0336	Dorsal Rim Buttress Plate, Left

Hand and Wrist Portfolio [continued]

Acu-Loc® 2 Wrist Plating System

Sterile Part Number	Non-Sterile Part Number	Description
70-0337-S	70-0337	Dorsal Lunate Plate, Right
70-0338-S	70-0338	Dorsal Lunate Plate, Left
70-0350-S	70-0350	Acu-Loc 2 Volar Distal Radius Proximal Plate, Standard, Left
70-0351-S	70-0351	Acu-Loc 2 Volar Distal Radius Proximal Plate, Standard, Right
70-0352-S	70-0352	Acu-Loc 2 Volar Distal Radius Proximal Plate, Narrow, Left
70-0353-S	70-0353	Acu-Loc 2 Volar Distal Radius Proximal Plate, Narrow, Right
70-0354-S	70-0354	Acu-Loc 2 Volar Distal Radius Proximal Plate, Wide, Left
70-0355-S	70-0355	Acu-Loc 2 Volar Distal Radius Proximal Plate, Wide, Right
70-0356-S	70-0356	Acu-Loc 2 Volar Distal Radius Plate Standard, Left
70-0357-S	70-0357	Acu-Loc 2 Volar Distal Radius Plate Standard, Right
70-0358-S	70-0358	Acu-Loc 2 Volar Distal Radius Plate, Narrow, Left
70-0359-S	70-0359	Acu-Loc 2 Volar Distal Radius Plate, Narrow, Right
70-0360-S	70-0360	Acu-Loc 2 Volar Distal Radius Plate, Wide, Left
70-0361-S	70-0361	Acu-Loc 2 Volar Distal Radius Plate, Wide, Right
70-0364-S	70-0364	Acu-Loc 2 Volar Distal Radius Extension Plate, Neutral
70-0365-S	70-0365	Acu-Loc 2 Volar Distal Radius Extension Plate, Long, Left
70-0366-S	70-0366	Acu-Loc 2 Volar Distal Radius Extension Plate, Long, Right
70-0368-S	70-0368	Acu-Loc 2 Volar Distal Radius Plate Standard Long, Left
70-0369-S	70-0369	Acu-Loc 2 Volar Distal Radius Plate Standard Long, Right
70-0370-S	70-0370	Acu-Loc 2 Volar Distal Radius Plate, Narrow Long, Left
70-0371-S	70-0371	Acu-Loc 2 Volar Distal Radius Plate, Narrow Long, Right
70-0372-S	70-0372	Acu-Loc 2 Volar Distal Radius Proximal Plate, Standard Long, Left
70-0373-S	70-0373	Acu-Loc 2 Volar Distal Radius Proximal Plate, Standard Long, Right
70-0382-S	70-0382	Acu-Loc 2 Volar Distal Radius Proximal Plate, Narrow Long, Left
70-0383-S	70-0383	Acu-Loc 2 Volar Distal Radius Proximal Plate, Narrow Long, Right
7005-08001-S	7005-08001	0.8 mm Avulsion Hook Plate
CO-3100-S	CO-3100	3.5 mm x 10.0 mm Cortical Screw
CO-3120-S	CO-3120	3.5 mm x 12.0 mm Cortical Screw
CO-3140-S	CO-3140	3.5 mm x 14.0 mm Cortical Screw

Hand and Wrist Portfolio [continued]

Acu-Loc® 2 Wrist Plating System

Sterile Part Number	Non-Sterile Part Number	Description
CO-3160-S	CO-3160	3.5 mm x 16.0 mm Cortical Screw
CO-3180-S	CO-3180	3.5 mm x 18.0 mm Cortical Screw
COL-3080-S	COL-3080	3.5 mm x 8.0 mm Locking Cortical Screw
COL-3100-S	COL-3100	3.5 mm x 10.0 mm Locking Cortical Screw
COL-3120-S	COL-3120	3.5 mm x 12.0 mm Locking Cortical Screw
COL-3140-S	COL-3140	3.5 mm x 14.0 mm Locking Cortical Screw
COL-3160-S	COL-3160	3.5 mm x 16.0 mm Locking Cortical Screw
COL-3180-S	COL-3180	3.5 mm x 18.0 mm Locking Cortical Screw
CO-N2308-S	CO-N2308	2.3 mm x 8 mm Non-Toggling Cortical Screw
CO-N2310-S	CO-N2310	2.3 mm x 10 mm Non-Toggling Cortical Screw
CO-N2312-S	CO-N2312	2.3 mm x 12 mm Non-Toggling Cortical Screw
CO-N2314-S	CO-N2314	2.3 mm x 14 mm Non-Toggling Cortical Screw
CO-N2316-S	CO-N2316	2.3 mm x 16 mm Non-Toggling Cortical Screw
CO-N2318-S	CO-N2318	2.3 mm x 18 mm Non-Toggling Cortical Screw
CO-N2320-S	CO-N2320	2.3 mm x 20 mm Non-Toggling Cortical Screw
CO-N2322-S	CO-N2322	2.3 mm x 22 mm Non-Toggling Cortical Screw
CO-N2324-S	CO-N2324	2.3 mm x 24 mm Non-Toggling Cortical Screw
CO-N2326-S	CO-N2326	2.3 mm x 26 mm Non-Toggling Cortical Screw
CO-N2328-S	CO-N2328	2.3 mm x 28 mm Non-Toggling Cortical Screw
CO-N2330-S	CO-N2330	2.3 mm x 30 mm Non-Toggling Cortical Screw
CO-N2332-S	CO-N2332	2.3 mm x 32 mm Non-Toggling Cortical Screw
CO-N2334-S	CO-N2334	2.3 mm x 34 mm Non-Toggling Cortical Screw
CO-N2336-S	CO-N2336	2.3 mm x 36 mm Non-Toggling Cortical Screw
CO-N2338-S	CO-N2338	2.3 mm x 38 mm Non-Toggling Cortical Screw
CO-N2340-S	CO-N2340	2.3 mm x 40 mm Non-Toggling Cortical Screw
CO-N2342-S	CO-N2342	2.3 mm x 42 mm Non-Toggling Cortical Screw
CO-N2344-S	CO-N2344	2.3 mm x 44 mm Non-Toggling Cortical Screw
CO-N2346-S	CO-N2346	2.3 mm x 46 mm Non-Toggling Cortical Screw
CO-S2308-S	CO-S2308	2.3 mm x 8 mm Locking Cortical Peg

Hand and Wrist Portfolio [continued]

Acu-Loc® 2 Wrist Plating System

Sterile Part Number	Non-Sterile Part Number	Description
CO-S2310-S	CO-S2310	2.3 mm x 10 mm Locking Cortical Peg
CO-S2312-S	CO-S2312	2.3 mm x 12 mm Locking Cortical Peg
CO-S2314-S	CO-S2314	2.3 mm x 14 mm Locking Cortical Peg
CO-S2316-S	CO-S2316	2.3 mm x 16 mm Locking Cortical Peg
CO-S2318-S	CO-S2318	2.3 mm x 18 mm Locking Cortical Peg
CO-S2320-S	CO-S2320	2.3 mm x 20 mm Locking Cortical Peg
CO-S2322-S	CO-S2322	2.3 mm x 22 mm Locking Cortical Peg
CO-S2324-S	CO-S2324	2.3 mm x 24 mm Locking Cortical Peg
CO-S2326-S	CO-S2326	2.3 mm x 26 mm Locking Cortical Peg
CO-S2328-S	CO-S2328	2.3 mm x 28 mm Locking Cortical Peg
CO-T2308-S	CO-T2308	2.3 mm x 8 mm Locking Cortical Screw
CO-T2310-S	CO-T2310	2.3 mm x 10 mm Locking Cortical Screw
CO-T2312-S	CO-T2312	2.3 mm x 12 mm Locking Cortical Screw
CO-T2314-S	CO-T2314	2.3 mm x 14 mm Locking Cortical Screw
CO-T2316-S	CO-T2316	2.3 mm x 16 mm Locking Cortical Screw
CO-T2318-S	CO-T2318	2.3 mm x 18 mm Locking Cortical Screw
CO-T2320-S	CO-T2320	2.3 mm x 20 mm Locking Cortical Screw
CO-T2322-S	CO-T2322	2.3 mm x 22 mm Locking Cortical Screw
CO-T2324-S	CO-T2324	2.3 mm x 24 mm Locking Cortical Screw
CO-T2326-S	CO-T2326	2.3 mm x 26 mm Locking Cortical Screw
CO-T2328-S	CO-T2328	2.3 mm x 28 mm Locking Cortical Screw
CO-T2330-S	CO-T2330	2.3 mm x 30 mm Locking Cortical Screw
CO-T2332-S	CO-T2332	2.3 mm x 32 mm Locking Cortical Screw
CO-T2334-S	CO-T2334	2.3 mm x 34 mm Locking Cortical Screw
CO-T2336-S	CO-T2336	2.3 mm x 36 mm Locking Cortical Screw
CO-T2338-S	CO-T2338	2.3 mm x 38 mm Locking Cortical Screw
CO-T2340-S	CO-T2340	2.3 mm x 40 mm Locking Cortical Screw
CO-T2342-S	CO-T2342	2.3 mm x 42 mm Locking Cortical Screw
CO-T2344-S	CO-T2344	2.3 mm x 44 mm Locking Cortical Screw
CO-T2346-S	CO-T2346	2.3 mm x 46 mm Locking Cortical Screw

Hand and Wrist Portfolio [continued]

Acu-Loc® Wrist Spanning Plate

Sterile Part Number	Non-Sterile Part Number	Description
30-0232-S	30-0232	3.5 mm x 8 mm Locking Hexalobe Screw
30-0233-S	30-0233	3.5 mm x 10 mm Locking Hexalobe Screw
30-0234-S	30-0234	3.5 mm x 12 mm Locking Hexalobe Screw
30-0235-S	30-0235	3.5 mm x 14 mm Locking Hexalobe Screw
30-0236-S	30-0236	3.5 mm x 16 mm Locking Hexalobe Screw
30-0237-S	30-0237	3.5 mm x 18 mm Locking Hexalobe Screw
30-0255-S	30-0255	3.5 mm x 8 mm Nonlocking Hexalobe Screw
30-0256-S	30-0256	3.5 mm x 10 mm Nonlocking Hexalobe Screw
30-0257-S	30-0257	3.5 mm x 12 mm Nonlocking Hexalobe Screw
30-0258-S	30-0258	3.5 mm x 14 mm Nonlocking Hexalobe Screw
30-0259-S	30-0259	3.5 mm x 16 mm Nonlocking Hexalobe Screw
30-0260-S	30-0260	3.5 mm x 18 mm Nonlocking Hexalobe Screw
30-0324-S	30-0324	2.7 mm x 8 mm Locking Hexalobe Screw
30-0325-S	30-0325	2.7 mm x 10 mm Locking Hexalobe Screw
30-0326-S	30-0326	2.7 mm x 12 mm Locking Hexalobe Screw
30-0327-S	30-0327	2.7 mm x 14 mm Locking Hexalobe Screw
30-0328-S	30-0328	2.7 mm x 16 mm Locking Hexalobe Screw
30-0329-S	30-0329	2.7 mm x 18 mm Locking Hexalobe Screw
30-0330-S	30-0330	2.7 mm x 20 mm Locking Hexalobe Screw
30-0331-S	30-0331	2.7 mm x 22 mm Locking Hexalobe Screw
30-0343-S	30-0343	2.7 mm x 8 mm Nonlocking Hexalobe Screw
30-0344-S	30-0344	2.7 mm x 10 mm Nonlocking Hexalobe Screw
30-0345-S	30-0345	2.7 mm x 12 mm Nonlocking Hexalobe Screw
30-0346-S	30-0346	2.7 mm x 14 mm Nonlocking Hexalobe Screw
30-0347-S	30-0347	2.7 mm x 16 mm Nonlocking Hexalobe Screw
30-0348-S	30-0348	2.7 mm x 18 mm Nonlocking Hexalobe Screw
30-0349-S	30-0349	2.7 mm x 20 mm Nonlocking Hexalobe Screw
30-0350-S	30-0350	2.7 mm x 22 mm Nonlocking Hexalobe Screw
7006-1170N-S	N/A	Acu-Loc Wrist Spanning Plate, Short

Hand and Wrist Portfolio [continued]

Acu-Loc® Wrist Spanning Plate

Sterile Part Number	Non-Sterile Part Number	Description
7006-1190N-S	N/A	Acu-Loc Wrist Spanning Plate, Long
N/A	3060-27008	2.7 mm x 8 mm Hexalobe Locking CoCr Screw
N/A	3060-27010	2.7 mm x 10 mm Hexalobe Locking CoCr Screw
N/A	3060-27012	2.7 mm x 12 mm Hexalobe Locking CoCr Screw
N/A	3060-27014	2.7 mm x 14 mm Hexalobe Locking CoCr Screw
N/A	3060-27016	2.7 mm x 16 mm Hexalobe Locking CoCr Screw
N/A	3060-27018	2.7 mm x 18 mm Hexalobe Locking CoCr Screw
N/A	3060-27020	2.7 mm x 20 mm Hexalobe Locking CoCr Screw
N/A	3060-27022	2.7 mm x 22 mm Hexalobe Locking CoCr Screw
N/A	3060-27024	2.7 mm x 24 mm Hexalobe Locking CoCr Screw
N/A	3060-27026	2.7 mm x 26 mm Hexalobe Locking CoCr Screw
N/A	3060-27028	2.7 mm x 28 mm Hexalobe Locking CoCr Screw
N/A	3060-27030	2.7 mm x 30 mm Hexalobe Locking CoCr Screw
N/A	3060-27032	2.7 mm x 32 mm Hexalobe Locking CoCr Screw
N/A	3061-27008	2.7 mm x 8 mm Hexalobe Non-Locking CoCr Screw
N/A	3061-27010	2.7 mm x 10 mm Hexalobe Non-Locking CoCr Screw
N/A	3061-27012	2.7 mm x 12 mm Hexalobe Non-Locking CoCr Screw
N/A	3061-27014	2.7 mm x 14 mm Hexalobe Non-Locking CoCr Screw
N/A	3061-27016	2.7 mm x 16 mm Hexalobe Non-Locking CoCr Screw
N/A	3061-27018	2.7 mm x 18 mm Hexalobe Non-Locking CoCr Screw
N/A	3061-27020	2.7 mm x 20 mm Hexalobe Non-Locking CoCr Screw
N/A	3061-27022	2.7 mm x 22 mm Hexalobe Non-Locking CoCr Screw
N/A	3061-27024	2.7 mm x 24 mm Hexalobe Non-Locking CoCr Screw
N/A	3061-27026	2.7 mm x 26 mm Hexalobe Non-Locking CoCr Screw
N/A	3061-27028	2.7 mm x 28 mm Hexalobe Non-Locking CoCr Screw
N/A	3061-27030	2.7 mm x 30 mm Hexalobe Non-Locking CoCr Screw
N/A	3061-27032	2.7 mm x 32 mm Hexalobe Non-Locking CoCr Screw
N/A	3062-35008	3.5 mm x 8 mm Hexalobe Locking CoCr Screw
N/A	3062-35010	3.5 mm x 10 mm Hexalobe Locking CoCr Screw

Hand and Wrist Portfolio [continued]

Acu-Loc® Wrist Spanning Plate

Sterile Part Number	Non-Sterile Part Number	Description
N/A	3062-35012	3.5 mm x 12 mm Hexalobe Locking CoCr Screw
N/A	3062-35014	3.5 mm x 14 mm Hexalobe Locking CoCr Screw
N/A	3062-35016	3.5 mm x 16 mm Hexalobe Locking CoCr Screw
N/A	3062-35018	3.5 mm x 18 mm Hexalobe Locking CoCr Screw
N/A	3062-35020	3.5 mm x 20 mm Hexalobe Locking CoCr Screw
N/A	3062-35022	3.5 mm x 22 mm Hexalobe Locking CoCr Screw
N/A	3062-35024	3.5 mm x 24 mm Hexalobe Locking CoCr Screw
N/A	3062-35026	3.5 mm x 26 mm Hexalobe Locking CoCr Screw
N/A	3062-35028	3.5 mm x 28 mm Hexalobe Locking CoCr Screw
N/A	3062-35030	3.5 mm x 30 mm Hexalobe Locking CoCr Screw
N/A	3062-35032	3.5 mm x 32 mm Hexalobe Locking CoCr Screw
N/A	3062-35034	3.5 mm x 34 mm Hexalobe Locking CoCr Screw
N/A	3062-35036	3.5 mm x 36 mm Hexalobe Locking CoCr Screw
N/A	3062-35038	3.5 mm x 38 mm Hexalobe Locking CoCr Screw
N/A	3062-35040	3.5 mm x 40 mm Hexalobe Locking CoCr Screw
N/A	3062-35045	3.5 mm x 45 mm Hexalobe Locking CoCr Screw
N/A	3062-35050	3.5 mm x 50 mm Hexalobe Locking CoCr Screw
N/A	3062-35055	3.5 mm x 55 mm Hexalobe Locking CoCr Screw
N/A	3062-35060	3.5 mm x 60 mm Hexalobe Locking CoCr Screw
N/A	3063-35008	3.5 mm x 8 mm Hexalobe Non-Locking CoCr Screw
N/A	3063-35010	3.5 mm x 10 mm Hexalobe Non-Locking CoCr Screw
N/A	3063-35012	3.5 mm x 12 mm Hexalobe Non-Locking CoCr Screw
N/A	3063-35014	3.5 mm x 14 mm Hexalobe Non-Locking CoCr Screw
N/A	3063-35016	3.5 mm x 16 mm Hexalobe Non-Locking CoCr Screw
N/A	3063-35018	3.5 mm x 18 mm Hexalobe Non-Locking CoCr Screw
N/A	3063-35020	3.5 mm x 20 mm Hexalobe Non-Locking CoCr Screw
N/A	3063-35022	3.5 mm x 22 mm Hexalobe Non-Locking CoCr Screw
N/A	3063-35024	3.5 mm x 24 mm Hexalobe Non-Locking CoCr Screw
N/A	3063-35026	3.5 mm x 26 mm Hexalobe Non-Locking CoCr Screw
N/A	3063-35028	3.5 mm x 28 mm Hexalobe Non-Locking CoCr Screw

Hand and Wrist Portfolio [continued]

Acu-Loc® Wrist Spanning Plate

Sterile Part Number	Non-Sterile Part Number	Description
N/A	3063-35030	3.5 mm x 30 mm Hexalobe Non-Locking CoCr Screw
N/A	3063-35032	3.5 mm x 32 mm Hexalobe Non-Locking CoCr Screw
N/A	3063-35034	3.5 mm x 34 mm Hexalobe Non-Locking CoCr Screw
N/A	3063-35036	3.5 mm x 36 mm Hexalobe Non-Locking CoCr Screw
N/A	3063-35038	3.5 mm x 38 mm Hexalobe Non-Locking CoCr Screw
N/A	3063-35040	3.5 mm x 40 mm Hexalobe Non-Locking CoCr Screw
N/A	3063-35045	3.5 mm x 45 mm Hexalobe Non-Locking CoCr Screw
N/A	3063-35050	3.5 mm x 50 mm Hexalobe Non-Locking CoCr Screw
N/A	3063-35055	3.5 mm x 55 mm Hexalobe Non-Locking CoCr Screw
N/A	3063-35060	3.5 mm x 60 mm Hexalobe Non-Locking CoCr Screw
N/A	3063-35065	3.5 mm x 65 mm Hexalobe Non-Locking CoCr Screw

Forearm Fracture Solutions (plates)

Sterile Part Number	Non-Sterile Part Number	Description
30-0232-S	30-0232	3.5 mm x 8 mm Locking Hexalobe Screw
30-0233-S	30-0233	3.5 mm x 10 mm Locking Hexalobe Screw
30-0234-S	30-0234	3.5 mm x 12 mm Locking Hexalobe Screw
30-0235-S	30-0235	3.5 mm x 14 mm Locking Hexalobe Screw
30-0236-S	30-0236	3.5 mm x 16 mm Locking Hexalobe Screw
30-0237-S	30-0237	3.5 mm x 18 mm Locking Hexalobe Screw
30-0238-S	30-0238	3.5 mm x 20 mm Locking Hexalobe Screw
30-0239-S	30-0239	3.5 mm x 22 mm Locking Hexalobe Screw
30-0256-S	30-0256	3.5 mm x 10 mm Nonlocking Hexalobe Screw
30-0257-S	30-0257	3.5 mm x 12 mm Nonlocking Hexalobe Screw
30-0258-S	30-0258	3.5 mm x 14 mm Nonlocking Hexalobe Screw
30-0259-S	30-0259	3.5 mm x 16 mm Nonlocking Hexalobe Screw
30-0260-S	30-0260	3.5 mm x 18 mm Nonlocking Hexalobe Screw
30-0261-S	30-0261	3.5 mm x 20 mm Nonlocking Hexalobe Screw
30-0262-S	30-0262	3.5 mm x 22 mm Nonlocking Hexalobe Screw

Hand and Wrist Portfolio [continued]

Forearm Fracture Solutions (plates)

Sterile Part Number	Non-Sterile Part Number	Description
30-0302-S	30-0302	3.0 mm x 10 mm Nonlocking Hexalobe Screw
30-0303-S	30-0303	3.0 mm x 12 mm Nonlocking Hexalobe Screw
30-0304-S	30-0304	3.0 mm x 14 mm Nonlocking Hexalobe Screw
30-0305-S	30-0305	3.0 mm x 16 mm Nonlocking Hexalobe Screw
30-0306-S	30-0306	3.0 mm x 18 mm Nonlocking Hexalobe Screw
30-0307-S	30-0307	3.0 mm x 20 mm Nonlocking Hexalobe Screw
30-0308-S	30-0308	3.0 mm x 22 mm Nonlocking Hexalobe Screw
70-0066-S	70-0066	6-Hole Volar Midshaft Radius Plate
70-0067-S	70-0067	8-Hole Volar Midshaft Radius Plate
70-0068-S	70-0068	10-Hole Volar Midshaft Radius Plate
70-0069-S	70-0069	12-Hole Volar Midshaft Radius Plate
70-0070-S	70-0070	6-Hole Midshaft Ulna Plate
70-0071-S	70-0071	8-Hole Midshaft Ulna Plate
70-0072-S	70-0072	10-Hole Midshaft Ulna Plate
70-0073-S	70-0073	12-Hole Midshaft Ulna Plate
70-0074-S	70-0074	6-Hole Dorsolateral Midshaft Radius Plate
70-0075-S	70-0075	8-Hole Dorsolateral Midshaft Radius Plate
70-0076-S	70-0076	10-Hole Dorsolateral Midshaft Radius Plate
70-0077-S	70-0077	12-Hole Dorsolateral Midshaft Radius Plate
70-0463-S	N/A	14-Hole Midshaft Ulna Plate
70-0464-S	N/A	16-Hole Midshaft Ulna Plate
70-0466-S	N/A	14-Hole Midshaft Dorsolateral Radius Plate
70-0467-S	N/A	16-Hole Midshaft Dorsolateral Radius Plate
70-0469-S	N/A	14-Hole Volar Midshaft Radius Plate
70-0470-S	N/A	16-Hole Volar Midshaft Radius Plate
CO-3125-S	N/A	3.5 mm x 12.5 mm Cortical Screw
CO-3150-S	N/A	3.5 mm x 15.0 mm Cortical Screw
CO-3175-S	N/A	3.5 mm x 17.5 mm Cortical Screw
CO-3200-S	CO-3200	3.5 mm x 20.0 mm Cortical Screw

Hand and Wrist Portfolio [continued]

Forearm Fracture Solutions (plates)

Sterile Part Number	Non-Sterile Part Number	Description
CO-3225-S	CO-3225	3.5 mm x 22.5 mm Cortical Screw
CO-3250-S	CO-3250	3.5 mm x 25.0 mm Cortical Screw
CO-3275-S	CO-3275	3.5 mm x 27.5 mm Cortical Screw
CO-3300-S	CO-3300	3.5 mm x 30.0 mm Cortical Screw
CO-3325-S	CO-3325	3.5 mm x 32.5 mm Cortical Screw
N/A	CO-T2326	2.3 mm x 26 mm Locking Cortical Screw
N/A	CO-T2328	2.3 mm x 28 mm Locking Cortical Screw

Hand Fracture System

Sterile Part Number	Non-Sterile Part Number	Description
3004-15005-S	3004-15005	1.5 mm x 5 mm Hexalobe MultiScrew
3004-15006-S	3004-15006	1.5 mm x 6 mm Hexalobe MultiScrew
3004-15007-S	3004-15007	1.5 mm x 7 mm Hexalobe MultiScrew
3004-15008-S	3004-15008	1.5 mm x 8 mm Hexalobe MultiScrew
3004-15009-S	3004-15009	1.5 mm x 9 mm Hexalobe MultiScrew
3004-15010-S	3004-15010	1.5 mm x 10 mm Hexalobe MultiScrew
3004-15011-S	3004-15011	1.5 mm x 11 mm Hexalobe MultiScrew
3004-15012-S	3004-15012	1.5 mm x 12 mm Hexalobe MultiScrew
3004-15013-S	3004-15013	1.5 mm x 13 mm Hexalobe MultiScrew
3004-15014-S	3004-15014	1.5 mm x 14 mm Hexalobe MultiScrew
3004-15016-S	3004-15016	1.5 mm x 16 mm Hexalobe MultiScrew
3004-15018-S	3004-15018	1.5 mm x 18 mm Hexalobe MultiScrew
3004-15020-S	3004-15020	1.5 mm x 20 mm Hexalobe MultiScrew
3004-23005-S	3004-23005	2.3 mm x 5 mm Hexalobe MultiScrew
3004-23006-S	3004-23006	2.3 mm x 6 mm Hexalobe MultiScrew
3004-23007-S	3004-23007	2.3 mm x 7 mm Hexalobe MultiScrew
3004-23008-S	3004-23008	2.3 mm x 8 mm Hexalobe MultiScrew
3004-23009-S	3004-23009	2.3 mm x 9 mm Hexalobe MultiScrew

Hand and Wrist Portfolio [continued]

Hand Fracture System

Sterile Part Number	Non-Sterile Part Number	Description
3004-23010-S	3004-23010	2.3 mm x 10 mm Hexalobe MultiScrew
3004-23011-S	3004-23011	2.3 mm x 11 mm Hexalobe MultiScrew
3004-23012-S	3004-23012	2.3 mm x 12 mm Hexalobe MultiScrew
3004-23013-S	3004-23013	2.3 mm x 13 mm Hexalobe MultiScrew
3004-23014-S	3004-23014	2.3 mm x 14 mm Hexalobe MultiScrew
3004-23016-S	3004-23016	2.3 mm x 16 mm Hexalobe MultiScrew
3004-23018-S	3004-23018	2.3 mm x 18 mm Hexalobe MultiScrew
3004-23020-S	3004-23020	2.3 mm x 20 mm Hexalobe MultiScrew
3012-15005-S	3012-15005	1.5 mm x 5 mm Hexalobe Lag Screw
3012-15006-S	3012-15006	1.5 mm x 6 mm Hexalobe Lag Screw
3012-15007-S	3012-15007	1.5 mm x 7 mm Hexalobe Lag Screw
3012-15008-S	3012-15008	1.5 mm x 8 mm Hexalobe Lag Screw
3012-15009-S	3012-15009	1.5 mm x 9 mm Hexalobe Lag Screw
3012-15010-S	3012-15010	1.5 mm x 10 mm Hexalobe Lag Screw
3012-15011-S	3012-15011	1.5 mm x 11 mm Hexalobe Lag Screw
3012-15012-S	3012-15012	1.5 mm x 12 mm Hexalobe Lag Screw
3012-15013-S	3012-15013	1.5 mm x 13 mm Hexalobe Lag Screw
3012-15014-S	3012-15014	1.5 mm x 14 mm Hexalobe Lag Screw
3012-15016-S	3012-15016	1.5 mm x 16 mm Hexalobe Lag Screw
3012-15018-S	3012-15018	1.5 mm x 18 mm Hexalobe Lag Screw
3012-15020-S	3012-15020	1.5 mm x 20 mm Hexalobe Lag Screw
3012-23005-S	3012-23005	2.3 mm x 5 mm Hexalobe Lag Screw
3012-23006-S	3012-23006	2.3 mm x 6 mm Hexalobe Lag Screw
3012-23007-S	3012-23007	2.3 mm x 7 mm Hexalobe Lag Screw
3012-23008-S	3012-23008	2.3 mm x 8 mm Hexalobe Lag Screw
3012-23009-S	3012-23009	2.3 mm x 9 mm Hexalobe Lag Screw
3012-23010-S	3012-23010	2.3 mm x 10 mm Hexalobe Lag Screw
3012-23011-S	3012-23011	2.3 mm x 11 mm Hexalobe Lag Screw
3012-23012-S	3012-23012	2.3 mm x 12 mm Hexalobe Lag Screw

Hand and Wrist Portfolio [continued]

Hand Fracture System

Sterile Part Number	Non-Sterile Part Number	Description
3012-23013-S	3012-23013	2.3 mm x 13 mm Hexalobe Lag Screw
3012-23014-S	3012-23014	2.3 mm x 14 mm Hexalobe Lag Screw
3012-23016-S	3012-23016	2.3 mm x 16 mm Hexalobe Lag Screw
3012-23018-S	3012-23018	2.3 mm x 18 mm Hexalobe Lag Screw
3012-23020-S	3012-23020	2.3 mm x 20 mm Hexalobe Lag Screw
7005-08001-S	7005-08001	0.8 mm Avulsion Hook Plate
7005-08003-S	7005-08003	0.8 mm T-Plate
7005-08004-S	7005-08004	0.8 mm Offset Plate
7005-08006-S	7005-08006	0.8 mm Compression Plate, 6-Hole
7005-08007-S	7005-08007	0.8 mm Curved Medial/Lateral Plate
7005-08010-S	7005-08010	0.8 mm Straight Plate, 10-Hole
7005-13002-S	7005-13002	1.3 mm T-Plate
7005-13004-S	7005-13004	1.3 mm Rolando Fracture Hook Plate
7005-13005-S	7005-13005	1.3 mm Rotational Correction Plate
7005-13006-S	7005-13006	1.3 mm Compression Plate, 6-Hole
7005-13010-S	7005-13010	1.3 mm Straight Plate, 10-Hole
7005-1303L-S	7005-1303L	1.3 mm Metacarpal Neck Plate, Left
7005-1303R-S	7005-1303R	1.3 mm Metacarpal Neck Plate, Right

Hand and Wrist Portfolio [continued]

Modular Hand System

Sterile Part Number	Non-Sterile Part Number	Description
CO-F2106-S	CO-F2106	2.1 mm x 6.0 mm Cruciform Screw
CO-F2108-S	CO-F2108	2.1 mm x 8.0 mm Cruciform Screw
CO-F2110-S	CO-F2110	2.1 mm x 10.0 mm Cruciform Screw
CO-F2112-S	CO-F2112	2.1 mm x 12.0 mm Cruciform Screw
CO-F2114-S	CO-F2114	2.1 mm x 14.0 mm Cruciform Screw
CO-F2116-S	CO-F2116	2.1 mm x 16.0 mm Cruciform Screw
CO-F2706-S	CO-F2706	2.7 mm x 6.0 mm Cruciform Screw
CO-F2708-S	CO-F2708	2.7 mm x 8.0 mm Cruciform Screw
CO-F2710-S	CO-F2710	2.7 mm x 10.0 mm Cruciform Screw
CO-F2712-S	CO-F2712	2.7 mm x 12.0 mm Cruciform Screw
CO-F2714-S	CO-F2714	2.7 mm x 14.0 mm Cruciform Screw
CO-F2716-S	CO-F2716	2.7 mm x 16.0 mm Cruciform Screw
PL-MCPL-S	PL-MCPL	1st Metacarpophalangeal Fusion Plate, Left
PL-MCPR-S	PL-MCPR	1st Metacarpophalangeal Fusion Plate, Right
PL-WF33-S	PL-WF33	Hub Cap® Wrist Fusion Plate, Mini, 3 corner
PL-WF40-S	PL-WF40	Hub Cap Wrist Fusion Plate
PL-WF44-S	PL-WF44	Hub Cap Wrist Fusion Plate, Mini, 4 corner
PL-WF60-S	PL-WF60	Hub Cap Wrist Fusion Plate, Screw Cover

Hand and Wrist Portfolio [continued]

Osteotomy System

Sterile Part Number	Non-Sterile Part Number	Description
30-0232-S	30-0232	3.5 mm x 8 mm Locking Hexalobe Screw
30-0233-S	30-0233	3.5 mm x 10 mm Locking Hexalobe Screw
30-0234-S	30-0234	3.5 mm x 12 mm Locking Hexalobe Screw
30-0235-S	30-0235	3.5 mm x 14 mm Locking Hexalobe Screw
30-0236-S	30-0236	3.5 mm x 16 mm Locking Hexalobe Screw
30-0237-S	30-0237	3.5 mm x 18 mm Locking Hexalobe Screw
30-0238-S	30-0238	3.5 mm x 20 mm Locking Hexalobe Screw
30-0255-S	30-0255	3.5 mm x 8 mm Nonlocking Hexalobe Screw
30-0256-S	30-0256	3.5 mm x 10 mm Nonlocking Hexalobe Screw
30-0257-S	30-0257	3.5 mm x 12 mm Nonlocking Hexalobe Screw
30-0258-S	30-0258	3.5 mm x 14 mm Nonlocking Hexalobe Screw
30-0259-S	30-0259	3.5 mm x 16 mm Nonlocking Hexalobe Screw
30-0260-S	30-0260	3.5 mm x 18 mm Nonlocking Hexalobe Screw
30-0261-S	30-0261	3.5 mm x 20 mm Nonlocking Hexalobe Screw
CO-3080-S	CO-3080	3.5 mm x 8.0 mm Cortical Screw
CO-3100-S	CO-3100	3.5 mm x 10.0 mm Cortical Screw
CO-3120-S	CO-3120	3.5 mm x 12.0 mm Cortical Screw
CO-3140-S	CO-3140	3.5 mm x 14.0 mm Cortical Screw
CO-3160-S	CO-3160	3.5 mm x 16.0 mm Cortical Screw
CO-3180-S	CO-3180	3.5 mm x 18.0 mm Cortical Screw
CO-3200-S	CO-3200	3.5 mm x 20.0 mm Cortical Screw
COL-3080-S	COL-3080	3.5 mm x 8.0 mm Locking Cortical Screw
COL-3100-S	COL-3100	3.5 mm x 10.0 mm Locking Cortical Screw
COL-3120-S	COL-3120	3.5 mm x 12.0 mm Locking Cortical Screw
COL-3140-S	COL-3140	3.5 mm x 14.0 mm Locking Cortical Screw
COL-3160-S	COL-3160	3.5 mm x 16.0 mm Locking Cortical Screw
COL-3180-S	COL-3180	3.5 mm x 18.0 mm Locking Cortical Screw
COL-3200-S	COL-3200	3.5 mm x 20.0 mm Locking Cortical Screw
PL-UL06-S	PL-UL06	6-Hole Ulna Shortening Plate

Hand and Wrist Portfolio [continued]

Total Wrist Fusion Plating System

Sterile Part Number	Non-Sterile Part Number	Description
30-0232-S	30-0232	3.5 mm x 8 mm Locking Hexalobe Screw
30-0233-S	30-0233	3.5 mm x 10 mm Locking Hexalobe Screw
30-0234-S	30-0234	3.5 mm x 12 mm Locking Hexalobe Screw
30-0235-S	30-0235	3.5 mm x 14 mm Locking Hexalobe Screw
30-0236-S	30-0236	3.5 mm x 16 mm Locking Hexalobe Screw
30-0237-S	30-0237	3.5 mm x 18 mm Locking Hexalobe Screw
30-0238-S	30-0238	3.5 mm x 20 mm Locking Hexalobe Screw
30-0239-S	30-0239	3.5 mm x 22 mm Locking Hexalobe Screw
30-0240-S	30-0240	3.5 mm x 24 mm Locking Hexalobe Screw
30-0241-S	30-0241	3.5 mm x 26 mm Locking Hexalobe Screw
30-0255-S	30-0255	3.5 mm x 8 mm Nonlocking Hexalobe Screw
30-0256-S	30-0256	3.5 mm x 10 mm Nonlocking Hexalobe Screw
30-0257-S	30-0257	3.5 mm x 12 mm Nonlocking Hexalobe Screw
30-0258-S	30-0258	3.5 mm x 14 mm Nonlocking Hexalobe Screw
30-0259-S	30-0259	3.5 mm x 16 mm Nonlocking Hexalobe Screw
30-0260-S	30-0260	3.5 mm x 18 mm Nonlocking Hexalobe Screw
30-0261-S	30-0261	3.5 mm x 20 mm Nonlocking Hexalobe Screw
30-0262-S	30-0262	3.5 mm x 22 mm Nonlocking Hexalobe Screw
30-0263-S	30-0263	3.5 mm x 24 mm Nonlocking Hexalobe Screw
30-0264-S	30-0264	3.5 mm x 26 mm Nonlocking Hexalobe Screw
3004-23005-S	3004-23005	2.3 mm x 5 mm Hexalobe MultiScrew
3004-23006-S	3004-23006	2.3 mm x 6 mm Hexalobe MultiScrew
3004-23007-S	3004-23007	2.3 mm x 7 mm Hexalobe MultiScrew
3004-23008-S	3004-23008	2.3 mm x 8 mm Hexalobe MultiScrew
3004-23009-S	3004-23009	2.3 mm x 9 mm Hexalobe MultiScrew
3004-23010-S	3004-23010	2.3 mm x 10 mm Hexalobe MultiScrew
3004-23011-S	3004-23011	2.3 mm x 11 mm Hexalobe MultiScrew
3004-23012-S	3004-23012	2.3 mm x 12 mm Hexalobe MultiScrew
3004-23013-S	3004-23013	2.3 mm x 13 mm Hexalobe MultiScrew

Hand and Wrist Portfolio [continued]

Total Wrist Fusion Plating System

Sterile Part Number	Non-Sterile Part Number	Description
3004-23014-S	3004-23014	2.3 mm x 14 mm Hexalobe MultiScrew
3004-23016-S	3004-23016	2.3 mm x 16 mm Hexalobe MultiScrew
3004-23018-S	3004-23018	2.3 mm x 18 mm Hexalobe MultiScrew
3004-23020-S	3004-23020	2.3 mm x 20 mm Hexalobe MultiScrew
3012-23008-S	3012-23008	2.3 mm x 8 mm Hexalobe Lag Screw
3012-23009-S	3012-23009	2.3 mm x 9 mm Hexalobe Lag Screw
3012-23010-S	3012-23010	2.3 mm x 10 mm Hexalobe Lag Screw
3012-23011-S	3012-23011	2.3 mm x 11 mm Hexalobe Lag Screw
3012-23012-S	3012-23012	2.3 mm x 12 mm Hexalobe Lag Screw
3012-23014-S	3012-23014	2.3 mm x 14 mm Hexalobe Lag Screw
3012-23016-S	3012-23016	2.3 mm x 16 mm Hexalobe Lag Screw
3012-23018-S	3012-23018	2.3 mm x 18 mm Hexalobe Lag Screw
70-0325-S	70-0325	Total Wrist Fusion Plate, Standard, Left
70-0326-S	70-0326	Total Wrist Fusion Plate, Standard, Right
70-0327-S	70-0327	Total Wrist Fusion Plate, Small, Left
70-0328-S	70-0328	Total Wrist Fusion Plate, Small, Right
70-0362-S	70-0362	Total Wrist Fusion Plate, Neutral
CO-3080-S	CO-3080	3.5 mm x 8.0 mm Cortical Screw
CO-3100-S	CO-3100	3.5 mm x 10.0 mm Cortical Screw
CO-3120-S	CO-3120	3.5 mm x 12.0 mm Cortical Screw
CO-3140-S	CO-3140	3.5 mm x 14.0 mm Cortical Screw
CO-3160-S	CO-3160	3.5 mm x 16.0 mm Cortical Screw
CO-3180-S	CO-3180	3.5 mm x 18.0 mm Cortical Screw
CO-3200-S	CO-3200	3.5 mm x 20.0 mm Cortical Screw
CO-3220-S	CO-3220	3.5 mm x 22.0 mm Cortical Screw
CO-3240-S	CO-3240	3.5 mm x 24.0 mm Cortical Screw
CO-3260-S	CO-3260	3.5 mm x 26.0 mm Cortical Screw
COL-3080-S	COL-3080	3.5 mm x 8.0 mm Locking Cortical Screw
COL-3100-S	COL-3100	3.5 mm x 10.0 mm Locking Cortical Screw

Hand and Wrist Portfolio [continued]

Total Wrist Fusion Plating System

Sterile Part Number	Non-Sterile Part Number	Description
COL-3120-S	COL-3120	3.5 mm x 12.0 mm Locking Cortical Screw
COL-3140-S	COL-3140	3.5 mm x 14.0 mm Locking Cortical Screw
COL-3160-S	COL-3160	3.5 mm x 16.0 mm Locking Cortical Screw
COL-3180-S	COL-3180	3.5 mm x 18.0 mm Locking Cortical Screw
COL-3200-S	COL-3200	3.5 mm x 20.0 mm Locking Cortical Screw
COL-3220-S	COL-3220	3.5 mm x 22.0 mm Locking Cortical Screw
COL-3240-S	COL-3240	3.5 mm x 24.0 mm Locking Cortical Screw
COL-3260-S	COL-3260	3.5 mm x 26.0 mm Locking Cortical Screw

Conditions for safe scanning in the MRI

A patient with the peg, plate, screw, sleeve, and screw cover implants from the Acu-Loc® Wrist Plating, Acu-Loc 2 Wrist Plating, Acu-Loc Wrist Spanning Plate, Anatomic Midshaft Forearm Plating, Hand Fracture, Modular Hand, Osteotomy, and Total Wrist Fusion Plating Systems may be safely scanned under the following conditions. Failure to follow these conditions may result in injury to the patient.

Nominal values(s) of Static Magnetic Field	1.5 T or 3 T
Maximum Spatial Field Gradient	30 T/m (3000 gauss/cm)
RF Excitation	Circularly Polarized (CP)
RF Transmit Coil Type	Whole body transmit coil, Head RF transmit-receive coil
Maximum Whole Body SAR at 1.5 T and 3.0 T	1.0 W/kg*
Limits on Scan Duration	1.0 W/kg whole body average SAR for 15 minutes of continuous RF (a sequence or back to back series/scan without breaks) with 5 minutes of cooling for an hour-long scanning session
MR Image Artifact	The presence of this implant may produce an image artifact of 37 mm.
If more information about a specific parameter is not included, there are no conditions associated with that parameter.	

*Whole body SAR of 1.0 W/kg may limit the availability and diagnostic quality of some MR procedures

Screw and Pin Portfolio

- ▶ Acutrak® Headless Compression Screw System
- ▶ Acutrak 2® Headless Compression Screw System
- ▶ Acutrak 3 Headless Compression Screw System
- ▶ AcuTwist® Acutrak Compression Screw System
- ▶ Biotrak® Headless Resorbable Compression System
- ▶ Cannulated Screw System
- ▶ Extremity Screw System
- ▶ Tension Band Pin System/Tension Band Pin System 2

Acutrak® Headless Compression Screw System

Sterile Part Number	Non-Sterile Part Number	Description
30-0033-S	30-0033	30 mm Acutrak Fusion Device 6 mm Nose
AM-0025-S	AM-0025	25.0 mm Acutrak 4/5 Bone Screw
AM-0030-S	AM-0030	30.0 mm Acutrak 4/5 Bone Screw
AM-0035-S	AM-0035	35.0 mm Acutrak 4/5 Bone Screw
AM-0040-S	AM-0040	40.0 mm Acutrak 4/5 Bone Screw
AM-0045-S	AM-0045	45.0 mm Acutrak 4/5 Bone Screw
AM-0050-S	AM-0050	50.0 mm Acutrak 4/5 Bone Screw
AP-0035-S	AP-0035	35.0 mm Acutrak Plus Fixation Screw
AP-0040-S	AP-0040	40.0 mm Acutrak Plus Fixation Screw
AP-0045-S	AP-0045	45.0 mm Acutrak Plus Fixation Screw
AP-0050-S	AP-0050	50.0 mm Acutrak Plus Fixation Screw
AP-0055-S	AP-0055	55.0 mm Acutrak Plus Fixation Screw
AP-0060-S	AP-0060	60.0 mm Acutrak Plus Fixation Screw
AP-0065-S	AP-0065	65.0 mm Acutrak Plus Fixation Screw
AP-0070-S	AP-0070	70.0 mm Acutrak Plus Fixation Screw
AP-0075-S	AP-0075	75.0 mm Acutrak Plus Fixation Screw
AP-0080-S	AP-0080	80.0 mm Acutrak Plus Fixation Screw
AP-6740-S	N/A	40 mm Acutrak 6/7 Fixation Screw
AP-6745-S	N/A	45 mm Acutrak 6/7 Fixation Screw
AP-6750-S	AP-6750	50 mm Acutrak 6/7 Fixation Screw
AP-6755-S	AP-6755	55 mm Acutrak 6/7 Fixation Screw
AP-6760-S	AP-6760	60 mm Acutrak 6/7 Fixation Screw
AP-6765-S	AP-6765	65 mm Acutrak 6/7 Fixation Screw
AP-6770-S	AP-6770	70 mm Acutrak 6/7 Fixation Screw
AP-6775-S	AP-6775	75 mm Acutrak 6/7 Fixation Screw

Screw and Pin Portfolio [continued]

Acutrak® Headless Compression Screw System

Sterile Part Number	Non-Sterile Part Number	Description
AP-6780-S	AP-6780	80 mm Acutrak 6/7 Fixation Screw
AP-6785-S	AP-6785	85 mm Acutrak 6/7 Fixation Screw
AP-6790-S	AP-6790	90 mm Acutrak 6/7 Fixation Screw
AP-6795-S	AP-6795	95 mm Acutrak 6/7 Fixation Screw
AP-67100-S	AP-67100	100 mm Acutrak 6/7 Fixation Screw
AT-1125-S	AT-1125	12.5 mm Standard Acutrak Fixation Screw
AT-1150-S	AT-1150	15.0 mm Standard Acutrak Fixation Screw
AT-1175-S	AT-1175	17.5 mm Standard Acutrak Fixation Screw
AT-1200-S	AT-1200	20.0 mm Standard Acutrak Fixation Screw
AT-1225-S	AT-1225	22.5 mm Standard Acutrak Fixation Screw
AT-1250-S	AT-1250	25.0 mm Standard Acutrak Fixation Screw
AT-1275-S	AT-1275	27.5 mm Standard Acutrak Fixation Screw
AT-1300-S	AT-1300	30.0 mm Standard Acutrak Fixation Screw
ATF-140-S	ATF-140	14.0 mm Acutrak Fusion Device
ATF-160-S	ATF-160	16.0 mm Acutrak Fusion Device
ATF-180-S	ATF-180	18.0 mm Acutrak Fusion Device
ATF-200-S	ATF-200	20.0 mm Acutrak Fusion Device
ATF-220-S	ATF-220	22.0 mm Acutrak Fusion Device
ATF-240-S	ATF-240	24.0 mm Acutrak Fusion Device
ATF-270-S	ATF-270	27.0 mm Acutrak Fusion Device
ATF-320-S	ATF-320	32.0 mm Acutrak Fusion Device
ATF-370-S	ATF-370	37.0 mm Acutrak Fusion Device
ATM-008-S	ATM-008	8.0 mm Mini Acutrak Fixation Screw
ATM-100-S	ATM-100	10.0 mm Mini Acutrak Fixation Screw
ATM-120-S	ATM-120	12.0 mm Mini Acutrak Fixation Screw
ATM-140-S	ATM-140	14.0 mm Mini Acutrak Fixation Screw
ATM-160-S	ATM-160	16.0 mm Mini Acutrak Fixation Screw
ATM-180-S	ATM-180	18.0 mm Mini Acutrak Fixation Screw
ATM-200-S	ATM-200	20.0 mm Mini Acutrak Fixation Screw

Screw and Pin Portfolio [continued]

Acutrak® Headless Compression Screw System

Sterile Part Number	Non-Sterile Part Number	Description
ATM-220-S	ATM-220	22.0 mm Mini Acutrak Fixation Screw
ATM-240-S	ATM-240	24.0 mm Mini Acutrak Fixation Screw
ATM-260-S	ATM-260	26.0 mm Mini Acutrak Fixation Screw

Acutrak 2® Headless Compression Screw System

Sterile Part Number	Non-Sterile Part Number	Description
30-0021-S	30-0021	25.0 mm Acutrak 2 - 5.5 mm Screw
30-0023-S	30-0023	30.0 mm Acutrak 2 - 5.5 mm Screw
30-0025-S	30-0025	35.0 mm Acutrak 2 - 5.5 mm Screw
30-0027-S	30-0027	40.0 mm Acutrak 2 - 5.5 mm Screw
30-0029-S	30-0029	45.0 mm Acutrak 2 - 5.5 mm Screw
30-0031-S	30-0031	50.0 mm Acutrak 2 - 5.5 mm Screw
30-0084-S	30-0084	55.0 mm Acutrak 2 - 5.5 mm Screw
30-0085-S	30-0085	60.0 mm Acutrak 2 - 5.5 mm Screw
30-0620-S	30-0620	20.0 mm Acutrak 2 - 4.7 mm Screw
30-0622-S	30-0622	22.0 mm Acutrak 2 - 4.7 mm Screw
30-0624-S	30-0624	24.0 mm Acutrak 2 - 4.7 mm Screw
30-0626-S	30-0626	26.0 mm Acutrak 2 - 4.7 mm Screw
30-0628-S	30-0628	28.0 mm Acutrak 2 - 4.7 mm Screw
30-0630-S	30-0630	30.0 mm Acutrak 2 - 4.7 mm Screw
30-0635-S	30-0635	35.0 mm Acutrak 2 - 4.7 mm Screw
30-0640-S	30-0640	40.0 mm Acutrak 2 - 4.7 mm Screw
30-0645-S	30-0645	45.0 mm Acutrak 2 - 4.7 mm Screw
30-0650-S	30-0650	50.0 mm Acutrak 2 - 4.7 mm Screw
30-0740-S	30-0740	40.0 mm Acutrak 2 - 7.5 mm Screw
30-0745-S	30-0745	45.0 mm Acutrak 2 - 7.5 mm Screw
30-0750-S	30-0750	50.0 mm Acutrak 2 - 7.5 mm Screw
30-0755-S	30-0755	55.0 mm Acutrak 2 - 7.5 mm Screw
30-0760-S	30-0760	60.0 mm Acutrak 2 - 7.5 mm Screw

Screw and Pin Portfolio [continued]

Acutrak 2® Headless Compression Screw System

Sterile Part Number	Non-Sterile Part Number	Description
30-0765-S	30-0765	65.0 mm Acutrak 2 - 7.5 mm Screw
30-0770-S	30-0770	70.0 mm Acutrak 2 - 7.5 mm Screw
30-0775-S	30-0775	75.0 mm Acutrak 2 - 7.5 mm Screw
30-0780-S	30-0780	80.0 mm Acutrak 2 - 7.5 mm Screw
30-0785-S	30-0785	85.0 mm Acutrak 2 - 7.5 mm Screw
30-0790-S	30-0790	90.0 mm Acutrak 2 - 7.5 mm Screw
30-0795-S	30-0795	95.0 mm Acutrak 2 - 7.5 mm Screw
30-0800-S	30-0800	100.0 mm Acutrak 2 - 7.5 mm Screw
30-0805-S	30-0805	105.0 mm Acutrak 2 - 7.5 mm Screw
30-0810-S	30-0810	110.0 mm Acutrak 2 - 7.5 mm Screw
30-0815-S	30-0815	115.0 mm Acutrak 2 - 7.5 mm Screw
30-0820-S	30-0820	120.0 mm Acutrak 2 - 7.5 mm Screw
AT2-C08-S	AT2-C08	8.0 mm, Micro Acutrak 2 Bone Screw
AT2-C09-S	AT2-C09	9.0 mm, Micro Acutrak 2 Bone Screw
AT2-C10-S	AT2-C10	10.0 mm, Micro Acutrak 2 Bone Screw
AT2-C11-S	AT2-C11	11.0 mm, Micro Acutrak 2 Bone Screw
AT2-C12-S	AT2-C12	12.0 mm, Micro Acutrak 2 Bone Screw
AT2-C13-S	AT2-C13	13.0 mm, Micro Acutrak 2 Bone Screw
AT2-C14-S	AT2-C14	14.0 mm, Micro Acutrak 2 Bone Screw
AT2-C16-S	AT2-C16	16.0 mm, Micro Acutrak 2 Bone Screw
AT2-C18-S	AT2-C18	18.0 mm, Micro Acutrak 2 Bone Screw
AT2-C20-S	AT2-C20	20.0 mm, Micro Acutrak 2 Bone Screw
AT2-C22-S	AT2-C22	22.0 mm, Micro Acutrak 2 Bone Screw
AT2-C24-S	AT2-C24	24.0 mm, Micro Acutrak 2 Bone Screw
AT2-C26-S	AT2-C26	26.0 mm, Micro Acutrak 2 Bone Screw
AT2-C28-S	AT2-C28	28.0 mm, Micro Acutrak 2 Bone Screw
AT2-C30-S	AT2-C30	30.0 mm, Micro Acutrak 2 Bone Screw
AT2-M16-S	AT2-M16	16.0 mm, Mini Acutrak 2 Bone Screw
AT2-M18-S	AT2-M18	18.0 mm, Mini Acutrak 2 Bone Screw

Screw and Pin Portfolio [continued]

Acutrak 2® Headless Compression Screw System

Sterile Part Number	Non-Sterile Part Number	Description
AT2-M20-S	AT2-M20	20.0 mm, Mini Acutrak 2 Bone Screw
AT2-M22-S	AT2-M22	22.0 mm, Mini Acutrak 2 Bone Screw
AT2-M24-S	AT2-M24	24.0 mm, Mini Acutrak 2 Bone Screw
AT2-M26-S	AT2-M26	26.0 mm, Mini Acutrak 2 Bone Screw
AT2-M28-S	AT2-M28	28.0 mm, Mini Acutrak 2 Bone Screw
AT2-M30-S	AT2-M30	30.0 mm, Mini Acutrak 2 Bone Screw
AT2-S16-S	AT2-S16	16.0 mm, Standard Acutrak 2 Bone Screw
AT2-S18-S	AT2-S18	18.0 mm, Standard Acutrak 2 Bone Screw
AT2-S20-S	AT2-S20	20.0 mm, Standard Acutrak 2 Bone Screw
AT2-S22-S	AT2-S22	22.0 mm, Standard Acutrak 2 Bone Screw
AT2-S24-S	AT2-S24	24.0 mm, Standard Acutrak 2 Bone Screw
AT2-S26-S	AT2-S26	26.0 mm, Standard Acutrak 2 Bone Screw
AT2-S28-S	AT2-S28	28.0 mm, Standard Acutrak 2 Bone Screw
AT2-S30-S	AT2-S30	30.0 mm, Standard Acutrak 2 Bone Screw
AT2-S32-S	AT2-S32	32.0 mm, Standard Acutrak 2 Bone Screw
AT2-S34-S	AT2-S34	34.0 mm, Standard Acutrak 2 Bone Screw

Acutrak® 3 Headless Compression Screw System

Sterile Part Number	Non-Sterile Part Number	Description
3050-20008-S	3050-20008	8.0 mm, Nano Acutrak 3 Bone Screw
3050-20009-S	3050-20009	9.0 mm, Nano Acutrak 3 Bone Screw
3050-20010-S	3050-20010	10.0 mm, Nano Acutrak 3 Bone Screw
3050-20011-S	3050-20011	11.0 mm, Nano Acutrak 3 Bone Screw
3050-20012-S	3050-20012	12.0 mm, Nano Acutrak 3 Bone Screw
3050-20013-S	3050-20013	13.0 mm, Nano Acutrak 3 Bone Screw
3050-20014-S	3050-20014	14.0 mm, Nano Acutrak 3 Bone Screw
3050-20016-S	3050-20016	16.0 mm, Nano Acutrak 3 Bone Screw
3050-20018-S	3050-20018	18.0 mm, Nano Acutrak 3 Bone Screw
3050-20020-S	3050-20020	20.0 mm, Nano Acutrak 3 Bone Screw

Screw and Pin Portfolio [continued]

Acutrak® 3 Headless Compression Screw System

Sterile Part Number	Non-Sterile Part Number	Description
3050-20022-S	3050-20022	22.0 mm, Nano Acutrak 3 Bone Screw
3050-20024-S	3050-20024	24.0 mm, Nano Acutrak 3 Bone Screw
3050-20026-S	3050-20026	26.0 mm, Nano Acutrak 3 Bone Screw
3050-20028-S	3050-20028	28.0 mm, Nano Acutrak 3 Bone Screw
3050-20030-S	3050-20030	30.0 mm, Nano Acutrak 3 Bone Screw
3050-20032-S	3050-20032	32.0 mm, Nano Acutrak 3 Bone Screw
3050-20034-S	3050-20034	34.0 mm, Nano Acutrak 3 Bone Screw
3050-20036-S	3050-20036	36.0 mm, Nano Acutrak 3 Bone Screw
3050-20038-S	3050-20038	38.0 mm, Nano Acutrak 3 Bone Screw
3050-20040-S	3050-20040	40.0 mm, Nano Acutrak 3 Bone Screw
3051-25008-S	3051-25008	8.0 mm, Micro Acutrak 3 Bone Screw
3051-25009-S	3051-25009	9.0 mm, Micro Acutrak 3 Bone Screw
3051-25010-S	3051-25010	10.0 mm, Micro Acutrak 3 Bone Screw
3051-25011-S	3051-25011	11.0 mm, Micro Acutrak 3 Bone Screw
3051-25012-S	3051-25012	12.0 mm, Micro Acutrak 3 Bone Screw
3051-25013-S	3051-25013	13.0 mm, Micro Acutrak 3 Bone Screw
3051-25014-S	3051-25014	14.0 mm, Micro Acutrak 3 Bone Screw
3051-25016-S	3051-25016	16.0 mm, Micro Acutrak 3 Bone Screw
3051-25018-S	3051-25018	18.0 mm, Micro Acutrak 3 Bone Screw
3051-25020-S	3051-25020	20.0 mm, Micro Acutrak 3 Bone Screw
3051-25022-S	3051-25022	22.0 mm, Micro Acutrak 3 Bone Screw
3051-25024-S	3051-25024	24.0 mm, Micro Acutrak 3 Bone Screw
3051-25026-S	3051-25026	26.0 mm, Micro Acutrak 3 Bone Screw
3051-25028-S	3051-25028	28.0 mm, Micro Acutrak 3 Bone Screw
3051-25030-S	3051-25030	30.0 mm, Micro Acutrak 3 Bone Screw
3051-25032-S	3051-25032	32.0 mm, Micro Acutrak 3 Bone Screw
3051-25034-S	3051-25034	34.0 mm, Micro Acutrak 3 Bone Screw
3051-25036-S	3051-25036	36.0 mm, Micro Acutrak 3 Bone Screw
3051-25038-S	3051-25038	38.0 mm, Micro Acutrak 3 Bone Screw

Screw and Pin Portfolio [continued]

Acutrak® 3 Headless Compression Screw System

Sterile Part Number	Non-Sterile Part Number	Description
3051-25040-S	3051-25040	40.0 mm, Micro Acutrak 3 Bone Screw
3051-25042-S	3051-25042	42.0 mm, Micro Acutrak 3 Bone Screw
3051-25044-S	3051-25044	44.0 mm, Micro Acutrak 3 Bone Screw
3051-25046-S	3051-25046	46.0 mm, Micro Acutrak 3 Bone Screw
3051-25048-S	3051-25048	48.0 mm, Micro Acutrak 3 Bone Screw
3051-25050-S	3051-25050	50.0 mm, Micro Acutrak 3 Bone Screw
3052-35012-S	3052-35012	12.0 mm, Mini Acutrak 3 Bone Screw
3052-35014-S	3052-35014	14.0 mm, Mini Acutrak 3 Bone Screw
3052-35016-S	3052-35016	16.0 mm, Mini Acutrak 3 Bone Screw
3052-35018-S	3052-35018	18.0 mm, Mini Acutrak 3 Bone Screw
3052-35020-S	3052-35020	20.0 mm, Mini Acutrak 3 Bone Screw
3052-35022-S	3052-35022	22.0 mm, Mini Acutrak 3 Bone Screw
3052-35024-S	3052-35024	24.0 mm, Mini Acutrak 3 Bone Screw
3052-35026-S	3052-35026	26.0 mm, Mini Acutrak 3 Bone Screw
3052-35028-S	3052-35028	28.0 mm, Mini Acutrak 3 Bone Screw
3052-35030-S	3052-35030	30.0 mm, Mini Acutrak 3 Bone Screw
3052-35032-S	3052-35032	32.0 mm, Mini Acutrak 3 Bone Screw
3052-35034-S	3052-35034	34.0 mm, Mini Acutrak 3 Bone Screw
3052-35036-S	3052-35036	36.0 mm, Mini Acutrak 3 Bone Screw
3052-35038-S	3052-35038	38.0 mm, Mini Acutrak 3 Bone Screw
3052-35040-S	3052-35040	40.0 mm, Mini Acutrak 3 Bone Screw
3052-35042-S	3052-35042	42.0 mm, Mini Acutrak 3 Bone Screw
3052-35044-S	3052-35044	44.0 mm, Mini Acutrak 3 Bone Screw
3052-35046-S	3052-35046	46.0 mm, Mini Acutrak 3 Bone Screw
3052-35048-S	3052-35048	48.0 mm, Mini Acutrak 3 Bone Screw
3052-35050-S	3052-35050	50.0 mm, Mini Acutrak 3 Bone Screw
3052-35052-S	3052-35052	52.0 mm, Mini Acutrak 3 Bone Screw
3052-35054-S	3052-35054	54.0 mm, Mini Acutrak 3 Bone Screw
3052-35056-S	3052-35056	56.0 mm, Mini Acutrak 3 Bone Screw
3052-35058-S	3052-35058	58.0 mm, Mini Acutrak 3 Bone Screw

Screw and Pin Portfolio [continued]

Acutrak® 3 Headless Compression Screw System

Sterile Part Number	Non-Sterile Part Number	Description
3052-35060-S	3052-35060	60.0 mm, Mini Acutrak 3 Bone Screw
3053-40016-S	3053-40016	16.0 mm, Standard Acutrak 3 Bone Screw
3053-40018-S	3053-40018	18.0 mm, Standard Acutrak 3 Bone Screw
3053-40020-S	3053-40020	20.0 mm, Standard Acutrak 3 Bone Screw
3053-40022-S	3053-40022	22.0 mm, Standard Acutrak 3 Bone Screw
3053-40024-S	3053-40024	24.0 mm, Standard Acutrak 3 Bone Screw
3053-40026-S	3053-40026	26.0 mm, Standard Acutrak 3 Bone Screw
3053-40028-S	3053-40028	28.0 mm, Standard Acutrak 3 Bone Screw
3053-40030-S	3053-40030	30.0 mm, Standard Acutrak 3 Bone Screw
3053-40032-S	3053-40032	32.0 mm, Standard Acutrak 3 Bone Screw
3053-40034-S	3053-40034	34.0 mm, Standard Acutrak 3 Bone Screw
3053-40036-S	3053-40036	36.0 mm, Standard Acutrak 3 Bone Screw
3053-40038-S	3053-40038	38.0 mm, Standard Acutrak 3 Bone Screw
3053-40040-S	3053-40040	40.0 mm, Standard Acutrak 3 Bone Screw
3053-40042-S	3053-40042	42.0 mm, Standard Acutrak 3 Bone Screw
3053-40044-S	3053-40044	44.0 mm, Standard Acutrak 3 Bone Screw
3053-40046-S	3053-40046	46.0 mm, Standard Acutrak 3 Bone Screw
3053-40048-S	3053-40048	48.0 mm, Standard Acutrak 3 Bone Screw
3053-40050-S	3053-40050	50.0 mm, Standard Acutrak 3 Bone Screw
3053-40052-S	3053-40052	52.0 mm, Standard Acutrak 3 Bone Screw
3053-40054-S	3053-40054	54.0 mm, Standard Acutrak 3 Bone Screw
3053-40056-S	3053-40056	56.0 mm, Standard Acutrak 3 Bone Screw
3053-40058-S	3053-40058	58.0 mm, Standard Acutrak 3 Bone Screw
3053-40060-S	3053-40060	60.0 mm, Standard Acutrak 3 Bone Screw

Screw and Pin Portfolio [continued]

AcuTwist® Acutrak Compression Screw System

Sterile Part Number	Non-Sterile Part Number	Description
AI-0010-S	AI-0010	10 mm AcuTwist Acutrak Compression Screw
AI-0012-S	AI-0012	12 mm AcuTwist Acutrak Compression Screw
AI-0014-S	AI-0014	14 mm AcuTwist Acutrak Compression Screw
AI-0016-S	AI-0016	16 mm AcuTwist Acutrak Compression Screw
AI-0018-S	AI-0018	18 mm AcuTwist Acutrak Compression Screw
AI-0020-S	AI-0020	20 mm AcuTwist Acutrak Compression Screw
AI-0022-S	AI-0022	22 mm AcuTwist Acutrak Compression Screw
AI-0024-S	AI-0024	24 mm AcuTwist Acutrak Compression Screw
AI-0026-S	AI-0026	26 mm AcuTwist Acutrak Compression Screw
AI-0028-S	AI-0028	28 mm AcuTwist Acutrak Compression Screw
AI-0030-S	AI-0030	30 mm AcuTwist Acutrak Compression Screw

Screw and Pin Portfolio [continued]

Biotrak® Headless Resorbable Compression Screw System

Sterile Part Number	Non-Sterile Part Number	Description
30170007-S	N/A	Biotrak 16 mm Screw
30170008-S	N/A	Biotrak 18 mm Screw
30170009-S	N/A	Biotrak 20 mm Screw
30170010-S	N/A	Biotrak 22 mm Screw
30170011-S	N/A	Biotrak 24 mm Screw
30170057-S	N/A	Biotrak Mini 16 mm Screw
30170058-S	N/A	Biotrak Mini 18 mm Screw
30170059-S	N/A	Biotrak Mini 20 mm Screw
30170060-S	N/A	Biotrak Mini 22 mm Screw
30170061-S	N/A	Biotrak Mini 24 mm Screw
30170320-S	N/A	Biotrak 2.0 mm x 20 mm Pin
30170330-S	N/A	Biotrak 2.0 mm x 30 mm Pin
30170340-S	N/A	Biotrak 2.0 mm x 40 mm Pin
55-0001-S	55-0001	Biotrak 20 mm Helical Nail
55-0002-S	55-0002	Biotrak 30 mm Helical Nail
55-0003-S	55-0003	Biotrak 40 mm Helical Nail

Cannulated Screw System

Sterile Part Number	Non-Sterile Part Number	Description
N/A	3005-40010	4.0 mm x 10 mm Cannulated Screw, Short Thread
N/A	3005-40012	4.0 mm x 12 mm Cannulated Screw, Short Thread
N/A	3005-40014	4.0 mm x 14 mm Cannulated Screw, Short Thread
N/A	3005-40016	4.0 mm x 16 mm Cannulated Screw, Short Thread
N/A	3005-40018	4.0 mm x 18 mm Cannulated Screw, Short Thread
N/A	3005-40020	4.0 mm x 20 mm Cannulated Screw, Short Thread
N/A	3005-40022	4.0 mm x 22 mm Cannulated Screw, Short Thread
N/A	3005-40024	4.0 mm x 24 mm Cannulated Screw, Short Thread
N/A	3005-40026	4.0 mm x 26 mm Cannulated Screw, Short Thread
N/A	3005-40028	4.0 mm x 28 mm Cannulated Screw, Short Thread

Screw and Pin Portfolio [continued]

Cannulated Screw System

Sterile Part Number	Non-Sterile Part Number	Description
N/A	3005-40030	4.0 mm x 30 mm Cannulated Screw, Short Thread
N/A	3005-40032	4.0 mm x 32 mm Cannulated Screw, Short Thread
N/A	3005-40034	4.0 mm x 34 mm Cannulated Screw, Short Thread
N/A	3005-40036	4.0 mm x 36 mm Cannulated Screw, Short Thread
N/A	3005-40038	4.0 mm x 38 mm Cannulated Screw, Short Thread
N/A	3005-40040	4.0 mm x 40 mm Cannulated Screw, Short Thread
N/A	3005-40042	4.0 mm x 42 mm Cannulated Screw, Short Thread
N/A	3005-40044	4.0 mm x 44 mm Cannulated Screw, Short Thread
N/A	3005-40046	4.0 mm x 46 mm Cannulated Screw, Short Thread
N/A	3005-40048	4.0 mm x 48 mm Cannulated Screw, Short Thread
N/A	3005-40050	4.0 mm x 50 mm Cannulated Screw, Short Thread
N/A	3005-40052	4.0 mm x 52 mm Cannulated Screw, Short Thread
N/A	3005-40054	4.0 mm x 54 mm Cannulated Screw, Short Thread
N/A	3005-40056	4.0 mm x 56 mm Cannulated Screw, Short Thread
N/A	3005-40058	4.0 mm x 58 mm Cannulated Screw, Short Thread
N/A	3005-40060	4.0 mm x 60 mm Cannulated Screw, Short Thread
N/A	3005-40064	4.0 mm x 64 mm Cannulated Screw, Short Thread
N/A	3005-40068	4.0 mm x 68 mm Cannulated Screw, Short Thread
N/A	3005-40072	4.0 mm x 72 mm Cannulated Screw, Short Thread
N/A	3005-65030	6.5 mm x 30 mm Cannulated Screw, 16 mm Thread
N/A	3005-65035	6.5 mm x 35 mm Cannulated Screw, 16 mm Thread
N/A	3005-65040	6.5 mm x 40 mm Cannulated Screw, 16 mm Thread
N/A	3005-65045	6.5 mm x 45 mm Cannulated Screw, 16 mm Thread
N/A	3005-65050	6.5 mm x 50 mm Cannulated Screw, 16 mm Thread
N/A	3005-65055	6.5 mm x 55 mm Cannulated Screw, 16 mm Thread
N/A	3005-65060	6.5 mm x 60 mm Cannulated Screw, 16 mm Thread
N/A	3005-65065	6.5 mm x 65 mm Cannulated Screw, 16 mm Thread
N/A	3005-65070	6.5 mm x 70 mm Cannulated Screw, 16 mm Thread
N/A	3005-65075	6.5 mm x 75 mm Cannulated Screw, 16 mm Thread

Screw and Pin Portfolio [continued]

Cannulated Screw System

Sterile Part Number	Non-Sterile Part Number	Description
N/A	3005-65080	6.5 mm x 80 mm Cannulated Screw, 16 mm Thread
N/A	3005-65085	6.5 mm x 85 mm Cannulated Screw, 16 mm Thread
N/A	3005-65090	6.5 mm x 90 mm Cannulated Screw, 16 mm Thread
N/A	3005-65095	6.5 mm x 95 mm Cannulated Screw, 16 mm Thread
N/A	3005-65100	6.5 mm x 100 mm Cannulated Screw, 16 mm Thread
N/A	3005-65105	6.5 mm x 105 mm Cannulated Screw, 16 mm Thread
N/A	3005-65110	6.5 mm x 110 mm Cannulated Screw, 16 mm Thread
N/A	3005-65115	6.5 mm x 115 mm Cannulated Screw, 16 mm Thread
N/A	3005-65120	6.5 mm x 120 mm Cannulated Screw, 16 mm Thread
N/A	3005-65125	6.5 mm x 125 mm Cannulated Screw, 16 mm Thread
N/A	3005-65130	6.5 mm x 130 mm Cannulated Screw, 16 mm Thread
N/A	3005-65135	6.5 mm x 135 mm Cannulated Screw, 16 mm Thread
N/A	3005-65140	6.5 mm x 140 mm Cannulated Screw, 16 mm Thread
N/A	3005-65145	6.5 mm x 145 mm Cannulated Screw, 16 mm Thread
N/A	3005-65150	6.5 mm x 150 mm Cannulated Screw, 16 mm Thread
N/A	3005-73030	7.3 mm x 30 mm Cannulated Screw, 16 mm Thread
N/A	3005-73035	7.3 mm x 35 mm Cannulated Screw, 16 mm Thread
N/A	3005-73040	7.3 mm x 40 mm Cannulated Screw, 16 mm Thread
N/A	3005-73045	7.3 mm x 45 mm Cannulated Screw, 16 mm Thread
N/A	3005-73050	7.3 mm x 50 mm Cannulated Screw, 16 mm Thread
N/A	3005-73055	7.3 mm x 55 mm Cannulated Screw, 16 mm Thread
N/A	3005-73060	7.3 mm x 60 mm Cannulated Screw, 16 mm Thread
N/A	3005-73065	7.3 mm x 65 mm Cannulated Screw, 16 mm Thread
N/A	3005-73070	7.3 mm x 70 mm Cannulated Screw, 16 mm Thread
N/A	3005-73075	7.3 mm x 75 mm Cannulated Screw, 16 mm Thread
N/A	3005-73080	7.3 mm x 80 mm Cannulated Screw, 16 mm Thread
N/A	3005-73085	7.3 mm x 85 mm Cannulated Screw, 16 mm Thread
N/A	3005-73090	7.3 mm x 90 mm Cannulated Screw, 16 mm Thread
N/A	3005-73095	7.3 mm x 95 mm Cannulated Screw, 16 mm Thread

Screw and Pin Portfolio [continued]

Cannulated Screw System

Sterile Part Number	Non-Sterile Part Number	Description
N/A	3005-73100	7.3 mm x 100 mm Cannulated Screw, 16 mm Thread
N/A	3005-73105	7.3 mm x 105 mm Cannulated Screw, 16 mm Thread
N/A	3005-73110	7.3 mm x 110 mm Cannulated Screw, 16 mm Thread
N/A	3005-73115	7.3 mm x 115 mm Cannulated Screw, 16 mm Thread
N/A	3005-73120	7.3 mm x 120 mm Cannulated Screw, 16 mm Thread
N/A	3005-73125	7.3 mm x 125 mm Cannulated Screw, 16 mm Thread
N/A	3005-73130	7.3 mm x 130 mm Cannulated Screw, 16 mm Thread
N/A	3005-73135	7.3 mm x 135 mm Cannulated Screw, 16 mm Thread
N/A	3005-73140	7.3 mm x 140 mm Cannulated Screw, 16 mm Thread
N/A	3005-73145	7.3 mm x 145 mm Cannulated Screw, 16 mm Thread
N/A	3005-73150	7.3 mm x 150 mm Cannulated Screw, 16 mm Thread
N/A	3006-40016	4.0 mm x 16 mm Cannulated Screw, Long Thread
N/A	3006-40018	4.0 mm x 18 mm Cannulated Screw, Long Thread
N/A	3006-40020	4.0 mm x 20 mm Cannulated Screw, Long Thread
N/A	3006-40022	4.0 mm x 22 mm Cannulated Screw, Long Thread
N/A	3006-40024	4.0 mm x 24 mm Cannulated Screw, Long Thread
N/A	3006-40026	4.0 mm x 26 mm Cannulated Screw, Long Thread
N/A	3006-40028	4.0 mm x 28 mm Cannulated Screw, Long Thread
N/A	3006-40030	4.0 mm x 30 mm Cannulated Screw, Long Thread
N/A	3006-40032	4.0 mm x 32 mm Cannulated Screw, Long Thread
N/A	3006-40034	4.0 mm x 34 mm Cannulated Screw, Long Thread
3006-40036-S	3006-40036	4.0 mm x 36 mm Cannulated Screw, Long Thread
N/A	3006-40038	4.0 mm x 38 mm Cannulated Screw, Long Thread
N/A	3006-40040	4.0 mm x 40 mm Cannulated Screw, Long Thread
3006-40042-S	3006-40042	4.0 mm x 42 mm Cannulated Screw, Long Thread
N/A	3006-40044	4.0 mm x 44 mm Cannulated Screw, Long Thread
N/A	3006-40046	4.0 mm x 46 mm Cannulated Screw, Long Thread
3006-40048-S	3006-40048	4.0 mm x 48 mm Cannulated Screw, Long Thread
N/A	3006-40050	4.0 mm x 50 mm Cannulated Screw, Long Thread

Screw and Pin Portfolio [continued]

Cannulated Screw System

Sterile Part Number	Non-Sterile Part Number	Description
N/A	3006-40052	4.0 mm x 52 mm Cannulated Screw, Long Thread
N/A	3006-40054	4.0 mm x 54 mm Cannulated Screw, Long Thread
N/A	3006-40056	4.0 mm x 56 mm Cannulated Screw, Long Thread
N/A	3006-40058	4.0 mm x 58 mm Cannulated Screw, Long Thread
N/A	3006-40060	4.0 mm x 60 mm Cannulated Screw, Long Thread
N/A	3006-40064	4.0 mm x 64 mm Cannulated Screw, Long Thread
N/A	3006-40068	4.0 mm x 68 mm Cannulated Screw, Long Thread
N/A	3006-40072	4.0 mm x 72 mm Cannulated Screw, Long Thread
N/A	3006-65045	6.5 mm x 45 mm Cannulated Screw, 32 mm Thread
N/A	3006-65050	6.5 mm x 50 mm Cannulated Screw, 32 mm Thread
N/A	3006-65055	6.5 mm x 55 mm Cannulated Screw, 32 mm Thread
N/A	3006-65060	6.5 mm x 60 mm Cannulated Screw, 32 mm Thread
N/A	3006-65065	6.5 mm x 65 mm Cannulated Screw, 32 mm Thread
N/A	3006-65070	6.5 mm x 70 mm Cannulated Screw, 32 mm Thread
N/A	3006-65075	6.5 mm x 75 mm Cannulated Screw, 32 mm Thread
N/A	3006-65080	6.5 mm x 80 mm Cannulated Screw, 32 mm Thread
N/A	3006-65085	6.5 mm x 85 mm Cannulated Screw, 32 mm Thread
N/A	3006-65090	6.5 mm x 90 mm Cannulated Screw, 32 mm Thread
N/A	3006-65095	6.5 mm x 95 mm Cannulated Screw, 32 mm Thread
N/A	3006-65100	6.5 mm x 100 mm Cannulated Screw, 32 mm Thread
N/A	3006-65105	6.5 mm x 105 mm Cannulated Screw, 32 mm Thread
N/A	3006-65110	6.5 mm x 110 mm Cannulated Screw, 32 mm Thread
N/A	3006-65115	6.5 mm x 115 mm Cannulated Screw, 32 mm Thread
N/A	3006-65120	6.5 mm x 120 mm Cannulated Screw, 32 mm Thread
N/A	3006-65125	6.5 mm x 125 mm Cannulated Screw, 32 mm Thread
N/A	3006-65130	6.5 mm x 130 mm Cannulated Screw, 32 mm Thread
N/A	3006-65135	6.5 mm x 135 mm Cannulated Screw, 32 mm Thread
N/A	3006-65140	6.5 mm x 140 mm Cannulated Screw, 32 mm Thread
N/A	3006-65145	6.5 mm x 145 mm Cannulated Screw, 32 mm Thread

Screw and Pin Portfolio [continued]

Cannulated Screw System

Sterile Part Number	Non-Sterile Part Number	Description
N/A	3006-65150	6.5 mm x 150 mm Cannulated Screw, 32 mm Thread
N/A	3006-73045	7.3 mm x 45 mm Cannulated Screw, 32 mm Thread
N/A	3006-73050	7.3 mm x 50 mm Cannulated Screw, 32 mm Thread
N/A	3006-73055	7.3 mm x 55 mm Cannulated Screw, 32 mm Thread
N/A	3006-73060	7.3 mm x 60 mm Cannulated Screw, 32 mm Thread
N/A	3006-73065	7.3 mm x 65 mm Cannulated Screw, 32 mm Thread
N/A	3006-73070	7.3 mm x 70 mm Cannulated Screw, 32 mm Thread
N/A	3006-73075	7.3 mm x 75 mm Cannulated Screw, 32 mm Thread
N/A	3006-73080	7.3 mm x 80 mm Cannulated Screw, 32 mm Thread
N/A	3006-73085	7.3 mm x 85 mm Cannulated Screw, 32 mm Thread
N/A	3006-73090	7.3 mm x 90 mm Cannulated Screw, 32 mm Thread
N/A	3006-73095	7.3 mm x 95 mm Cannulated Screw, 32 mm Thread
N/A	3006-73100	7.3 mm x 100 mm Cannulated Screw, 32 mm Thread
N/A	3006-73105	7.3 mm x 105 mm Cannulated Screw, 32 mm Thread
N/A	3006-73110	7.3 mm x 110 mm Cannulated Screw, 32 mm Thread
N/A	3006-73115	7.3 mm x 115 mm Cannulated Screw, 32 mm Thread
N/A	3006-73120	7.3 mm x 120 mm Cannulated Screw, 32 mm Thread
N/A	3006-73125	7.3 mm x 125 mm Cannulated Screw, 32 mm Thread
N/A	3006-73130	7.3 mm x 130 mm Cannulated Screw, 32 mm Thread
N/A	3006-73135	7.3 mm x 135 mm Cannulated Screw, 32 mm Thread
N/A	3006-73140	7.3 mm x 140 mm Cannulated Screw, 32 mm Thread
N/A	3006-73145	7.3 mm x 145 mm Cannulated Screw, 32 mm Thread
N/A	3006-73150	7.3 mm x 150 mm Cannulated Screw, 32 mm Thread
N/A	3007-65030	6.5 mm x 30 mm Cannulated Screw, Full Thread
N/A	3007-65035	6.5 mm x 35 mm Cannulated Screw, Full Thread
N/A	3007-65040	6.5 mm x 40 mm Cannulated Screw, Full Thread
N/A	3007-65045	6.5 mm x 45 mm Cannulated Screw, Full Thread
N/A	3007-65050	6.5 mm x 50 mm Cannulated Screw, Full Thread
N/A	3007-65055	6.5 mm x 55 mm Cannulated Screw, Full Thread

Screw and Pin Portfolio [continued]

Cannulated Screw System

Sterile Part Number	Non-Sterile Part Number	Description
N/A	3007-65060	6.5 mm x 60 mm Cannulated Screw, Full Thread
N/A	3007-65065	6.5 mm x 65 mm Cannulated Screw, Full Thread
N/A	3007-65070	6.5 mm x 70 mm Cannulated Screw, Full Thread
N/A	3007-65075	6.5 mm x 75 mm Cannulated Screw, Full Thread
N/A	3007-65080	6.5 mm x 80 mm Cannulated Screw, Full Thread
N/A	3007-65085	6.5 mm x 85 mm Cannulated Screw, Full Thread
N/A	3007-65090	6.5 mm x 90 mm Cannulated Screw, Full Thread
N/A	3007-65095	6.5 mm x 95 mm Cannulated Screw, Full Thread
N/A	3007-65100	6.5 mm x 100 mm Cannulated Screw, Full Thread
N/A	3007-65105	6.5 mm x 105 mm Cannulated Screw, Full Thread
N/A	3007-65110	6.5 mm x 110 mm Cannulated Screw, Full Thread
N/A	3007-65115	6.5 mm x 115 mm Cannulated Screw, Full Thread
N/A	3007-65120	6.5 mm x 120 mm Cannulated Screw, Full Thread
N/A	3007-65125	6.5 mm x 125 mm Cannulated Screw, Full Thread
N/A	3007-65130	6.5 mm x 130 mm Cannulated Screw, Full Thread
N/A	3007-65135	6.5 mm x 135 mm Cannulated Screw, Full Thread
N/A	3007-65140	6.5 mm x 140 mm Cannulated Screw, Full Thread
N/A	3007-65145	6.5 mm x 145 mm Cannulated Screw, Full Thread
N/A	3007-65150	6.5 mm x 150 mm Cannulated Screw, Full Thread
N/A	3007-73030	7.3 mm x 30 mm Cannulated Screw, Full Thread
N/A	3007-73035	7.3 mm x 35 mm Cannulated Screw, Full Thread
N/A	3007-73040	7.3 mm x 40 mm Cannulated Screw, Full Thread
N/A	3007-73045	7.3 mm x 45 mm Cannulated Screw, Full Thread
N/A	3007-73050	7.3 mm x 50 mm Cannulated Screw, Full Thread
N/A	3007-73055	7.3 mm x 55 mm Cannulated Screw, Full Thread
N/A	3007-73060	7.3 mm x 60 mm Cannulated Screw, Full Thread
N/A	3007-73065	7.3 mm x 65 mm Cannulated Screw, Full Thread
N/A	3007-73070	7.3 mm x 70 mm Cannulated Screw, Full Thread
N/A	3007-73075	7.3 mm x 75 mm Cannulated Screw, Full Thread

Screw and Pin Portfolio [continued]

Cannulated Screw System

Sterile Part Number	Non-Sterile Part Number	Description
N/A	3007-73080	7.3 mm x 80 mm Cannulated Screw, Full Thread
N/A	3007-73085	7.3 mm x 85 mm Cannulated Screw, Full Thread
N/A	3007-73090	7.3 mm x 90 mm Cannulated Screw, Full Thread
N/A	3007-73095	7.3 mm x 95 mm Cannulated Screw, Full Thread
N/A	3007-73100	7.3 mm x 100 mm Cannulated Screw, Full Thread
N/A	3007-73105	7.3 mm x 105 mm Cannulated Screw, Full Thread
N/A	3007-73110	7.3 mm x 110 mm Cannulated Screw, Full Thread
N/A	3007-73115	7.3 mm x 115 mm Cannulated Screw, Full Thread
N/A	3007-73120	7.3 mm x 120 mm Cannulated Screw, Full Thread
N/A	3007-73125	7.3 mm x 125 mm Cannulated Screw, Full Thread
N/A	3007-73130	7.3 mm x 130 mm Cannulated Screw, Full Thread
N/A	3007-73135	7.3 mm x 135 mm Cannulated Screw, Full Thread
N/A	3007-73140	7.3 mm x 140 mm Cannulated Screw, Full Thread
N/A	3007-73145	7.3 mm x 145 mm Cannulated Screw, Full Thread
N/A	3007-73150	7.3 mm x 150 mm Cannulated Screw, Full Thread
7003-07036-S	7003-07036	Cannulated Screw Washer 7.0 mm OD x 3.6 mm ID
N/A	7003-13067	Cannulated Screw Washer 13.0 mm OD x 6.7 mm ID

Extremity Screw System

Sterile Part Number	Non-Sterile Part Number	Description
27-0120-S	27-0120	2.7 mm x 12.0 mm Cannulated Extremity Screw
27-0140-S	27-0140	2.7 mm x 14.0 mm Cannulated Extremity Screw
27-0160-S	27-0160	2.7 mm x 16.0 mm Cannulated Extremity Screw
27-0180-S	27-0180	2.7 mm x 18.0 mm Cannulated Extremity Screw
27-0200-S	27-0200	2.7 mm x 20.0 mm Cannulated Extremity Screw
27-0220-S	27-0220	2.7 mm x 22.0 mm Cannulated Extremity Screw
27-0240-S	27-0240	2.7 mm x 24.0 mm Cannulated Extremity Screw
FA-0140-S	FA-0140	3.5 mm x 14.0 mm Cannulated Extremity Screw
FA-0160-S	FA-0160	3.5 mm x 16.0 mm Cannulated Extremity Screw

Screw and Pin Portfolio [continued]

Extremity Screw System

Sterile Part Number	Non-Sterile Part Number	Description
FA-0180-S	FA-0180	3.5 mm x 18.0 mm Cannulated Extremity Screw
FA-0200-S	FA-0200	3.5 mm x 20.0 mm Cannulated Extremity Screw
FA-0220-S	FA-0220	3.5 mm x 22.0 mm Cannulated Extremity Screw
FA-0240-S	FA-0240	3.5 mm x 24.0 mm Cannulated Extremity Screw
FA-0260-S	FA-0260	3.5 mm x 26.0 mm Cannulated Extremity Screw
FA-0280-S	FA-0280	3.5 mm x 28.0 mm Cannulated Extremity Screw
FA-0300-S	FA-0300	3.5 mm x 30.0 mm Cannulated Extremity Screw
FA-0320-S	FA-0320	3.5 mm x 32.0 mm Cannulated Extremity Screw
FA-0340-S	FA-0340	3.5 mm x 34.0 mm Cannulated Extremity Screw
FA-0360-S	FA-0360	3.5 mm x 36.0 mm Cannulated Extremity Screw
FA-0380-S	FA-0380	3.5 mm x 38.0 mm Cannulated Extremity Screw
FA-0400-S	FA-0400	3.5 mm x 40.0 mm Cannulated Extremity Screw
FA-0420-S	FA-0420	3.5 mm x 42.0 mm Cannulated Extremity Screw
FA-4014-S	FA-4014	4.0 mm x 14.0 mm Cannulated Extremity Screw
FA-4016-S	FA-4016	4.0 mm x 16.0 mm Cannulated Extremity Screw
FA-4018-S	FA-4018	4.0 mm x 18.0 mm Cannulated Extremity Screw
FA-4020-S	FA-4020	4.0 mm x 20.0 mm Cannulated Extremity Screw
FA-4022-S	FA-4022	4.0 mm x 22.0 mm Cannulated Extremity Screw
FA-4024-S	FA-4024	4.0 mm x 24.0 mm Cannulated Extremity Screw
FA-4026-S	FA-4026	4.0 mm x 26.0 mm Cannulated Extremity Screw
FA-4028-S	FA-4028	4.0 mm x 28.0 mm Cannulated Extremity Screw
FA-4030-S	FA-4030	4.0 mm x 30.0 mm Cannulated Extremity Screw
FA-4032-S	FA-4032	4.0 mm x 32.0 mm Cannulated Extremity Screw
FA-4034-S	FA-4034	4.0 mm x 34.0 mm Cannulated Extremity Screw
FA-4036-S	FA-4036	4.0 mm x 36.0 mm Cannulated Extremity Screw
FA-4038-S	FA-4038	4.0 mm x 38.0 mm Cannulated Extremity Screw
FA-4040-S	FA-4040	4.0 mm x 40.0 mm Cannulated Extremity Screw
FA-4042-S	FA-4042	4.0 mm x 42.0 mm Cannulated Extremity Screw

Screw and Pin Portfolio [continued]

Tension Band Pin System / Tension Band Pin System 2

Sterile Part Number	Non-Sterile Part Number	Description
N/A	30-0097	50.0 mm Tension Band Pin
N/A	30-0098	70.0 mm Tension Band Pin
N/A	30-0099	90.0 mm Tension Band Pin
N/A	30-0166	2 mm x 50 mm Tension Band Pin
N/A	30-0167	2 mm x 70 mm Tension Band Pin
N/A	30-0168	2 mm x 90 mm Tension Band Pin

Conditions for safe scanning in the MRI

A patient with the screw, fusion device, helical nail, pin, and washer implants from the Acutrak® Headless Compression Screw, Acutrak 2® Headless Compression Screw, Acutrak 3 Headless Compression Screw, AcuTwist® Acutrak Compression Screw, Biotrak® Headless Resorbable Compression, Cannulated Screw, Tension Band Pin/Tension Band Pin 2, and Extremity Screw Systems may be safely scanned under the following conditions. Failure to follow these conditions may result in injury to the patient.

Nominal value(s) of Static Magnetic Field	1.5 T or 3 T
Maximum Spatial Field Gradient	30 T/m (3000 gauss/cm)
RF Excitation	Circularly Polarized (CP)
RF Transmit Coil Type	Whole body transmit coil, Head RF transmit-receive coil
Maximum Whole Body SAR	1.0 W/kg*
Limits on Scan Duration	1.0 W/kg whole body average SAR for 10 minutes of continuous RF (a sequence or back to back series/scan without breaks) followed by a 5 minute cooling period for a one hour session.
MR Image Artifact	The presence of this implant may produce an image artifact of 35 mm.
If more information about a specific parameter is not included, there are no conditions associated with that parameter.	

* Whole body SAR of 1.0 W/kg may limit the availability and diagnostic quality of some MR procedures

Shoulder and Elbow Portfolio

- ▶ Elbow Plating System
- ▶ Congruent Elbow Plating System
- ▶ Radial Head Plating System
- ▶ Clavicle Plating System
- ▶ Clavicle Hook Plating System
- ▶ Locking Clavicle Plating System
- ▶ Acu-Sinch® Repair System
- ▶ Dual-Trak Clavicle Screw System
- ▶ Polarus® 3 Solution (plates)
- ▶ Polarus Proximal Humeral Plating System
- ▶ Scapula Plating System

Elbow Plating System

Sterile Part Number	Non-Sterile Part Number	Description
N/A	30-0098	70.0 mm Tension Band Pin
N/A	30-0099	90.0 mm Tension Band Pin
30-0232-S	30-0232	3.5 mm x 8 mm Locking Hexalobe Screw
30-0233-S	30-0233	3.5 mm x 10 mm Locking Hexalobe Screw
30-0234-S	30-0234	3.5 mm x 12 mm Locking Hexalobe Screw
30-0235-S	30-0235	3.5 mm x 14 mm Locking Hexalobe Screw
30-0236-S	30-0236	3.5 mm x 16 mm Locking Hexalobe Screw
30-0237-S	30-0237	3.5 mm x 18 mm Locking Hexalobe Screw
30-0238-S	30-0238	3.5 mm x 20 mm Locking Hexalobe Screw
30-0239-S	30-0239	3.5 mm x 22 mm Locking Hexalobe Screw
30-0240-S	30-0240	3.5 mm x 24 mm Locking Hexalobe Screw
30-0241-S	30-0241	3.5 mm x 26 mm Locking Hexalobe Screw
30-0242-S	30-0242	3.5 mm x 28 mm Locking Hexalobe Screw
30-0243-S	30-0243	3.5 mm x 30 mm Locking Hexalobe Screw
30-0244-S	30-0244	3.5 mm x 32 mm Locking Hexalobe Screw
30-0245-S	30-0245	3.5 mm x 34 mm Locking Hexalobe Screw
30-0246-S	30-0246	3.5 mm x 36 mm Locking Hexalobe Screw
30-0247-S	30-0247	3.5 mm x 38 mm Locking Hexalobe Screw
30-0248-S	30-0248	3.5 mm x 40 mm Locking Hexalobe Screw
30-0249-S	30-0249	3.5 mm x 45 mm Locking Hexalobe Screw
30-0250-S	30-0250	3.5 mm x 50 mm Locking Hexalobe Screw
30-0251-S	30-0251	3.5 mm x 55 mm Locking Hexalobe Screw
30-0252-S	30-0252	3.5 mm x 60 mm Locking Hexalobe Screw
30-0255-S	30-0255	3.5 mm x 8 mm Nonlocking Hexalobe Screw
30-0256-S	30-0256	3.5 mm x 10 mm Nonlocking Hexalobe Screw

Shoulder and Elbow Portfolio [continued]

Elbow Plating System

Sterile Part Number	Non-Sterile Part Number	Description
30-0257-S	30-0257	3.5 mm x 12 mm Nonlocking Hexalobe Screw
30-0258-S	30-0258	3.5 mm x 14 mm Nonlocking Hexalobe Screw
30-0259-S	30-0259	3.5 mm x 16 mm Nonlocking Hexalobe Screw
30-0260-S	30-0260	3.5 mm x 18 mm Nonlocking Hexalobe Screw
30-0261-S	30-0261	3.5 mm x 20 mm Nonlocking Hexalobe Screw
30-0262-S	30-0262	3.5 mm x 22 mm Nonlocking Hexalobe Screw
30-0263-S	30-0263	3.5 mm x 24 mm Nonlocking Hexalobe Screw
30-0264-S	30-0264	3.5 mm x 26 mm Nonlocking Hexalobe Screw
30-0265-S	30-0265	3.5 mm x 28 mm Nonlocking Hexalobe Screw
30-0266-S	30-0266	3.5 mm x 30 mm Nonlocking Hexalobe Screw
30-0267-S	30-0267	3.5 mm x 32 mm Nonlocking Hexalobe Screw
30-0268-S	30-0268	3.5 mm x 34 mm Nonlocking Hexalobe Screw
30-0269-S	30-0269	3.5 mm x 36 mm Nonlocking Hexalobe Screw
30-0270-S	30-0270	3.5 mm x 38 mm Nonlocking Hexalobe Screw
30-0271-S	30-0271	3.5 mm x 40 mm Nonlocking Hexalobe Screw
30-0272-S	30-0272	3.5 mm x 45 mm Nonlocking Hexalobe Screw
30-0273-S	30-0273	3.5 mm x 50 mm Nonlocking Hexalobe Screw
30-0274-S	30-0274	3.5 mm x 55 mm Nonlocking Hexalobe Screw
30-0275-S	30-0275	3.5 mm x 60 mm Nonlocking Hexalobe Screw
30-0276-S	30-0276	3.5 mm x 65 mm Nonlocking Hexalobe Screw
30-0278-S	30-0278	3.0 mm x 8 mm Locking Hexalobe Screw
30-0279-S	30-0279	3.0 mm x 10 mm Locking Hexalobe Screw
30-0280-S	30-0280	3.0 mm x 12 mm Locking Hexalobe Screw
30-0281-S	30-0281	3.0 mm x 14 mm Locking Hexalobe Screw
30-0282-S	30-0282	3.0 mm x 16 mm Locking Hexalobe Screw
30-0283-S	30-0283	3.0 mm x 18 mm Locking Hexalobe Screw
30-0284-S	30-0284	3.0 mm x 20 mm Locking Hexalobe Screw
30-0285-S	30-0285	3.0 mm x 22 mm Locking Hexalobe Screw
30-0286-S	30-0286	3.0 mm x 24 mm Locking Hexalobe Screw

Shoulder and Elbow Portfolio [continued]

Elbow Plating System

Sterile Part Number	Non-Sterile Part Number	Description
30-0287-S	30-0287	3.0 mm x 26 mm Locking Hexalobe Screw
30-0288-S	30-0288	3.0 mm x 28 mm Locking Hexalobe Screw
30-0289-S	30-0289	3.0 mm x 30 mm Locking Hexalobe Screw
30-0290-S	30-0290	3.0 mm x 32 mm Locking Hexalobe Screw
30-0291-S	30-0291	3.0 mm x 34 mm Locking Hexalobe Screw
30-0292-S	30-0292	3.0 mm x 36 mm Locking Hexalobe Screw
30-0293-S	30-0293	3.0 mm x 38 mm Locking Hexalobe Screw
30-0294-S	30-0294	3.0 mm x 40 mm Locking Hexalobe Screw
30-0295-S	30-0295	3.0 mm x 45 mm Locking Hexalobe Screw
30-0296-S	30-0296	3.0 mm x 50 mm Locking Hexalobe Screw
30-0297-S	30-0297	3.0 mm x 55 mm Locking Hexalobe Screw
30-0298-S	30-0298	3.0 mm x 60 mm Locking Hexalobe Screw
30-0301-S	30-0301	3.0 mm x 8 mm Nonlocking Hexalobe Screw
30-0302-S	30-0302	3.0 mm x 10 mm Nonlocking Hexalobe Screw
30-0303-S	30-0303	3.0 mm x 12 mm Nonlocking Hexalobe Screw
30-0304-S	30-0304	3.0 mm x 14 mm Nonlocking Hexalobe Screw
30-0305-S	30-0305	3.0 mm x 16 mm Nonlocking Hexalobe Screw
30-0306-S	30-0306	3.0 mm x 18 mm Nonlocking Hexalobe Screw
30-0307-S	30-0307	3.0 mm x 20 mm Nonlocking Hexalobe Screw
30-0308-S	30-0308	3.0 mm x 22 mm Nonlocking Hexalobe Screw
30-0309-S	30-0309	3.0 mm x 24 mm Nonlocking Hexalobe Screw
30-0310-S	30-0310	3.0 mm x 26 mm Nonlocking Hexalobe Screw
30-0311-S	30-0311	3.0 mm x 28 mm Nonlocking Hexalobe Screw
30-0312-S	30-0312	3.0 mm x 30 mm Nonlocking Hexalobe Screw
30-0313-S	30-0313	3.0 mm x 32 mm Nonlocking Hexalobe Screw
30-0314-S	30-0314	3.0 mm x 34 mm Nonlocking Hexalobe Screw
30-0315-S	30-0315	3.0 mm x 36 mm Nonlocking Hexalobe Screw
30-0316-S	30-0316	3.0 mm x 38 mm Nonlocking Hexalobe Screw
30-0317-S	30-0317	3.0 mm x 40 mm Nonlocking Hexalobe Screw

Shoulder and Elbow Portfolio [continued]

Elbow Plating System

Sterile Part Number	Non-Sterile Part Number	Description
30-0318-S	30-0318	3.0 mm x 45 mm Nonlocking Hexalobe Screw
30-0319-S	30-0319	3.0 mm x 50 mm Nonlocking Hexalobe Screw
30-0320-S	30-0320	3.0 mm x 55 mm Nonlocking Hexalobe Screw
30-0321-S	30-0321	3.0 mm x 60 mm Nonlocking Hexalobe Screw
30-0322-S	30-0322	3.0 mm x 65 mm Nonlocking Hexalobe Screw
30-0324-S	30-0324	2.7 mm x 8 mm Locking Hexalobe Screw
30-0325-S	30-0325	2.7 mm x 10 mm Locking Hexalobe Screw
30-0326-S	30-0326	2.7 mm x 12 mm Locking Hexalobe Screw
30-0327-S	30-0327	2.7 mm x 14 mm Locking Hexalobe Screw
30-0328-S	30-0328	2.7 mm x 16 mm Locking Hexalobe Screw
30-0329-S	30-0329	2.7 mm x 18 mm Locking Hexalobe Screw
30-0330-S	30-0330	2.7 mm x 20 mm Locking Hexalobe Screw
30-0331-S	30-0331	2.7 mm x 22 mm Locking Hexalobe Screw
30-0332-S	30-0332	2.7 mm x 24 mm Locking Hexalobe Screw
30-0333-S	30-0333	2.7 mm x 26 mm Locking Hexalobe Screw
30-0334-S	30-0334	2.7 mm x 28 mm Locking Hexalobe Screw
30-0335-S	30-0335	2.7 mm x 30 mm Locking Hexalobe Screw
30-0336-S	30-0336	2.7 mm x 32 mm Locking Hexalobe Screw
30-0343-S	30-0343	2.7 mm x 8 mm Nonlocking Hexalobe Screw
30-0344-S	30-0344	2.7 mm x 10 mm Nonlocking Hexalobe Screw
30-0345-S	30-0345	2.7 mm x 12 mm Nonlocking Hexalobe Screw
30-0346-S	30-0346	2.7 mm x 14 mm Nonlocking Hexalobe Screw
30-0347-S	30-0347	2.7 mm x 16 mm Nonlocking Hexalobe Screw
30-0348-S	30-0348	2.7 mm x 18 mm Nonlocking Hexalobe Screw
30-0349-S	30-0349	2.7 mm x 20 mm Nonlocking Hexalobe Screw
30-0350-S	30-0350	2.7 mm x 22 mm Nonlocking Hexalobe Screw
30-0351-S	30-0351	2.7 mm x 24 mm Nonlocking Hexalobe Screw
30-0352-S	30-0352	2.7 mm x 26 mm Nonlocking Hexalobe Screw
30-0353-S	30-0353	2.7 mm x 28 mm Nonlocking Hexalobe Screw

Shoulder and Elbow Portfolio [continued]

Elbow Plating System

Sterile Part Number	Non-Sterile Part Number	Description
30-0354-S	30-0354	2.7 mm x 30 mm Nonlocking Hexalobe Screw
30-0355-S	30-0355	2.7 mm x 32 mm Nonlocking Hexalobe Screw
N/A	3060-27008	2.7 mm x 8 mm Hexalobe Locking CoCr Screw
N/A	3060-27010	2.7 mm x 10 mm Hexalobe Locking CoCr Screw
N/A	3060-27012	2.7 mm x 12 mm Hexalobe Locking CoCr Screw
N/A	3060-27014	2.7 mm x 14 mm Hexalobe Locking CoCr Screw
N/A	3060-27016	2.7 mm x 16 mm Hexalobe Locking CoCr Screw
N/A	3060-27018	2.7 mm x 18 mm Hexalobe Locking CoCr Screw
N/A	3060-27020	2.7 mm x 20 mm Hexalobe Locking CoCr Screw
N/A	3060-27022	2.7 mm x 22 mm Hexalobe Locking CoCr Screw
N/A	3060-27024	2.7 mm x 24 mm Hexalobe Locking CoCr Screw
N/A	3060-27026	2.7 mm x 26 mm Hexalobe Locking CoCr Screw
N/A	3060-27028	2.7 mm x 28 mm Hexalobe Locking CoCr Screw
N/A	3060-27030	2.7 mm x 30 mm Hexalobe Locking CoCr Screw
N/A	3060-27032	2.7 mm x 32 mm Hexalobe Locking CoCr Screw
N/A	3061-27008	2.7 mm x 8 mm Hexalobe Non-Locking CoCr Screw
N/A	3061-27010	2.7 mm x 10 mm Hexalobe Non-Locking CoCr Screw
N/A	3061-27012	2.7 mm x 12 mm Hexalobe Non-Locking CoCr Screw
N/A	3061-27014	2.7 mm x 14 mm Hexalobe Non-Locking CoCr Screw
N/A	3061-27016	2.7 mm x 16 mm Hexalobe Non-Locking CoCr Screw
N/A	3061-27018	2.7 mm x 18 mm Hexalobe Non-Locking CoCr Screw
N/A	3061-27020	2.7 mm x 20 mm Hexalobe Non-Locking CoCr Screw
N/A	3061-27022	2.7 mm x 22 mm Hexalobe Non-Locking CoCr Screw
N/A	3061-27024	2.7 mm x 24 mm Hexalobe Non-Locking CoCr Screw
N/A	3061-27026	2.7 mm x 26 mm Hexalobe Non-Locking CoCr Screw
N/A	3061-27028	2.7 mm x 28 mm Hexalobe Non-Locking CoCr Screw
N/A	3061-27030	2.7 mm x 30 mm Hexalobe Non-Locking CoCr Screw
N/A	3061-27032	2.7 mm x 32 mm Hexalobe Non-Locking CoCr Screw
N/A	3062-35008	3.5 mm x 8 mm Hexalobe Locking CoCr Screw

Shoulder and Elbow Portfolio [continued]

Elbow Plating System

Sterile Part Number	Non-Sterile Part Number	Description
N/A	3062-35010	3.5 mm x 10 mm Hexalobe Locking CoCr Screw
N/A	3062-35012	3.5 mm x 12 mm Hexalobe Locking CoCr Screw
N/A	3062-35014	3.5 mm x 14 mm Hexalobe Locking CoCr Screw
N/A	3062-35016	3.5 mm x 16 mm Hexalobe Locking CoCr Screw
N/A	3062-35018	3.5 mm x 18 mm Hexalobe Locking CoCr Screw
N/A	3062-35020	3.5 mm x 20 mm Hexalobe Locking CoCr Screw
N/A	3062-35022	3.5 mm x 22 mm Hexalobe Locking CoCr Screw
N/A	3062-35024	3.5 mm x 24 mm Hexalobe Locking CoCr Screw
N/A	3062-35026	3.5 mm x 26 mm Hexalobe Locking CoCr Screw
N/A	3062-35028	3.5 mm x 28 mm Hexalobe Locking CoCr Screw
N/A	3062-35030	3.5 mm x 30 mm Hexalobe Locking CoCr Screw
N/A	3062-35032	3.5 mm x 32 mm Hexalobe Locking CoCr Screw
N/A	3062-35034	3.5 mm x 34 mm Hexalobe Locking CoCr Screw
N/A	3062-35036	3.5 mm x 36 mm Hexalobe Locking CoCr Screw
N/A	3062-35038	3.5 mm x 38 mm Hexalobe Locking CoCr Screw
N/A	3062-35040	3.5 mm x 40 mm Hexalobe Locking CoCr Screw
N/A	3062-35045	3.5 mm x 45 mm Hexalobe Locking CoCr Screw
N/A	3062-35050	3.5 mm x 50 mm Hexalobe Locking CoCr Screw
N/A	3062-35055	3.5 mm x 55 mm Hexalobe Locking CoCr Screw
N/A	3062-35060	3.5 mm x 60 mm Hexalobe Locking CoCr Screw
N/A	3063-35008	3.5 mm x 8 mm Hexalobe Non-Locking CoCr Screw
N/A	3063-35010	3.5 mm x 10 mm Hexalobe Non-Locking CoCr Screw
N/A	3063-35012	3.5 mm x 12 mm Hexalobe Non-Locking CoCr Screw
N/A	3063-35014	3.5 mm x 14 mm Hexalobe Non-Locking CoCr Screw
N/A	3063-35016	3.5 mm x 16 mm Hexalobe Non-Locking CoCr Screw
N/A	3063-35018	3.5 mm x 18 mm Hexalobe Non-Locking CoCr Screw
N/A	3063-35020	3.5 mm x 20 mm Hexalobe Non-Locking CoCr Screw
N/A	3063-35022	3.5 mm x 22 mm Hexalobe Non-Locking CoCr Screw

Shoulder and Elbow Portfolio [continued]

Elbow Plating System

Sterile Part Number	Non-Sterile Part Number	Description
N/A	3063-35024	3.5 mm x 24 mm Hexalobe Non-Locking CoCr Screw
N/A	3063-35026	3.5 mm x 26 mm Hexalobe Non-Locking CoCr Screw
N/A	3063-35028	3.5 mm x 28 mm Hexalobe Non-Locking CoCr Screw
N/A	3063-35030	3.5 mm x 30 mm Hexalobe Non-Locking CoCr Screw
N/A	3063-35032	3.5 mm x 32 mm Hexalobe Non-Locking CoCr Screw
N/A	3063-35034	3.5 mm x 34 mm Hexalobe Non-Locking CoCr Screw
N/A	3063-35036	3.5 mm x 36 mm Hexalobe Non-Locking CoCr Screw
N/A	3063-35038	3.5 mm x 38 mm Hexalobe Non-Locking CoCr Screw
N/A	3063-35040	3.5 mm x 40 mm Hexalobe Non-Locking CoCr Screw
N/A	3063-35045	3.5 mm x 45 mm Hexalobe Non-Locking CoCr Screw
N/A	3063-35050	3.5 mm x 50 mm Hexalobe Non-Locking CoCr Screw
N/A	3063-35055	3.5 mm x 55 mm Hexalobe Non-Locking CoCr Screw
N/A	3063-35060	3.5 mm x 60 mm Hexalobe Non-Locking CoCr Screw
N/A	3063-35065	3.5 mm x 65 mm Hexalobe Non-Locking CoCr Screw
70-0302-S	70-0302	Olecranon Plate, Standard 3-Hole, Left
70-0303-S	70-0303	Olecranon Plate, Standard 3-Hole, Right
70-0304-S	70-0304	Olecranon Plate, Standard 5-Hole, Left
70-0305-S	70-0305	Olecranon Plate, Standard 5-Hole, Right
70-0306-S	70-0306	Olecranon Plate, Standard 7-Hole, Left
70-0307-S	70-0307	Olecranon Plate, Standard 7-Hole, Right
70-0308-S	70-0308	Olecranon Plate, Standard 11-Hole, Left
70-0309-S	70-0309	Olecranon Plate, Standard 11-Hole, Right
70-0310-S	70-0310	Olecranon Plate, Standard 15-Hole, Left
70-0311-S	70-0311	Olecranon Plate, Standard 15-Hole, Right
70-0312-S	70-0312	Olecranon Plate, Extended 5-Hole, Left
70-0313-S	70-0313	Olecranon Plate, Extended 5-Hole, Right
70-0314-S	70-0314	Olecranon Plate, Extended 9-Hole, Left
70-0315-S	70-0315	Olecranon Plate, Extended 9-Hole, Right
70-0316-S	70-0316	Olecranon Plate, Narrow 5-Hole, Left

Shoulder and Elbow Portfolio [continued]

Elbow Plating System

Sterile Part Number	Non-Sterile Part Number	Description
70-0317-S	70-0317	Olecranon Plate, Narrow 5-Hole, Right
70-0374-S	70-0374	Posterolateral Distal Humerus Plate 5-Hole, Left
70-0375-S	70-0375	Posterolateral Distal Humerus Plate 5-Hole, Right
70-0376-S	70-0376	Posterolateral Distal Humerus Plate 7-Hole, Left
70-0377-S	70-0377	Posterolateral Distal Humerus Plate 7-Hole, Right
70-0378-S	70-0378	Posterolateral Distal Humerus Plate 11-Hole, Left
70-0379-S	70-0379	Posterolateral Distal Humerus Plate 11-Hole, Right
70-0380-S	70-0380	Posterolateral Distal Humerus Plate 15-Hole, Left
70-0381-S	70-0381	Posterolateral Distal Humerus Plate 15-Hole, Right
70-0413-S	70-0413	Coronoid Plate, Standard, Left
70-0414-S	70-0414	Coronoid Plate, Standard, Right
70-0415-S	70-0415	Coronoid Plate, Sm, Left
70-0416-S	70-0416	Coronoid Plate, Sm, Right
PL-LEL10L-S	PL-LEL10L	Left 10-Hole Locking Lateral Plate
PL-LEL10R-S	PL-LEL10R	Right 10-Hole Locking Lateral Plate
PL-LEL14L-S	PL-LEL14L	Left 14-Hole Locking Lateral Plate
PL-LEL14R-S	PL-LEL14R	Right 14-Hole Locking Lateral Plate
PL-LEL20L-S	PL-LEL20L	Left 20-Hole Locking Lateral Plate
PL-LEL20R-S	PL-LEL20R	Right 20-Hole Locking Lateral Plate
PL-LEL6L-S	PL-LEL6L	Left 6-Hole Locking Lateral Plate
PL-LEL6R-S	PL-LEL6R	Right 6-Hole Locking Lateral Plate
PL-LEM12-S	PL-LEM12	12-Hole Locking Medial Plate
PL-LEM16-S	PL-LEM16	16-Hole Locking Medial Plate
PL-LEM7-S	PL-LEM7	7-Hole Locking Medial Plate
PL-LEM8-S	PL-LEM8	8-Hole Locking Medial Plate
PL-LEM9L-S	PL-LEM9L	Long 9-Hole Locking Medial Plate
PL-LEM9S-S	PL-LEM9S	Short 9-Hole Locking Medial Plate

Shoulder and Elbow Portfolio [continued]

Congruent Elbow Plating System

Sterile Part Number	Non-Sterile Part Number	Description
CA-4120-S	CA-4120	4.0 mm x 12.0 mm Cancellous Screw
CA-4140-S	CA-4140	4.0 mm x 14.0 mm Cancellous Screw
CA-4160-S	CA-4160	4.0 mm x 16.0 mm Cancellous Screw
CA-4180-S	CA-4180	4.0 mm x 18.0 mm Cancellous Screw
CA-4200-S	CA-4200	4.0 mm x 20.0 mm Cancellous Screw
CA-4220-S	CA-4220	4.0 mm x 22.0 mm Cancellous Screw
CA-4240-S	CA-4240	4.0 mm x 24.0 mm Cancellous Screw
CA-4260-S	CA-4260	4.0 mm x 26.0 mm Cancellous Screw
CA-4280-S	CA-4280	4.0 mm x 28.0 mm Cancellous Screw
CA-4300-S	CA-4300	4.0 mm x 30.0 mm Cancellous Screw
CA-4350-S	CA-4350	4.0 mm x 35.0 mm Cancellous Screw
CA-4400-S	CA-4400	4.0 mm x 40.0 mm Cancellous Screw
CA-4450-S	CA-4450	4.0 mm x 45.0 mm Cancellous Screw
CA-4500-S	CA-4500	4.0 mm x 50.0 mm Cancellous Screw
CA-4550-S	CA-4550	4.0 mm x 55.0 mm Cancellous Screw
CA-4600-S	CA-4600	4.0 mm x 60.0 mm Cancellous Screw
CO-2712-S	CO-2712	2.7 mm x 12.0 mm Cortical Screw
CO-2714-S	CO-2714	2.7 mm x 14.0 mm Cortical Screw
CO-2716-S	CO-2716	2.7 mm x 16.0 mm Cortical Screw
CO-2718-S	CO-2718	2.7 mm x 18.0 mm Cortical Screw
CO-2720-S	CO-2720	2.7 mm x 20.0 mm Cortical Screw
CO-2722-S	CO-2722	2.7 mm x 22.0 mm Cortical Screw
CO-2724-S	CO-2724	2.7 mm x 24.0 mm Cortical Screw
CO-2726-S	CO-2726	2.7 mm x 26.0 mm Cortical Screw
CO-2728-S	CO-2728	2.7 mm x 28.0 mm Cortical Screw
CO-2730-S	CO-2730	2.7 mm x 30.0 mm Cortical Screw
CO-2732-S	CO-2732	2.7 mm x 32.0 mm Cortical Screw
CO-2734-S	CO-2734	2.7 mm x 34.0 mm Cortical Screw
CO-2736-S	CO-2736	2.7 mm x 36.0 mm Cortical Screw
CO-2738-S	CO-2738	2.7 mm x 38.0 mm Cortical Screw

Shoulder and Elbow Portfolio [continued]

Congruent Elbow Plating System

Sterile Part Number	Non-Sterile Part Number	Description
CO-2740-S	CO-2740	2.7 mm x 40.0 mm Cortical Screw
CO-2745-S	CO-2745	2.7 mm x 45.0 mm Cortical Screw
CO-2750-S	CO-2750	2.7 mm x 50.0 mm Cortical Screw
CO-2755-S	CO-2755	2.7 mm x 55.0 mm Cortical Screw
CO-2760-S	CO-2760	2.7 mm x 60.0 mm Cortical Screw
CO-2765-S	CO-2765	2.7 mm x 65.0 mm Cortical Screw
CO-3120-S	CO-3120	3.5 mm x 12.0 mm Cortical Screw
CO-3140-S	CO-3140	3.5 mm x 14.0 mm Cortical Screw
CO-3160-S	CO-3160	3.5 mm x 16.0 mm Cortical Screw
CO-3180-S	CO-3180	3.5 mm x 18.0 mm Cortical Screw
CO-3200-S	CO-3200	3.5 mm x 20.0 mm Cortical Screw
CO-3220-S	CO-3220	3.5 mm x 22.0 mm Cortical Screw
CO-3240-S	CO-3240	3.5 mm x 24.0 mm Cortical Screw
CO-3260-S	CO-3260	3.5 mm x 26.0 mm Cortical Screw
CO-3280-S	CO-3280	3.5 mm x 28.0 mm Cortical Screw
CO-3300-S	CO-3300	3.5 mm x 30.0 mm Cortical Screw
CO-3320-S	CO-3320	3.5 mm x 32.0 mm Cortical Screw
CO-3340-S	CO-3340	3.5 mm x 34.0 mm Cortical Screw
CO-3360-S	CO-3360	3.5 mm x 36.0 mm Cortical Screw
CO-3380-S	CO-3380	3.5 mm x 38.0 mm Cortical Screw
CO-3400-S	CO-3400	3.5 mm x 40.0 mm Cortical Screw
CO-3450-S	CO-3450	3.5 mm x 45.0 mm Cortical Screw
CO-3500-S	CO-3500	3.5 mm x 50.0 mm Cortical Screw
CO-3550-S	CO-3550	3.5 mm x 55.0 mm Cortical Screw
CO-3600-S	CO-3600	3.5 mm x 60.0 mm Cortical Screw
CO-3650-S	CO-3650	3.5 mm x 65.0 mm Cortical Screw
COL-2120-S	COL-2120	2.7 mm x 12.0 mm Locking Cortical Screw
COL-2140-S	COL-2140	2.7 mm x 14.0 mm Locking Cortical Screw
COL-2160-S	COL-2160	2.7 mm x 16.0 mm Locking Cortical Screw
COL-2180-S	COL-2180	2.7 mm x 18.0 mm Locking Cortical Screw

Shoulder and Elbow Portfolio [continued]

Congruent Elbow Plating System

Sterile Part Number	Non-Sterile Part Number	Description
COL-2200-S	COL-2200	2.7 mm x 20.0 mm Locking Cortical Screw
COL-2220-S	COL-2220	2.7 mm x 22.0 mm Locking Cortical Screw
COL-3120-S	COL-3120	3.5 mm x 12.0 mm Locking Cortical Screw
COL-3140-S	COL-3140	3.5 mm x 14.0 mm Locking Cortical Screw
COL-3160-S	COL-3160	3.5 mm x 16.0 mm Locking Cortical Screw
COL-3180-S	COL-3180	3.5 mm x 18.0 mm Locking Cortical Screw
COL-3200-S	COL-3200	3.5 mm x 20.0 mm Locking Cortical Screw
COL-3220-S	COL-3220	3.5 mm x 22.0 mm Locking Cortical Screw
COL-3240-S	COL-3240	3.5 mm x 24.0 mm Locking Cortical Screw
COL-3260-S	COL-3260	3.5 mm x 26.0 mm Locking Cortical Screw
COL-3280-S	COL-3280	3.5 mm x 28.0 mm Locking Cortical Screw
COL-3300-S	COL-3300	3.5 mm x 30.0 mm Locking Cortical Screw
COL-3320-S	COL-3320	3.5 mm x 32.0 mm Locking Cortical Screw
COL-3340-S	COL-3340	3.5 mm x 34.0 mm Locking Cortical Screw
COL-3360-S	COL-3360	3.5 mm x 36.0 mm Locking Cortical Screw
COL-3380-S	COL-3380	3.5 mm x 38.0 mm Locking Cortical Screw
COL-3400-S	COL-3400	3.5 mm x 40.0 mm Locking Cortical Screw
COL-3450-S	COL-3450	3.5 mm x 45.0 mm Locking Cortical Screw
COL-3500-S	COL-3500	3.5 mm x 50.0 mm Locking Cortical Screw
FA-CO2736-S	FA-CO2736	2.7 mm x 36 mm Tap-Loc® Screw
FA-CO2738-S	FA-CO2738	2.7 mm x 38 mm Tap-Loc Screw
FA-CO2740-S	FA-CO2740	2.7 mm x 40 mm Tap-Loc Screw
FA-CO2745-S	FA-CO2745	2.7 mm x 45 mm Tap-Loc Screw
FA-CO2750-S	FA-CO2750	2.7 mm x 50 mm Tap-Loc Screw
FA-CO2755-S	FA-CO2755	2.7 mm x 55 mm Tap-Loc Screw
FA-CO2760-S	FA-CO2760	2.7 mm x 60 mm Tap-Loc Screw
FA-CO3536-S	FA-CO3536	3.5 mm x 36 mm Tap-Loc Screw
FA-CO3538-S	FA-CO3538	3.5 mm x 38 mm Tap-Loc Screw
FA-CO3540-S	FA-CO3540	3.5 mm x 40 mm Tap-Loc Screw
FA-CO3545-S	FA-CO3545	3.5 mm x 45 mm Tap-Loc Screw

Shoulder and Elbow Portfolio [continued]

Congruent Elbow Plating System

Sterile Part Number	Non-Sterile Part Number	Description
FA-CO3550-S	FA-CO3550	3.5 mm x 50 mm Tap-Loc Screw
FA-CO3555-S	FA-CO3555	3.5 mm x 55 mm Tap-Loc Screw
FA-CO3560-S	FA-CO3560	3.5 mm x 60 mm Tap-Loc Screw
PL-ELCLL-S	PL-ELCLL	Coronoid Plate, Left, Extended
PL-ELCLR-S	PL-ELCLR	Coronoid Plate, Right, Extended
PL-ELCOL-S	PL-ELCOL	Coronoid Plate, Left
PL-ELCOR-S	PL-ELCOR	Coronoid Plate, Right
PL-ELPO-S	PL-ELPO	Elbow, Posterior Plate
PL-LEL10L-S	PL-LEL10L	Left 10-Hole Locking Lateral Plate
PL-LEL10R-S	PL-LEL10R	Right 10-Hole Locking Lateral Plate
PL-LEL14L-S	PL-LEL14L	Left 14-Hole Locking Lateral Plate
PL-LEL14R-S	PL-LEL14R	Right 14-Hole Locking Lateral Plate
PL-LEL20L-S	PL-LEL20L	Left 20-Hole Locking Lateral Plate
PL-LEL20R-S	PL-LEL20R	Right 20-Hole Locking Lateral Plate
PL-LEL6L-S	PL-LEL6L	Left 6-Hole Locking Lateral Plate
PL-LEL6R-S	PL-LEL6R	Right 6-Hole Locking Lateral Plate
PL-LEM12-S	PL-LEM12	12-Hole Locking Medial Plate
PL-LEM16-S	PL-LEM16	16-Hole Locking Medial Plate
PL-LEM7-S	PL-LEM7	7-Hole Locking Medial Plate
PL-LEM8-S	PL-LEM8	8-Hole Locking Medial Plate
PL-LEM9L-S	PL-LEM9L	Long 9-Hole Locking Medial Plate
PL-LEM9S-S	PL-LEM9S	Short 9-Hole Locking Medial Plate
PL-LEO11-S	PL-LEO11	11-Hole Locking Olecranon Plate
PL-LEO13-S	PL-LEO13	13-Hole Locking Olecranon Plate
PL-LEO13E-S	PL-LEO13E	Extended 13-Hole Locking Olecranon Plate
PL-LEO17L-S	PL-LEO17L	Left 17-Hole Locking Olecranon Plate
PL-LEO17R-S	PL-LEO17R	Right 17-Hole Locking Olecranon Plate
PL-LEO9-S	PL-LEO9	9-Hole Locking Olecranon Plate

Shoulder and Elbow Portfolio [continued]

Radial Head Plating System

Sterile Part Number	Non-Sterile Part Number	Description
70-0097-S	70-0097	Locking Radial Head Plate 3-Hole Standard Curvature
70-0098-S	70-0098	Locking Radial Head Plate 5-Hole Standard Curvature
70-0099-S	70-0099	Locking Radial Head Plate 3-Hole Small Curvature
70-0100-S	70-0100	Locking Radial Head Plate 5-Hole Small Curvature
CO-N2312-S	CO-N2312	2.3 mm x 12 mm Non-Toggling Cortical Screw
CO-N2314-S	CO-N2314	2.3 mm x 14 mm Non-Toggling Cortical Screw
CO-N2316-S	CO-N2316	2.3 mm x 16 mm Non-Toggling Cortical Screw
CO-N2318-S	CO-N2318	2.3 mm x 18 mm Non-Toggling Cortical Screw
CO-N2320-S	CO-N2320	2.3 mm x 20 mm Non-Toggling Cortical Screw
CO-N2322-S	CO-N2322	2.3 mm x 22 mm Non-Toggling Cortical Screw
CO-N2324-S	CO-N2324	2.3 mm x 24 mm Non-Toggling Cortical Screw
CO-N2326-S	CO-N2326	2.3 mm x 26 mm Non-Toggling Cortical Screw
CO-N2328-S	CO-N2328	2.3 mm x 28 mm Non-Toggling Cortical Screw
CO-T2312-S	CO-T2312	2.3 mm x 12 mm Locking Cortical Screw
CO-T2314-S	CO-T2314	2.3 mm x 14 mm Locking Cortical Screw
CO-T2316-S	CO-T2316	2.3 mm x 16 mm Locking Cortical Screw
CO-T2318-S	CO-T2318	2.3 mm x 18 mm Locking Cortical Screw
CO-T2320-S	CO-T2320	2.3 mm x 20 mm Locking Cortical Screw
CO-T2322-S	CO-T2322	2.3 mm x 22 mm Locking Cortical Screw
CO-T2324-S	CO-T2324	2.3 mm x 24 mm Locking Cortical Screw
CO-T2326-S	CO-T2326	2.3 mm x 26 mm Locking Cortical Screw
CO-T2328-S	CO-T2328	2.3 mm x 28 mm Locking Cortical Screw

Shoulder and Elbow Portfolio [continued]

Clavicle Plating System

Sterile Part Number	Non-Sterile Part Number	Description
30-0232-S	30-0232	3.5 mm x 8 mm Locking Hexalobe Screw
30-0233-S	30-0233	3.5 mm x 10 mm Locking Hexalobe Screw
30-0234-S	30-0234	3.5 mm x 12 mm Locking Hexalobe Screw
30-0235-S	30-0235	3.5 mm x 14 mm Locking Hexalobe Screw
30-0236-S	30-0236	3.5 mm x 16 mm Locking Hexalobe Screw
30-0237-S	30-0237	3.5 mm x 18 mm Locking Hexalobe Screw
30-0238-S	30-0238	3.5 mm x 20 mm Locking Hexalobe Screw
30-0239-S	30-0239	3.5 mm x 22 mm Locking Hexalobe Screw
30-0240-S	30-0240	3.5 mm x 24 mm Locking Hexalobe Screw
30-0241-S	30-0241	3.5 mm x 26 mm Locking Hexalobe Screw
30-0255-S	30-0255	3.5 mm x 8 mm Nonlocking Hexalobe Screw
30-0256-S	30-0256	3.5 mm x 10 mm Nonlocking Hexalobe Screw
30-0257-S	30-0257	3.5 mm x 12 mm Nonlocking Hexalobe Screw
30-0258-S	30-0258	3.5 mm x 14 mm Nonlocking Hexalobe Screw
30-0259-S	30-0259	3.5 mm x 16 mm Nonlocking Hexalobe Screw
30-0260-S	30-0260	3.5 mm x 18 mm Nonlocking Hexalobe Screw
30-0261-S	30-0261	3.5 mm x 20 mm Nonlocking Hexalobe Screw
30-0262-S	30-0262	3.5 mm x 22 mm Nonlocking Hexalobe Screw
30-0263-S	30-0263	3.5 mm x 24 mm Nonlocking Hexalobe Screw
30-0264-S	30-0264	3.5 mm x 26 mm Nonlocking Hexalobe Screw
30-0278-S	30-0278	3.0 mm x 8 mm Locking Hexalobe Screw
30-0279-S	30-0279	3.0 mm x 10 mm Locking Hexalobe Screw
30-0280-S	30-0280	3.0 mm x 12 mm Locking Hexalobe Screw
30-0281-S	30-0281	3.0 mm x 14 mm Locking Hexalobe Screw
30-0282-S	30-0282	3.0 mm x 16 mm Locking Hexalobe Screw
30-0283-S	30-0283	3.0 mm x 18 mm Locking Hexalobe Screw
30-0284-S	30-0284	3.0 mm x 20 mm Locking Hexalobe Screw
30-0285-S	30-0285	3.0 mm x 22 mm Locking Hexalobe Screw
30-0286-S	30-0286	3.0 mm x 24 mm Locking Hexalobe Screw

Shoulder and Elbow Portfolio [continued]

Clavicle Plating System

Sterile Part Number	Non-Sterile Part Number	Description
30-0287-S	30-0287	3.0 mm x 26 mm Locking Hexalobe Screw
30-0301-S	30-0301	3.0 mm x 8 mm Nonlocking Hexalobe Screw
30-0302-S	30-0302	3.0 mm x 10 mm Nonlocking Hexalobe Screw
30-0303-S	30-0303	3.0 mm x 12 mm Nonlocking Hexalobe Screw
30-0304-S	30-0304	3.0 mm x 14 mm Nonlocking Hexalobe Screw
30-0305-S	30-0305	3.0 mm x 16 mm Nonlocking Hexalobe Screw
30-0306-S	30-0306	3.0 mm x 18 mm Nonlocking Hexalobe Screw
30-0307-S	30-0307	3.0 mm x 20 mm Nonlocking Hexalobe Screw
30-0308-S	30-0308	3.0 mm x 22 mm Nonlocking Hexalobe Screw
30-0309-S	30-0309	3.0 mm x 24 mm Nonlocking Hexalobe Screw
30-0310-S	30-0310	3.0 mm x 26 mm Nonlocking Hexalobe Screw
70-0111-S	70-0111	Distal Clavicle Plate 3.5 mm 12-Hole, Right
70-0112-S	70-0112	Distal Clavicle Plate 3.5 mm 12-Hole, Left
70-0116-S	70-0116	Distal Clavicle Plate 3.5 mm 9-Hole, Right
70-0117-S	70-0117	Distal Clavicle Plate 3.5 mm 9-Hole, Left
70-0118-S	70-0118	8-Hole Lateral Anterior Clavicle Plate
70-0119-S	70-0119	8-Hole Medial Anterior Clavicle Plate
70-0120-S	70-0120	6-Hole Medial Anterior Clavicle Plate
70-0121-S	70-0121	10-Hole Anterior Clavicle Plate
70-0122-S	70-0122	6-Hole Lateral Anterior Clavicle Plate
70-0123-S	70-0123	Distal Clavicle Plate 2.3 mm 16-Hole, Right
70-0124-S	70-0124	Distal Clavicle Plate 2.3 mm 16-Hole, Left
70-0125-S	70-0125	Distal Clavicle Plate 2.3 mm 13-Hole, Right
70-0126-S	70-0126	Distal Clavicle Plate 2.3 mm 13-Hole, Left
70-0286-S	70-0286	Low-Profile Clavicle Plate, 8-Hole Straight Left
70-0287-S	70-0287	Low-Profile Clavicle Plate, 8-Hole Straight Right
70-0288-S	70-0288	Low-Profile Clavicle Plate, 8-Hole Large Left
70-0289-S	70-0289	Low-Profile Clavicle Plate, 8-Hole Large Right
70-0290-S	70-0290	Low-Profile Clavicle Plate, 8-Hole Medium Left

Shoulder and Elbow Portfolio [continued]

Clavicle Plating System

Sterile Part Number	Non-Sterile Part Number	Description
70-0291-S	70-0291	Low-Profile Clavicle Plate, 8-Hole Medium Right
70-0292-S	70-0292	Low-Profile Clavicle Plate, 8-Hole Small Left
70-0293-S	70-0293	Low-Profile Clavicle Plate, 8-Hole Small Right
70-0294-S	70-0294	Low-Profile Clavicle Plate, 10-Hole Left
70-0295-S	70-0295	Low-Profile Clavicle Plate, 10-Hole Right
70-0296-S	70-0296	Narrow Profile Clavicle Plate, 6-Hole Left
70-0297-S	70-0297	Narrow Profile Clavicle Plate, 6-Hole Right
70-0298-S	70-0298	Narrow Profile Clavicle Plate, 8-Hole Straight Left
70-0299-S	70-0299	Narrow Profile Clavicle Plate, 8-Hole Straight Right
70-0300-S	70-0300	Narrow Profile Clavicle Plate, 8-Hole Large Left
70-0301-S	70-0301	Narrow Profile Clavicle Plate, 8-Hole Large Right
70-0319-S	70-0319	Low-Profile Clavicle J-Plate, 8-Hole Left
70-0320-S	70-0320	Low-Profile Clavicle J-Plate, 8-Hole Right
7002-0416L-S	7002-0416L	Distal Clavicle Plate 3.5 mm 16-Hole Left
7002-0416R-S	7002-0416R	Distal Clavicle Plate 3.5 mm 16-Hole Right
CA-4120-S	CA-4120	4.0 mm x 12.0 mm Cancellous Screw
CA-4140-S	CA-4140	4.0 mm x 14.0 mm Cancellous Screw
CA-4160-S	CA-4160	4.0 mm x 16.0 mm Cancellous Screw
CA-4180-S	CA-4180	4.0 mm x 18.0 mm Cancellous Screw
CA-4200-S	CA-4200	4.0 mm x 20.0 mm Cancellous Screw
CA-4220-S	CA-4220	4.0 mm x 22.0 mm Cancellous Screw
CA-4240-S	CA-4240	4.0 mm x 24.0 mm Cancellous Screw
CA-4260-S	CA-4260	4.0 mm x 26.0 mm Cancellous Screw
CA-4280-S	CA-4280	4.0 mm x 28.0 mm Cancellous Screw
CA-4300-S	CA-4300	4.0 mm x 30.0 mm Cancellous Screw
CA-4350-S	CA-4350	4.0 mm x 35.0 mm Cancellous Screw
CA-4400-S	CA-4400	4.0 mm x 40.0 mm Cancellous Screw
CA-4450-S	CA-4450	4.0 mm x 45.0 mm Cancellous Screw
CA-4500-S	CA-4500	4.0 mm x 50.0 mm Cancellous Screw

Shoulder and Elbow Portfolio [continued]

Clavicle Plating System

Sterile Part Number	Non-Sterile Part Number	Description
CA-4550-S	CA-4550	4.0 mm x 55.0 mm Cancellous Screw
CA-4600-S	CA-4600	4.0 mm x 60.0 mm Cancellous Screw
CO-2708-S	CO-2708	2.7 mm x 8.0 mm Cortical Screw
CO-2710-S	CO-2710	2.7 mm x 10.0 mm Cortical Screw
CO-2712-S	CO-2712	2.7 mm x 12.0 mm Cortical Screw
CO-2714-S	CO-2714	2.7 mm x 14.0 mm Cortical Screw
CO-2716-S	CO-2716	2.7 mm x 16.0 mm Cortical Screw
CO-2718-S	CO-2718	2.7 mm x 18.0 mm Cortical Screw
CO-2720-S	CO-2720	2.7 mm x 20.0 mm Cortical Screw
CO-2722-S	CO-2722	2.7 mm x 22.0 mm Cortical Screw
CO-2724-S	CO-2724	2.7 mm x 24.0 mm Cortical Screw
CO-2726-S	CO-2726	2.7 mm x 26.0 mm Cortical Screw
CO-2728-S	CO-2728	2.7 mm x 28.0 mm Cortical Screw
CO-2730-S	CO-2730	2.7 mm x 30.0 mm Cortical Screw
CO-2732-S	CO-2732	2.7 mm x 32.0 mm Cortical Screw
CO-2734-S	CO-2734	2.7 mm x 34.0 mm Cortical Screw
CO-2736-S	CO-2736	2.7 mm x 36.0 mm Cortical Screw
CO-2738-S	CO-2738	2.7 mm x 38.0 mm Cortical Screw
CO-2740-S	CO-2740	2.7 mm x 40.0 mm Cortical Screw
CO-2745-S	CO-2745	2.7 mm x 45.0 mm Cortical Screw
CO-2750-S	CO-2750	2.7 mm x 50.0 mm Cortical Screw
CO-2755-S	CO-2755	2.7 mm x 55.0 mm Cortical Screw
CO-2760-S	CO-2760	2.7 mm x 60.0 mm Cortical Screw
CO-2765-S	CO-2765	2.7 mm x 65.0 mm Cortical Screw
CO-3060-S	CO-3060	3.5 mm x 6.0 mm Cortical Screw
CO-3080-S	CO-3080	3.5 mm x 8.0 mm Cortical Screw
CO-3100-S	CO-3100	3.5 mm x 10.0 mm Cortical Screw
CO-3120-S	CO-3120	3.5 mm x 12.0 mm Cortical Screw
CO-3140-S	CO-3140	3.5 mm x 14.0 mm Cortical Screw

Shoulder and Elbow Portfolio [continued]

Clavicle Plating System

Sterile Part Number	Non-Sterile Part Number	Description
CO-3160-S	CO-3160	3.5 mm x 16.0 mm Cortical Screw
CO-3180-S	CO-3180	3.5 mm x 18.0 mm Cortical Screw
CO-3200-S	CO-3200	3.5 mm x 20.0 mm Cortical Screw
CO-3220-S	CO-3220	3.5 mm x 22.0 mm Cortical Screw
CO-3240-S	CO-3240	3.5 mm x 24.0 mm Cortical Screw
CO-3260-S	CO-3260	3.5 mm x 26.0 mm Cortical Screw
CO-3280-S	CO-3280	3.5 mm x 28.0 mm Cortical Screw
CO-3300-S	CO-3300	3.5 mm x 30.0 mm Cortical Screw
CO-3320-S	CO-3320	3.5 mm x 32.0 mm Cortical Screw
CO-3340-S	CO-3340	3.5 mm x 34.0 mm Cortical Screw
CO-3360-S	CO-3360	3.5 mm x 36.0 mm Cortical Screw
CO-3380-S	CO-3380	3.5 mm x 38.0 mm Cortical Screw
CO-3400-S	CO-3400	3.5 mm x 40.0 mm Cortical Screw
CO-3450-S	CO-3450	3.5 mm x 45.0 mm Cortical Screw
CO-3500-S	CO-3500	3.5 mm x 50.0 mm Cortical Screw
CO-3550-S	CO-3550	3.5 mm x 55.0 mm Cortical Screw
CO-3600-S	CO-3600	3.5 mm x 60.0 mm Cortical Screw
CO-3650-S	CO-3650	3.5 mm x 65.0 mm Cortical Screw
CO-N2308-S	CO-N2308	2.3 mm x 8 mm Non-Toggling Cortical Screw
CO-N2310-S	CO-N2310	2.3 mm x 10 mm Non-Toggling Cortical Screw
CO-N2312-S	CO-N2312	2.3 mm x 12 mm Non-Toggling Cortical Screw
CO-N2314-S	CO-N2314	2.3 mm x 14 mm Non-Toggling Cortical Screw
CO-N2316-S	CO-N2316	2.3 mm x 16 mm Non-Toggling Cortical Screw
CO-N2318-S	CO-N2318	2.3 mm x 18 mm Non-Toggling Cortical Screw
CO-N2320-S	CO-N2320	2.3 mm x 20 mm Non-Toggling Cortical Screw
CO-N2322-S	CO-N2322	2.3 mm x 22 mm Non-Toggling Cortical Screw
CO-N2324-S	CO-N2324	2.3 mm x 24 mm Non-Toggling Cortical Screw

Shoulder and Elbow Portfolio [continued]

Clavicle Plating System

Sterile Part Number	Non-Sterile Part Number	Description
CO-N2326-S	CO-N2326	2.3 mm x 26 mm Non-Toggling Cortical Screw
CO-T2308-S	CO-T2308	2.3 mm x 8 mm Locking Cortical Screw
CO-T2310-S	CO-T2310	2.3 mm x 10 mm Locking Cortical Screw
CO-T2312-S	CO-T2312	2.3 mm x 12 mm Locking Cortical Screw
CO-T2314-S	CO-T2314	2.3 mm x 14 mm Locking Cortical Screw
CO-T2316-S	CO-T2316	2.3 mm x 16 mm Locking Cortical Screw
CO-T2318-S	CO-T2318	2.3 mm x 18 mm Locking Cortical Screw
CO-T2320-S	CO-T2320	2.3 mm x 20 mm Locking Cortical Screw
CO-T2322-S	CO-T2322	2.3 mm x 22 mm Locking Cortical Screw
CO-T2324-S	CO-T2324	2.3 mm x 24 mm Locking Cortical Screw
CO-T2326-S	CO-T2326	2.3 mm x 26 mm Locking Cortical Screw
COL-2080-S	COL-2080	2.7 mm x 8.0 mm Locking Cortical Screw
COL-2100-S	COL-2100	2.7 mm x 10.0 mm Locking Cortical Screw
COL-2120-S	COL-2120	2.7 mm x 12.0 mm Locking Cortical Screw
COL-2140-S	COL-2140	2.7 mm x 14.0 mm Locking Cortical Screw
COL-2160-S	COL-2160	2.7 mm x 16.0 mm Locking Cortical Screw
COL-2180-S	COL-2180	2.7 mm x 18.0 mm Locking Cortical Screw
COL-2200-S	COL-2200	2.7 mm x 20.0 mm Locking Cortical Screw
COL-2220-S	COL-2220	2.7 mm x 22.0 mm Locking Cortical Screw
COL-2240-S	COL-2240	2.7 mm x 24.0 mm Locking Cortical Screw
COL-2260-S	COL-2260	2.7 mm x 26.0 mm Locking Cortical Screw
COL-2280-S	COL-2280	2.7 mm x 28.0 mm Locking Cortical Screw
COL-2300-S	COL-2300	2.7 mm x 30.0 mm Locking Cortical Screw
COL-2320-S	COL-2320	2.7 mm x 32.0 mm Locking Cortical Screw
COL-2340-S	COL-2340	2.7 mm x 34.0 mm Locking Cortical Screw
COL-2360-S	COL-2360	2.7 mm x 36.0 mm Locking Cortical Screw
COL-2380-S	COL-2380	2.7 mm x 38.0 mm Locking Cortical Screw
COL-2400-S	COL-2400	2.7 mm x 40.0 mm Locking Cortical Screw
COL-2450-S	COL-2450	2.7 mm x 45.0 mm Locking Cortical Screw

Shoulder and Elbow Portfolio [continued]

Clavicle Plating System

Sterile Part Number	Non-Sterile Part Number	Description
COL-2500-S	COL-2500	2.7 mm x 50.0 mm Locking Cortical Screw
COL-2550-S	COL-2550	2.7 mm x 55.0 mm Locking Cortical Screw
COL-2600-S	COL-2600	2.7 mm x 60.0 mm Locking Cortical Screw
COL-2650-S	COL-2650	2.7 mm x 65.0 mm Locking Cortical Screw
COL-3060-S	COL-3060	3.5 mm x 6.0 mm Locking Cortical Screw
COL-3080-S	COL-3080	3.5 mm x 8.0 mm Locking Cortical Screw
COL-3100-S	COL-3100	3.5 mm x 10.0 mm Locking Cortical Screw
COL-3120-S	COL-3120	3.5 mm x 12.0 mm Locking Cortical Screw
COL-3140-S	COL-3140	3.5 mm x 14.0 mm Locking Cortical Screw
COL-3160-S	COL-3160	3.5 mm x 16.0 mm Locking Cortical Screw
COL-3180-S	COL-3180	3.5 mm x 18.0 mm Locking Cortical Screw
COL-3200-S	COL-3200	3.5 mm x 20.0 mm Locking Cortical Screw
COL-3220-S	COL-3220	3.5 mm x 22.0 mm Locking Cortical Screw
COL-3240-S	COL-3240	3.5 mm x 24.0 mm Locking Cortical Screw
COL-3260-S	COL-3260	3.5 mm x 26.0 mm Locking Cortical Screw
COL-3280-S	COL-3280	3.5 mm x 28.0 mm Locking Cortical Screw
COL-3300-S	COL-3300	3.5 mm x 30.0 mm Locking Cortical Screw
COL-3320-S	COL-3320	3.5 mm x 32.0 mm Locking Cortical Screw
COL-3340-S	COL-3340	3.5 mm x 34.0 mm Locking Cortical Screw
COL-3360-S	COL-3360	3.5 mm x 36.0 mm Locking Cortical Screw
COL-3380-S	COL-3380	3.5 mm x 38.0 mm Locking Cortical Screw
COL-3400-S	COL-3400	3.5 mm x 40.0 mm Locking Cortical Screw
COL-3450-S	COL-3450	3.5 mm x 45.0 mm Locking Cortical Screw
COL-3500-S	COL-3500	3.5 mm x 50.0 mm Locking Cortical Screw
COL-3550-S	COL-3550	3.5 mm x 55.0 mm Locking Cortical Screw
COL-3600-S	COL-3600	3.5 mm x 60.0 mm Locking Cortical Screw
COL-3650-S	COL-3650	3.5 mm x 65.0 mm Locking Cortical Screw

Shoulder and Elbow Portfolio [continued]

Clavicle Hook Plating System

Sterile Part Number	Non-Sterile Part Number	Description
7014-0512L-S	7014-0512L	Clavicle Hook Plate, 5-Hole 12 mm, Left
7014-0512R-S	7014-0512R	Clavicle Hook Plate, 5-Hole 12 mm, Right
7014-0516L-S	7014-0516L	Clavicle Hook Plate, 5-Hole 16 mm, Left
7014-0516R-S	7014-0516R	Clavicle Hook Plate, 5-Hole 16 mm, Right
7014-0520L-S	7014-0520L	Clavicle Hook Plate, 5-Hole 20 mm, Left
7014-0520R-S	7014-0520R	Clavicle Hook Plate, 5-Hole 20 mm, Right
7014-0612L-S	7014-0612L	Clavicle Hook Plate, 6-Hole 12 mm, Left
7014-0612R-S	7014-0612R	Clavicle Hook Plate, 6-Hole 12 mm, Right
7014-0616L-S	7014-0616L	Clavicle Hook Plate, 6-Hole 16 mm, Left
7014-0616R-S	7014-0616R	Clavicle Hook Plate, 6-Hole 16 mm, Right
7014-0620L-S	7014-0620L	Clavicle Hook Plate, 6-Hole 20 mm, Left
7014-0620R-S	7014-0620R	Clavicle Hook Plate, 6-Hole 20 mm, Right
7014-0712L-S	7014-0712L	Clavicle Hook Plate, 7-Hole 12 mm, Left
7014-0712R-S	7014-0712R	Clavicle Hook Plate, 7-Hole 12 mm, Right
7014-0716L-S	7014-0716L	Clavicle Hook Plate, 7-Hole 16 mm, Left
7014-0716R-S	7014-0716R	Clavicle Hook Plate, 7-Hole 16 mm, Right
7014-0720L-S	7014-0720L	Clavicle Hook Plate, 7-Hole 20 mm, Left
7014-0720R-S	7014-0720R	Clavicle Hook Plate, 7-Hole 20 mm, Right
7014-0912L-S	7014-0912L	Clavicle Hook Plate, 9-Hole 12 mm, Left
7014-0912R-S	7014-0912R	Clavicle Hook Plate, 9-Hole 12 mm, Right
7014-0916L-S	7014-0916L	Clavicle Hook Plate, 9-Hole 16 mm, Left
7014-0916R-S	7014-0916R	Clavicle Hook Plate, 9-Hole 16 mm, Right
7014-0920L-S	7014-0920L	Clavicle Hook Plate, 9-Hole 20 mm, Left
7014-0920R-S	7014-0920R	Clavicle Hook Plate, 9-Hole 20 mm, Right
30-0232-S	30-0232	3.5 mm x 8 mm Locking Hexalobe Screw
30-0233-S	30-0233	3.5 mm x 10 mm Locking Hexalobe Screw
30-0234-S	30-0234	3.5 mm x 12 mm Locking Hexalobe Screw
30-0235-S	30-0235	3.5 mm x 14 mm Locking Hexalobe Screw
30-0236-S	30-0236	3.5 mm x 16 mm Locking Hexalobe Screw

Shoulder and Elbow Portfolio [continued]

Clavicle Hook Plating System

Sterile Part Number	Non-Sterile Part Number	Description
30-0237-S	30-0237	3.5 mm x 18 mm Locking Hexalobe Screw
30-0238-S	30-0238	3.5 mm x 20 mm Locking Hexalobe Screw
30-0239-S	30-0239	3.5 mm x 22 mm Locking Hexalobe Screw
30-0240-S	30-0240	3.5 mm x 24 mm Locking Hexalobe Screw
30-0241-S	30-0241	3.5 mm x 26 mm Locking Hexalobe Screw
30-0255-S	30-0255	3.5 mm x 8 mm Non-Locking Hexalobe Screw
30-0256-S	30-0256	3.5 mm x 10 mm Non-Locking Hexalobe Screw
30-0257-S	30-0257	3.5 mm x 12 mm Non-Locking Hexalobe Screw
30-0258-S	30-0258	3.5 mm x 14 mm Non-Locking Hexalobe Screw
30-0259-S	30-0259	3.5 mm x 16 mm Non-Locking Hexalobe Screw
30-0260-S	30-0260	3.5 mm x 18 mm Non-Locking Hexalobe Screw
30-0261-S	30-0261	3.5 mm x 20 mm Non-Locking Hexalobe Screw
30-0262-S	30-0262	3.5 mm x 22 mm Non-Locking Hexalobe Screw
30-0263-S	30-0263	3.5 mm x 24 mm Non-Locking Hexalobe Screw
30-0264-S	30-0264	3.5 mm x 26 mm Non-Locking Hexalobe Screw
CO-3080-S	CO-3080	3.5 mm x 8.0 mm Cortical Screw
CO-3100-S	CO-3100	3.5 mm x 10.0 mm Cortical Screw
CO-3120-S	CO-3120	3.5 mm x 12.0 mm Cortical Screw
CO-3140-S	CO-3140	3.5 mm x 14.0 mm Cortical Screw
CO-3160-S	CO-3160	3.5 mm x 16.0 mm Cortical Screw
CO-3180-S	CO-3180	3.5 mm x 18.0 mm Cortical Screw
CO-3200-S	CO-3200	3.5 mm x 20.0 mm Cortical Screw
CO-3220-S	CO-3220	3.5 mm x 22.0 mm Cortical Screw
CO-3240-S	CO-3240	3.5 mm x 24.0 mm Cortical Screw
CO-3260-S	CO-3260	3.5 mm x 26.0 mm Cortical Screw
COL-3080-S	COL-3080	3.5 mm x 8.0 mm Locking Cortical Screw
COL-3100-S	COL-3100	3.5 mm x 10.0 mm Locking Cortical Screw
COL-3120-S	COL-3120	3.5 mm x 12.0 mm Locking Cortical Screw
COL-3140-S	COL-3140	3.5 mm x 14.0 mm Locking Cortical Screw

Shoulder and Elbow Portfolio [continued]

Clavicle Hook Plating System

Sterile Part Number	Non-Sterile Part Number	Description
COL-3160-S	COL-3160	3.5 mm x 16.0 mm Locking Cortical Screw
COL-3180-S	COL-3180	3.5 mm x 18.0 mm Locking Cortical Screw
COL-3200-S	COL-3200	3.5 mm x 20.0 mm Locking Cortical Screw
COL-3220-S	COL-3220	3.5 mm x 22.0 mm Locking Cortical Screw
COL-3240-S	COL-3240	3.5 mm x 24.0 mm Locking Cortical Screw
COL-3260-S	COL-3260	3.5 mm x 26.0 mm Locking Cortical Screw

Locking Clavicle Plating System

Sterile Part Number	Non-Sterile Part Number	Description
30-0232-S	30-0232	3.5 mm x 8 mm Locking Hexalobe Screw
30-0233-S	30-0233	3.5 mm x 10 mm Locking Hexalobe Screw
30-0234-S	30-0234	3.5 mm x 12 mm Locking Hexalobe Screw
30-0235-S	30-0235	3.5 mm x 14 mm Locking Hexalobe Screw
30-0236-S	30-0236	3.5 mm x 16 mm Locking Hexalobe Screw
30-0237-S	30-0237	3.5 mm x 18 mm Locking Hexalobe Screw
30-0238-S	30-0238	3.5 mm x 20 mm Locking Hexalobe Screw
30-0239-S	30-0239	3.5 mm x 22 mm Locking Hexalobe Screw
30-0240-S	30-0240	3.5 mm x 24 mm Locking Hexalobe Screw
30-0241-S	30-0241	3.5 mm x 26 mm Locking Hexalobe Screw
30-0255-S	30-0255	3.5 mm x 8 mm Nonlocking Hexalobe Screw
30-0232-S	30-0232	3.5 mm x 8 mm Locking Hexalobe Screw
30-0233-S	30-0233	3.5 mm x 10 mm Locking Hexalobe Screw
30-0234-S	30-0234	3.5 mm x 12 mm Locking Hexalobe Screw
30-0235-S	30-0235	3.5 mm x 14 mm Locking Hexalobe Screw
30-0236-S	30-0236	3.5 mm x 16 mm Locking Hexalobe Screw
30-0237-S	30-0237	3.5 mm x 18 mm Locking Hexalobe Screw
30-0238-S	30-0238	3.5 mm x 20 mm Locking Hexalobe Screw
30-0239-S	30-0239	3.5 mm x 22 mm Locking Hexalobe Screw

Shoulder and Elbow Portfolio [continued]

Locking Clavicle Plating System

Sterile Part Number	Non-Sterile Part Number	Description
30-0240-S	30-0240	3.5 mm x 24 mm Locking Hexalobe Screw
30-0241-S	30-0241	3.5 mm x 26 mm Locking Hexalobe Screw
30-0255-S	30-0255	3.5 mm x 8 mm Nonlocking Hexalobe Screw
30-0256-S	30-0256	3.5 mm x 10 mm Nonlocking Hexalobe Screw
30-0257-S	30-0257	3.5 mm x 12 mm Nonlocking Hexalobe Screw
30-0258-S	30-0258	3.5 mm x 14 mm Nonlocking Hexalobe Screw
30-0259-S	30-0259	3.5 mm x 16 mm Nonlocking Hexalobe Screw
30-0260-S	30-0260	3.5 mm x 18 mm Nonlocking Hexalobe Screw
30-0261-S	30-0261	3.5 mm x 20 mm Nonlocking Hexalobe Screw
30-0262-S	30-0262	3.5 mm x 22 mm Nonlocking Hexalobe Screw
30-0263-S	30-0263	3.5 mm x 24 mm Nonlocking Hexalobe Screw
30-0264-S	30-0264	3.5 mm x 26 mm Nonlocking Hexalobe Screw
30-0301-S	30-0301	3.0 mm x 8 mm Nonlocking Hexalobe Screw
30-0302-S	30-0302	3.0 mm x 10 mm Nonlocking Hexalobe Screw
30-0303-S	30-0303	3.0 mm x 12 mm Nonlocking Hexalobe Screw
30-0304-S	30-0304	3.0 mm x 14 mm Nonlocking Hexalobe Screw
30-0305-S	30-0305	3.0 mm x 16 mm Nonlocking Hexalobe Screw
30-0306-S	30-0306	3.0 mm x 18 mm Nonlocking Hexalobe Screw
30-0307-S	30-0307	3.0 mm x 20 mm Nonlocking Hexalobe Screw
30-0308-S	30-0308	3.0 mm x 22 mm Nonlocking Hexalobe Screw
30-0309-S	30-0309	3.0 mm x 24 mm Nonlocking Hexalobe Screw
30-0310-S	30-0310	3.0 mm x 26 mm Nonlocking Hexalobe Screw
70-0095-S	70-0095	Locking Clavicle Plate 8-Hole Straight Right
70-0096-S	70-0096	Locking Clavicle Plate 8-Hole Straight Left
CO-2708-S	CO-2708	2.7 mm x 8.0 mm Cortical Screw
CO-2710-S	CO-2710	2.7 mm x 10.0 mm Cortical Screw
CO-2712-S	CO-2712	2.7 mm x 12.0 mm Cortical Screw
CO-2714-S	CO-2714	2.7 mm x 14.0 mm Cortical Screw
CO-2716-S	CO-2716	2.7 mm x 16.0 mm Cortical Screw
CO-2718-S	CO-2718	2.7 mm x 18.0 mm Cortical Screw

Shoulder and Elbow Portfolio [continued]

Locking Clavicle Plating System

Sterile Part Number	Non-Sterile Part Number	Description
CO-2720-S	CO-2720	2.7 mm x 20.0 mm Cortical Screw
CO-2722-S	CO-2722	2.7 mm x 22.0 mm Cortical Screw
CO-2724-S	CO-2724	2.7 mm x 24.0 mm Cortical Screw
CO-2726-S	CO-2726	2.7 mm x 26.0 mm Cortical Screw
CO-2728-S	CO-2728	2.7 mm x 28.0 mm Cortical Screw
CO-2730-S	CO-2730	2.7 mm x 30.0 mm Cortical Screw
CO-2732-S	CO-2732	2.7 mm x 32.0 mm Cortical Screw
CO-2734-S	CO-2734	2.7 mm x 34.0 mm Cortical Screw
CO-2736-S	CO-2736	2.7 mm x 36.0 mm Cortical Screw
CO-2738-S	CO-2738	2.7 mm x 38.0 mm Cortical Screw
CO-2740-S	CO-2740	2.7 mm x 40.0 mm Cortical Screw
CO-2745-S	CO-2745	2.7 mm x 45.0 mm Cortical Screw
CO-2750-S	CO-2750	2.7 mm x 50.0 mm Cortical Screw
CO-2755-S	CO-2755	2.7 mm x 55.0 mm Cortical Screw
CO-2760-S	CO-2760	2.7 mm x 60.0 mm Cortical Screw
CO-2765-S	CO-2765	2.7 mm x 65.0 mm Cortical Screw
CO-3080-S	CO-3080	3.5 mm x 8.0 mm Cortical Screw
CO-3100-S	CO-3100	3.5 mm x 10.0 mm Cortical Screw
CO-3120-S	CO-3120	3.5 mm x 12.0 mm Cortical Screw
CO-3140-S	CO-3140	3.5 mm x 14.0 mm Cortical Screw
CO-3160-S	CO-3160	3.5 mm x 16.0 mm Cortical Screw
CO-3180-S	CO-3180	3.5 mm x 18.0 mm Cortical Screw
CO-3200-S	CO-3200	3.5 mm x 20.0 mm Cortical Screw
CO-3220-S	CO-3220	3.5 mm x 22.0 mm Cortical Screw
CO-3240-S	CO-3240	3.5 mm x 24.0 mm Cortical Screw
CO-3260-S	CO-3260	3.5 mm x 26.0 mm Cortical Screw
COL-3080-S	COL-3080	3.5 mm x 8.0 mm Locking Cortical Screw
COL-3100-S	COL-3100	3.5 mm x 10.0 mm Locking Cortical Screw
COL-3120-S	COL-3120	3.5 mm x 12.0 mm Locking Cortical Screw

Shoulder and Elbow Portfolio [continued]

Locking Clavicle Plating System

Sterile Part Number	Non-Sterile Part Number	Description
COL-3140-S	COL-3140	3.5 mm x 14.0 mm Locking Cortical Screw
COL-3160-S	COL-3160	3.5 mm x 16.0 mm Locking Cortical Screw
COL-3180-S	COL-3180	3.5 mm x 18.0 mm Locking Cortical Screw
COL-3200-S	COL-3200	3.5 mm x 20.0 mm Locking Cortical Screw
COL-3220-S	COL-3220	3.5 mm x 22.0 mm Locking Cortical Screw
COL-3240-S	COL-3240	3.5 mm x 24.0 mm Locking Cortical Screw
COL-3260-S	COL-3260	3.5 mm x 26.0 mm Locking Cortical Screw
PL-CL10LL-S	PL-CL10LL	Locking Clavicle Plate, 10-Hole Large Left
PL-CL10LR-S	PL-CL10LR	Locking Clavicle Plate, 10-Hole Large Right
PL-CL6SL-S	PL-CL6SL	Locking Clavicle Plate, 6-Hole Small Left
PL-CL6SR-S	PL-CL6SR	Locking Clavicle Plate, 6-Hole Small Right
PL-CL8JL-S	PL-CL8JL	Locking Clavicle J-Plate, 8-Hole, Left
PL-CL8JR-S	PL-CL8JR	Locking Clavicle J-Plate, 8-Hole, Right
PL-CL8LL-S	PL-CL8LL	Locking Clavicle Plate, 8-Hole Large Left
PL-CL8LR-S	PL-CL8LR	Locking Clavicle Plate, 8-Hole Large Right
PL-CL8ML-S	PL-CL8ML	Locking Clavicle Plate, 8-Hole Medium Left
PL-CL8MR-S	PL-CL8MR	Locking Clavicle Plate, 8-Hole Medium Right
PL-CL8SL-S	PL-CL8SL	Locking Clavicle Plate, 8-Hole Small Left
PL-CL8SR-S	PL-CL8SR	Locking Clavicle Plate, 8-Hole Small Right
PL-CL9JL-S	PL-CL9JL	Locking Clavicle J-Plate, 9-Hole, Left
PL-CL9JR-S	PL-CL9JR	Locking Clavicle J-Plate, 9-Hole, Right

Shoulder and Elbow Portfolio [continued]

Acu-Sinch® Repair System

Sterile Part Number	Non-Sterile Part Number	Description
N/A	55-0004	Acu-Sinch Anchor
N/A	55-0014	Suture, FlexBraid®, #5 x 28"
N/A	55-0005	Acu-Sinch Suture Retainer

Dual-Trak Clavicle Screw System

Sterile Part Number	Non-Sterile Part Number	Description
40-0136-S	40-0136	3.0 mm x 80 mm Dual-Trak Clavicle Screw
40-0137-S	40-0137	3.0 mm x 90 mm Dual-Trak Clavicle Screw
40-0138-S	40-0138	3.0 mm x 100 mm Dual-Trak Clavicle Screw
40-0139-S	40-0139	3.0 mm x 110 mm Dual-Trak Clavicle Screw
40-0140-S	40-0140	3.0 mm x 120 mm Dual-Trak Clavicle Screw
40-0143-S	40-0143	3.8 mm x 80 mm Dual-Trak Clavicle Screw
40-0144-S	40-0144	3.8 mm x 90 mm Dual-Trak Clavicle Screw
40-0145-S	40-0145	3.8 mm x 100 mm Dual-Trak Clavicle Screw
40-0146-S	40-0146	3.8 mm x 110 mm Dual-Trak Clavicle Screw
40-0147-S	40-0147	3.8 mm x 120 mm Dual-Trak Clavicle Screw

Shoulder and Elbow Portfolio [continued]

Polarus® 3 Solution (plates)

Sterile Part Number	Non-Sterile Part Number	Description
3011-43018-S	3011-43018	4.3 mm x 18 mm Low-Profile Hexalobe Screw
3011-43020-S	3011-43020	4.3 mm x 20 mm Low-Profile Hexalobe Screw
3011-43022-S	3011-43022	4.3 mm x 22 mm Low-Profile Hexalobe Screw
3011-43024-S	3011-43024	4.3 mm x 24 mm Low-Profile Hexalobe Screw
3011-43026-S	3011-43026	4.3 mm x 26 mm Low-Profile Hexalobe Screw
3011-43028-S	3011-43028	4.3 mm x 28 mm Low-Profile Hexalobe Screw
3011-43030-S	3011-43030	4.3 mm x 30 mm Low-Profile Hexalobe Screw
3011-43032-S	3011-43032	4.3 mm x 32 mm Low-Profile Hexalobe Screw
3011-43034-S	3011-43034	4.3 mm x 34 mm Low-Profile Hexalobe Screw
3011-43036-S	3011-43036	4.3 mm x 36 mm Low-Profile Hexalobe Screw
3011-43038-S	3011-43038	4.3 mm x 38 mm Low-Profile Hexalobe Screw
3011-43040-S	3011-43040	4.3 mm x 40 mm Low-Profile Hexalobe Screw
3011-43042-S	3011-43042	4.3 mm x 42 mm Low-Profile Hexalobe Screw
3011-43044-S	3011-43044	4.3 mm x 44 mm Low-Profile Hexalobe Screw
3011-43046-S	3011-43046	4.3 mm x 46 mm Low-Profile Hexalobe Screw
3011-43048-S	3011-43048	4.3 mm x 48 mm Low-Profile Hexalobe Screw
3011-43050-S	3011-43050	4.3 mm x 50 mm Low-Profile Hexalobe Screw
3011-43052-S	3011-43052	4.3 mm x 52 mm Low-Profile Hexalobe Screw
3011-43054-S	3011-43054	4.3 mm x 54 mm Low-Profile Hexalobe Screw
3011-43056-S	3011-43056	4.3 mm x 56 mm Low-Profile Hexalobe Screw
3011-43058-S	3011-43058	4.3 mm x 58 mm Low-Profile Hexalobe Screw
3011-43060-S	3011-43060	4.3 mm x 60 mm Low-Profile Hexalobe Screw
3011-43062-S	3011-43062	4.3 mm x 62 mm Low-Profile Hexalobe Screw
3011-43064-S	3011-43064	4.3 mm x 64 mm Low-Profile Hexalobe Screw
3025-35018-S	3025-35018	3.5 x 18 mm Nonlocking Low-Profile Hexalobe Screw
3025-35020-S	3025-35020	3.5 x 20 mm Nonlocking Low-Profile Hexalobe Screw
3025-35022-S	3025-35022	3.5 x 22 mm Nonlocking Low-Profile Hexalobe Screw
3025-35024-S	3025-35024	3.5 x 24 mm Nonlocking Low-Profile Hexalobe Screw
3025-35026-S	3025-35026	3.5 x 26 mm Nonlocking Low-Profile Hexalobe Screw

Shoulder and Elbow Portfolio [continued]

Polarus® 3 Solution (plates)

Sterile Part Number	Non-Sterile Part Number	Description
3025-35028-S	3025-35028	3.5 x 28 mm Nonlocking Low-Profile Hexalobe Screw
3025-35030-S	3025-35030	3.5 x 30 mm Nonlocking Low-Profile Hexalobe Screw
3025-35032-S	3025-35032	3.5 x 32 mm Nonlocking Low-Profile Hexalobe Screw
3025-35034-S	3025-35034	3.5 x 34 mm Nonlocking Low-Profile Hexalobe Screw
3025-35036-S	3025-35036	3.5 x 36 mm Nonlocking Low-Profile Hexalobe Screw
3025-35038-S	3025-35038	3.5 x 38 mm Nonlocking Low-Profile Hexalobe Screw
3025-35040-S	3025-35040	3.5 x 40 mm Nonlocking Low-Profile Hexalobe Screw
3025-35042-S	3025-35042	3.5 x 42 mm Nonlocking Low-Profile Hexalobe Screw
3025-35044-S	3025-35044	3.5 x 44 mm Nonlocking Low-Profile Hexalobe Screw
3025-35046-S	3025-35046	3.5 x 46 mm Nonlocking Low-Profile Hexalobe Screw
3025-35048-S	3025-35048	3.5 x 48 mm Nonlocking Low-Profile Hexalobe Screw
3025-35050-S	3025-35050	3.5 x 50 mm Nonlocking Low-Profile Hexalobe Screw
3025-35052-S	3025-35052	3.5 x 52 mm Nonlocking Low-Profile Hexalobe Screw
3025-35054-S	3025-35054	3.5 x 54 mm Nonlocking Low-Profile Hexalobe Screw
3025-35056-S	3025-35056	3.5 x 56 mm Nonlocking Low-Profile Hexalobe Screw
3025-35058-S	3025-35058	3.5 x 58 mm Nonlocking Low-Profile Hexalobe Screw
3025-35060-S	3025-35060	3.5 x 60 mm Nonlocking Low-Profile Hexalobe Screw
3025-35062-S	3025-35062	3.5 x 62 mm Nonlocking Low-Profile Hexalobe Screw
3025-35064-S	3025-35064	3.5 x 64 mm Nonlocking Low-Profile Hexalobe Screw
7001-0104L-S	7001-0104L	Polarus 3 Standard Plate 4-Hole Left
7001-0104R-S	7001-0104R	Polarus 3 Standard Plate 4-Hole Right
7001-0106L-S	7001-0106L	Polarus 3 Standard Plate 6-Hole Left
7001-0106R-S	7001-0106R	Polarus 3 Standard Plate 6-Hole Right
7001-0110L-S	7001-0110L	Polarus 3 Standard Plate 10-Hole Left
7001-0110R-S	7001-0110R	Polarus 3 Standard Plate 10-Hole Right
N/A	7001-0114L-S	Polarus 3 Standard Plate 14-Hole Left
N/A	7001-0114R-S	Polarus 3 Standard Plate 14-Hole Right
N/A	7001-0118L-S	Polarus 3 Standard Plate 18-Hole Left
N/A	7001-0118R-S	Polarus 3 Standard Plate 18-Hole Right

Shoulder and Elbow Portfolio [continued]

Polarus® 3 Solution (plates)

Sterile Part Number	Non-Sterile Part Number	Description
N/A	7001-0122L-S	Polarus 3 Standard Plate 22-Hole Left
N/A	7001-0122R-S	Polarus 3 Standard Plate 22-Hole Right
7001-0204L-S	7001-0204L	Polarus 3 Posterior Plate 4-Hole Left
7001-0204R-S	7001-0204R	Polarus 3 Posterior Plate 4-Hole Right
7001-0206L-S	7001-0206L	Polarus 3 Posterior Plate 6-Hole Left
7001-0206R-S	7001-0206R	Polarus 3 Posterior Plate 6-Hole Right
N/A	7001-03001-S	Polarus 3 8 mm Washer, Locking

Polarus® Proximal Humeral Plating System

Sterile Part Number	Non-Sterile Part Number	Description
30-0426-S	30-0426	5.7 x 26 mm Locking Screw
30-0428-S	30-0428	5.7 x 28 mm Locking Screw
30-0430-S	30-0430	5.7 x 30 mm Locking Screw
30-0432-S	30-0432	5.7 x 32 mm Locking Screw
30-0434-S	30-0434	5.7 x 34 mm Locking Screw
30-0436-S	30-0436	5.7 x 36 mm Locking Screw
30-0438-S	30-0438	5.7 x 38 mm Locking Screw
30-0440-S	30-0440	5.7 x 40 mm Locking Screw
30-0442-S	30-0442	5.7 x 42 mm Locking Screw
30-0444-S	30-0444	5.7 x 44 mm Locking Screw
30-0446-S	30-0446	5.7 x 46 mm Locking Screw
30-0448-S	30-0448	5.7 x 48 mm Locking Screw
30-0450-S	30-0450	5.7 x 50 mm Locking Screw
30-0452-S	30-0452	5.7 x 52 mm Locking Screw
30-0454-S	30-0454	5.7 x 54 mm Locking Screw
CA-PHB25-S	CA-PHB25	4.5 mm x 25 mm Locking Buttress Screw
CA-PHB275-S	CA-PHB275	4.5 mm x 27.5 mm Locking Buttress Screw
CA-PHB30-S	CA-PHB30	4.5 mm x 30 mm Locking Buttress Screw
CA-PHB325-S	CA-PHB325	4.5 mm x 32.5 mm Locking Buttress Screw

Shoulder and Elbow Portfolio [continued]

Polarus® Proximal Humeral Plating System

Sterile Part Number	Non-Sterile Part Number	Description
CA-PHB35-S	CA-PHB35	4.5 mm x 35 mm Locking Buttress Screw
CA-PHB375-S	CA-PHB375	4.5 mm x 37.5 mm Locking Buttress Screw
CA-PHB40-S	CA-PHB40	4.5 mm x 40 mm Locking Buttress Screw
CA-PHB45-S	CA-PHB45	4.5 mm x 45 mm Locking Buttress Screw
CA-PHB50-S	CA-PHB50	4.5 mm x 50 mm Locking Buttress Screw
CA-PHB55-S	CA-PHB55	4.5 mm x 55 mm Locking Buttress Screw
CO-3200-S	CO-3200	3.5 mm x 20.0 mm Cortical Screw
CO-3225-S	CO-3225	3.5 mm x 22.5 mm Cortical Screw
CO-3250-S	CO-3250	3.5 mm x 25.0 mm Cortical Screw
CO-3275-S	CO-3275	3.5 mm x 27.5 mm Cortical Screw
CO-3300-S	CO-3300	3.5 mm x 30.0 mm Cortical Screw
CO-3325-S	CO-3325	3.5 mm x 32.5 mm Cortical Screw
CO-3350-S	CO-3350	3.5 mm x 35.0 mm Cortical Screw
CO-3400-S	CO-3400	3.5 mm x 40.0 mm Cortical Screw
CO-3450-S	CO-3450	3.5 mm x 45.0 mm Cortical Screw
CO-3500-S	CO-3500	3.5 mm x 50.0 mm Cortical Screw
CO-3550-S	CO-3550	3.5 mm x 55.0 mm Cortical Screw
COL-3200-S	COL-3200	3.5 mm x 20.0 mm Locking Cortical Screw
COL-3225-S	COL-3225	3.5 mm x 22.5 mm Locking Cortical Screw
COL-3250-S	COL-3250	3.5 mm x 25.0 mm Locking Cortical Screw
COL-3275-S	COL-3275	3.5 mm x 27.5 mm Locking Cortical Screw
COL-3300-S	COL-3300	3.5 mm x 30.0 mm Locking Cortical Screw
N/A	HCA-5125	5.0 mm x 25.0 mm Cancellous Screw
N/A	HCA-5127	5.0 mm x 27.5 mm Cancellous Screw
N/A	HCA-5130	5.0 mm x 30.0 mm Cancellous Screw
N/A	HCA-5132	5.0 mm x 32.5 mm Cancellous Screw
N/A	HCA-5135	5.0 mm x 35.0 mm Cancellous Screw
N/A	HCA-5137	5.0 mm x 37.5 mm Cancellous Screw
N/A	HCA-5140	5.0 mm x 40.0 mm Cancellous Screw

Shoulder and Elbow Portfolio [continued]

Polarus® Proximal Humeral Plating System

Sterile Part Number	Non-Sterile Part Number	Description
N/A	HCA-5145	5.0 mm x 45.0 mm Cancellous Screw
N/A	HCA-5150	5.0 mm x 50.0 mm Cancellous Screw
N/A	HCA-5155	5.0 mm x 55.0 mm Cancellous Screw
PL-PHGL-S	PL-PHGL	Proximal Humeral Plate Large, Left
PL-PHGR-S	PL-PHGR	Proximal Humeral Plate Large, Right
PL-PHSL-S	PL-PHSL	Proximal Humeral Plate, Small, Left
PL-PHSR-S	PL-PHSR	Proximal Humeral Plate, Small, Right
PL-PHXGL-S	PL-PHXGL	Proximal Humeral Plate, Extra Long, Left
PL-PHXGR-S	PL-PHXGR	Proximal Humeral Plate, Extra Long, Right

Scapula Plating System

Sterile Part Number	Non-Sterile Part Number	Description
CA-4120-S	CA-4120	4.0 mm x 12.0 mm Cancellous Screw
CA-4140-S	CA-4140	4.0 mm x 14.0 mm Cancellous Screw
CA-4160-S	CA-4160	4.0 mm x 16.0 mm Cancellous Screw
CA-4180-S	CA-4180	4.0 mm x 18.0 mm Cancellous Screw
CA-4200-S	CA-4200	4.0 mm x 20.0 mm Cancellous Screw
CA-4220-S	CA-4220	4.0 mm x 22.0 mm Cancellous Screw
CA-4240-S	CA-4240	4.0 mm x 24.0 mm Cancellous Screw
CA-4260-S	CA-4260	4.0 mm x 26.0 mm Cancellous Screw
CA-4280-S	CA-4280	4.0 mm x 28.0 mm Cancellous Screw
CA-4300-S	CA-4300	4.0 mm x 30.0 mm Cancellous Screw
CA-4350-S	CA-4350	4.0 mm x 35.0 mm Cancellous Screw
CA-4400-S	CA-4400	4.0 mm x 40.0 mm Cancellous Screw
CA-4450-S	CA-4450	4.0 mm x 45.0 mm Cancellous Screw
CA-4500-S	CA-4500	4.0 mm x 50.0 mm Cancellous Screw
CA-4550-S	CA-4550	4.0 mm x 55.0 mm Cancellous Screw
CA-4600-S	CA-4600	4.0 mm x 60.0 mm Cancellous Screw
CO-2708-S	CO-2708	2.7 mm x 8.0 mm Cortical Screw

Shoulder and Elbow Portfolio [continued]

Scapula Plating System

Sterile Part Number	Non-Sterile Part Number	Description
CO-2710-S	CO-2710	2.7 mm x 10.0 mm Cortical Screw
CO-2712-S	CO-2712	2.7 mm x 12.0 mm Cortical Screw
CO-2714-S	CO-2714	2.7 mm x 14.0 mm Cortical Screw
CO-2716-S	CO-2716	2.7 mm x 16.0 mm Cortical Screw
CO-2718-S	CO-2718	2.7 mm x 18.0 mm Cortical Screw
CO-2720-S	CO-2720	2.7 mm x 20.0 mm Cortical Screw
CO-2722-S	CO-2722	2.7 mm x 22.0 mm Cortical Screw
CO-2724-S	CO-2724	2.7 mm x 24.0 mm Cortical Screw
CO-2726-S	CO-2726	2.7 mm x 26.0 mm Cortical Screw
CO-2728-S	CO-2728	2.7 mm x 28.0 mm Cortical Screw
CO-2730-S	CO-2730	2.7 mm x 30.0 mm Cortical Screw
CO-2732-S	CO-2732	2.7 mm x 32.0 mm Cortical Screw
CO-2734-S	CO-2734	2.7 mm x 34.0 mm Cortical Screw
CO-2736-S	CO-2736	2.7 mm x 36.0 mm Cortical Screw
CO-2738-S	CO-2738	2.7 mm x 38.0 mm Cortical Screw
CO-2740-S	CO-2740	2.7 mm x 40.0 mm Cortical Screw
CO-2745-S	CO-2745	2.7 mm x 45.0 mm Cortical Screw
CO-2750-S	CO-2750	2.7 mm x 50.0 mm Cortical Screw
CO-2755-S	CO-2755	2.7 mm x 55.0 mm Cortical Screw
CO-2760-S	CO-2760	2.7 mm x 60.0 mm Cortical Screw
CO-2765-S	CO-2765	2.7 mm x 65.0 mm Cortical Screw
CO-3060-S	CO-3060	3.5 mm x 6.0 mm Cortical Screw
CO-3080-S	CO-3080	3.5 mm x 8.0 mm Cortical Screw
CO-3100-S	CO-3100	3.5 mm x 10.0 mm Cortical Screw
CO-3120-S	CO-3120	3.5 mm x 12.0 mm Cortical Screw
CO-3140-S	CO-3140	3.5 mm x 14.0 mm Cortical Screw
CO-3160-S	CO-3160	3.5 mm x 16.0 mm Cortical Screw
CO-3180-S	CO-3180	3.5 mm x 18.0 mm Cortical Screw
CO-3200-S	CO-3200	3.5 mm x 20.0 mm Cortical Screw

Shoulder and Elbow Portfolio [continued]

Scapula Plating System

Sterile Part Number	Non-Sterile Part Number	Description
CO-3220-S	CO-3220	3.5 mm x 22.0 mm Cortical Screw
CO-3240-S	CO-3240	3.5 mm x 24.0 mm Cortical Screw
CO-3260-S	CO-3260	3.5 mm x 26.0 mm Cortical Screw
CO-3280-S	CO-3280	3.5 mm x 28.0 mm Cortical Screw
CO-3300-S	CO-3300	3.5 mm x 30.0 mm Cortical Screw
CO-3320-S	CO-3320	3.5 mm x 32.0 mm Cortical Screw
CO-3340-S	CO-3340	3.5 mm x 34.0 mm Cortical Screw
CO-3360-S	CO-3360	3.5 mm x 36.0 mm Cortical Screw
CO-3380-S	CO-3380	3.5 mm x 38.0 mm Cortical Screw
CO-3400-S	CO-3400	3.5 mm x 40.0 mm Cortical Screw
CO-3450-S	CO-3450	3.5 mm x 45.0 mm Cortical Screw
CO-3500-S	CO-3500	3.5 mm x 50.0 mm Cortical Screw
CO-3550-S	CO-3550	3.5 mm x 55.0 mm Cortical Screw
CO-3600-S	CO-3600	3.5 mm x 60.0 mm Cortical Screw
CO-3650-S	CO-3650	3.5 mm x 65.0 mm Cortical Screw
COL-2080-S	COL-2080	2.7 mm x 8.0 mm Locking Cortical Screw
COL-2100-S	COL-2100	2.7 mm x 10.0 mm Locking Cortical Screw
COL-2120-S	COL-2120	2.7 mm x 12.0 mm Locking Cortical Screw
COL-2140-S	COL-2140	2.7 mm x 14.0 mm Locking Cortical Screw
COL-2160-S	COL-2160	2.7 mm x 16.0 mm Locking Cortical Screw
COL-2180-S	COL-2180	2.7 mm x 18.0 mm Locking Cortical Screw
COL-2200-S	COL-2200	2.7 mm x 20.0 mm Locking Cortical Screw
COL-2220-S	COL-2220	2.7 mm x 22.0 mm Locking Cortical Screw
COL-2240-S	COL-2240	2.7 mm x 24.0 mm Locking Cortical Screw
COL-2260-S	COL-2260	2.7 mm x 26.0 mm Locking Cortical Screw
COL-2280-S	COL-2280	2.7 mm x 28.0 mm Locking Cortical Screw
COL-2300-S	COL-2300	2.7 mm x 30.0 mm Locking Cortical Screw
COL-2320-S	COL-2320	2.7 mm x 32.0 mm Locking Cortical Screw
COL-2340-S	COL-2340	2.7 mm x 34.0 mm Locking Cortical Screw

Shoulder and Elbow Portfolio [continued]

Scapula Plating System

Sterile Part Number	Non-Sterile Part Number	Description
COL-2360-S	COL-2360	2.7 mm x 36.0 mm Locking Cortical Screw
COL-2380-S	COL-2380	2.7 mm x 38.0 mm Locking Cortical Screw
COL-2400-S	COL-2400	2.7 mm x 40.0 mm Locking Cortical Screw
COL-2450-S	COL-2450	2.7 mm x 45.0 mm Locking Cortical Screw
COL-2500-S	COL-2500	2.7 mm x 50.0 mm Locking Cortical Screw
COL-2550-S	COL-2550	2.7 mm x 55.0 mm Locking Cortical Screw
COL-2600-S	COL-2600	2.7 mm x 60.0 mm Locking Cortical Screw
COL-2650-S	COL-2650	2.7 mm x 65.0 mm Locking Cortical Screw
COL-3060-S	COL-3060	3.5 mm x 6.0 mm Locking Cortical Screw
COL-3080-S	COL-3080	3.5 mm x 8.0 mm Locking Cortical Screw
COL-3100-S	COL-3100	3.5 mm x 10.0 mm Locking Cortical Screw
COL-3120-S	COL-3120	3.5 mm x 12.0 mm Locking Cortical Screw
COL-3140-S	COL-3140	3.5 mm x 14.0 mm Locking Cortical Screw
COL-3160-S	COL-3160	3.5 mm x 16.0 mm Locking Cortical Screw
COL-3180-S	COL-3180	3.5 mm x 18.0 mm Locking Cortical Screw
COL-3200-S	COL-3200	3.5 mm x 20.0 mm Locking Cortical Screw
COL-3220-S	COL-3220	3.5 mm x 22.0 mm Locking Cortical Screw
COL-3240-S	COL-3240	3.5 mm x 24.0 mm Locking Cortical Screw
COL-3260-S	COL-3260	3.5 mm x 26.0 mm Locking Cortical Screw
COL-3280-S	COL-3280	3.5 mm x 28.0 mm Locking Cortical Screw
COL-3300-S	COL-3300	3.5 mm x 30.0 mm Locking Cortical Screw
COL-3320-S	COL-3320	3.5 mm x 32.0 mm Locking Cortical Screw
COL-3340-S	COL-3340	3.5 mm x 34.0 mm Locking Cortical Screw
COL-3360-S	COL-3360	3.5 mm x 36.0 mm Locking Cortical Screw
COL-3380-S	COL-3380	3.5 mm x 38.0 mm Locking Cortical Screw
COL-3400-S	COL-3400	3.5 mm x 40.0 mm Locking Cortical Screw
COL-3450-S	COL-3450	3.5 mm x 45.0 mm Locking Cortical Screw
COL-3500-S	COL-3500	3.5 mm x 50.0 mm Locking Cortical Screw
COL-3550-S	COL-3550	3.5 mm x 55.0 mm Locking Cortical Screw

Shoulder and Elbow Portfolio [continued]

Scapula Plating System

Sterile Part Number	Non-Sterile Part Number	Description
COL-3600-S	COL-3600	3.5 mm x 60.0 mm Locking Cortical Screw
COL-3650-S	COL-3650	3.5 mm x 65.0 mm Locking Cortical Screw
PL-SA06L-S	PL-SA06L	6-Hole Scapula Acromion Plate, Left
PL-SA06R-S	PL-SA06R	6-Hole Scapula Acromion Plate, Right
PL-SA07L-S	PL-SA07L	7-Hole Scapula Acromion Plate, Left
PL-SA07R-S	PL-SA07R	7-Hole Scapula Acromion Plate, Right
PL-SG04L-S	PL-SG04L	4-Hole Scapula Glenoid Plate, Left
PL-SG04R-S	PL-SG04R	4-Hole Scapula Glenoid Plate, Right
PL-SLB10L-S	PL-SLB10L	10-Hole Scapula Lateral Border Plate, Left
PL-SLB10R-S	PL-SLB10R	10-Hole Scapula Lateral Border Plate, Right
PL-SMB09L-S	PL-SMB09L	9-Hole Scapula Medial Border Plate, Left
PL-SMB09R-S	PL-SMB09R	9-Hole Scapula Medial Border Plate, Right
PL-SMB13L-S	PL-SMB13L	13-Hole Scapula Medial Border Plate, Left
PL-SMB13R-S	PL-SMB13R	13-Hole Scapula Medial Border Plate, Right

Shoulder and Elbow Portfolio [continued]

Conditions for safe scanning in the MRI

A patient with the screw, plate, pin, suture anchor, suture retainer, suture, and washer implants from the Elbow Plating, Congruent Elbow Plating, Radial Head Plating, Clavicle Plating, Locking Clavicle Plating, Clavicle Hook Plating, Acu-Sinch® Repair, Dual-Trak Clavicle Screw, Polarus® 3 Solution, Polarus Proximal Humeral Plating, and Scapula Plating Systems may be safely scanned under the following conditions. Failure to follow these conditions may result in injury to the patient.

Nominal values(s) of Static Magnetic Field	1.5 T or 3 T
Maximum Spatial Field Gradient	30 T/m (3000 gauss/cm)
RF Excitation	Circularly Polarized (CP)
RF Transmit Coil Type	Whole body transmit coil, Head RF transmit-receive coil
Maximum Whole Body SAR	1.0 W/kg*
Limits on Scan Duration	1.0 W/kg whole body average SAR for 20 minutes of continuous RF (a sequence or back to back series/scan without breaks) followed by a 3 minute cooling period for a one hour session.
MR Image Artifact	The presence of this implant may produce an image artifact of 37 mm.
If more information about a specific parameter is not included, there are no conditions associated with that parameter.	

* Whole body SAR of 1.0 W/kg may limit the availability and diagnostic quality of some MR procedures

Foot and Ankle Portfolio

- ▶ Ankle Plating System
- ▶ Ankle Plating System 3
- ▶ Calcaneal Plating System
- ▶ Forefoot-Midfoot Plating System
- ▶ Small Fragment Base Set

Ankle Plating System

Sterile Part Number	Non-Sterile Part Number	Description
30-0232-S	30-0232	3.5 mm x 8 mm Locking Hexalobe Screw
30-0233-S	30-0233	3.5 mm x 10 mm Locking Hexalobe Screw
30-0234-S	30-0234	3.5 mm x 12 mm Locking Hexalobe Screw
30-0235-S	30-0235	3.5 mm x 14 mm Locking Hexalobe Screw
30-0236-S	30-0236	3.5 mm x 16 mm Locking Hexalobe Screw
30-0237-S	30-0237	3.5 mm x 18 mm Locking Hexalobe Screw
30-0238-S	30-0238	3.5 mm x 20 mm Locking Hexalobe Screw
30-0239-S	30-0239	3.5 mm x 22 mm Locking Hexalobe Screw
30-0240-S	30-0240	3.5 mm x 24 mm Locking Hexalobe Screw
30-0241-S	30-0241	3.5 mm x 26 mm Locking Hexalobe Screw
30-0242-S	30-0242	3.5 mm x 28 mm Locking Hexalobe Screw
30-0243-S	30-0243	3.5 mm x 30 mm Locking Hexalobe Screw
30-0244-S	30-0244	3.5 mm x 32 mm Locking Hexalobe Screw
30-0245-S	30-0245	3.5 mm x 34 mm Locking Hexalobe Screw
30-0246-S	30-0246	3.5 mm x 36 mm Locking Hexalobe Screw
30-0247-S	30-0247	3.5 mm x 38 mm Locking Hexalobe Screw
30-0248-S	30-0248	3.5 mm x 40 mm Locking Hexalobe Screw
30-0249-S	30-0249	3.5 mm x 45 mm Locking Hexalobe Screw
30-0250-S	30-0250	3.5 mm x 50 mm Locking Hexalobe Screw
30-0251-S	30-0251	3.5 mm x 55 mm Locking Hexalobe Screw
30-0252-S	30-0252	3.5 mm x 60 mm Locking Hexalobe Screw
30-0255-S	30-0255	3.5 mm x 8 mm Nonlocking Hexalobe Screw
30-0256-S	30-0256	3.5 mm x 10 mm Nonlocking Hexalobe Screw
30-0257-S	30-0257	3.5 mm x 12 mm Nonlocking Hexalobe Screw
30-0258-S	30-0258	3.5 mm x 14 mm Nonlocking Hexalobe Screw
30-0259-S	30-0259	3.5 mm x 16 mm Nonlocking Hexalobe Screw
30-0260-S	30-0260	3.5 mm x 18 mm Nonlocking Hexalobe Screw

Foot and Ankle Portfolio [continued]

Ankle Plating System

Sterile Part Number	Non-Sterile Part Number	Description
30-0261-S	30-0261	3.5 mm x 20 mm Nonlocking Hexalobe Screw
30-0262-S	30-0262	3.5 mm x 22 mm Nonlocking Hexalobe Screw
30-0263-S	30-0263	3.5 mm x 24 mm Nonlocking Hexalobe Screw
30-0264-S	30-0264	3.5 mm x 26 mm Nonlocking Hexalobe Screw
30-0265-S	30-0265	3.5 mm x 28 mm Nonlocking Hexalobe Screw
30-0266-S	30-0266	3.5 mm x 30 mm Nonlocking Hexalobe Screw
30-0267-S	30-0267	3.5 mm x 32 mm Nonlocking Hexalobe Screw
30-0268-S	30-0268	3.5 mm x 34 mm Nonlocking Hexalobe Screw
30-0269-S	30-0269	3.5 mm x 36 mm Nonlocking Hexalobe Screw
30-0270-S	30-0270	3.5 mm x 38 mm Nonlocking Hexalobe Screw
30-0271-S	30-0271	3.5 mm x 40 mm Nonlocking Hexalobe Screw
30-0272-S	30-0272	3.5 mm x 45 mm Nonlocking Hexalobe Screw
30-0273-S	30-0273	3.5 mm x 50 mm Nonlocking Hexalobe Screw
30-0274-S	30-0274	3.5 mm x 55 mm Nonlocking Hexalobe Screw
30-0275-S	30-0275	3.5 mm x 60 mm Nonlocking Hexalobe Screw
30-0278-S	30-0278	3.0 mm x 8 mm Locking Hexalobe Screw
30-0279-S	30-0279	3.0 mm x 10 mm Locking Hexalobe Screw
30-0280-S	30-0280	3.0 mm x 12 mm Locking Hexalobe Screw
30-0281-S	30-0281	3.0 mm x 14 mm Locking Hexalobe Screw
30-0282-S	30-0282	3.0 mm x 16 mm Locking Hexalobe Screw
30-0283-S	30-0283	3.0 mm x 18 mm Locking Hexalobe Screw
30-0284-S	30-0284	3.0 mm x 20 mm Locking Hexalobe Screw
30-0285-S	30-0285	3.0 mm x 22 mm Locking Hexalobe Screw
30-0286-S	30-0286	3.0 mm x 24 mm Locking Hexalobe Screw
30-0287-S	30-0287	3.0 mm x 26 mm Locking Hexalobe Screw
30-0288-S	30-0288	3.0 mm x 28 mm Locking Hexalobe Screw
30-0289-S	30-0289	3.0 mm x 30 mm Locking Hexalobe Screw
30-0290-S	30-0290	3.0 mm x 32 mm Locking Hexalobe Screw
30-0291-S	30-0291	3.0 mm x 34 mm Locking Hexalobe Screw

Foot and Ankle Portfolio [continued]

Ankle Plating System

Sterile Part Number	Non-Sterile Part Number	Description
30-0292-S	30-0292	3.0 mm x 36 mm Locking Hexalobe Screw
30-0293-S	30-0293	3.0 mm x 38 mm Locking Hexalobe Screw
30-0294-S	30-0294	3.0 mm x 40 mm Locking Hexalobe Screw
30-0295-S	30-0295	3.0 mm x 45 mm Locking Hexalobe Screw
30-0296-S	30-0296	3.0 mm x 50 mm Locking Hexalobe Screw
30-0297-S	30-0297	3.0 mm x 55 mm Locking Hexalobe Screw
30-0301-S	30-0301	3.0 mm x 8 mm Nonlocking Hexalobe Screw
30-0302-S	30-0302	3.0 mm x 10 mm Nonlocking Hexalobe Screw
30-0303-S	30-0303	3.0 mm x 12 mm Nonlocking Hexalobe Screw
30-0304-S	30-0304	3.0 mm x 14 mm Nonlocking Hexalobe Screw
30-0305-S	30-0305	3.0 mm x 16 mm Nonlocking Hexalobe Screw
30-0306-S	30-0306	3.0 mm x 18 mm Nonlocking Hexalobe Screw
30-0307-S	30-0307	3.0 mm x 20 mm Nonlocking Hexalobe Screw
30-0308-S	30-0308	3.0 mm x 22 mm Nonlocking Hexalobe Screw
30-0309-S	30-0309	3.0 mm x 24 mm Nonlocking Hexalobe Screw
30-0310-S	30-0310	3.0 mm x 26 mm Nonlocking Hexalobe Screw
30-0311-S	30-0311	3.0 mm x 28 mm Nonlocking Hexalobe Screw
30-0312-S	30-0312	3.0 mm x 30 mm Nonlocking Hexalobe Screw
30-0313-S	30-0313	3.0 mm x 32 mm Nonlocking Hexalobe Screw
30-0314-S	30-0314	3.0 mm x 34 mm Nonlocking Hexalobe Screw
30-0315-S	30-0315	3.0 mm x 36 mm Nonlocking Hexalobe Screw
30-0316-S	30-0316	3.0 mm x 38 mm Nonlocking Hexalobe Screw
30-0317-S	30-0317	3.0 mm x 40 mm Nonlocking Hexalobe Screw
30-0318-S	30-0318	3.0 mm x 45 mm Nonlocking Hexalobe Screw
30-0319-S	30-0319	3.0 mm x 50 mm Nonlocking Hexalobe Screw
30-0320-S	30-0320	3.0 mm x 55 mm Nonlocking Hexalobe Screw
70-0145-S	70-0145	Low-profile Locking Lateral Fibula Plate 5-Hole
70-0147-S	70-0147	Low-profile Locking Lateral Fibula Plate 7-Hole
70-0149-S	70-0149	Low-profile Locking Lateral Fibula Plate 9-Hole

Foot and Ankle Portfolio [continued]

Ankle Plating System

Sterile Part Number	Non-Sterile Part Number	Description
70-0151-S	70-0151	Low-profile Locking Lateral Fibula Plate 11-Hole
70-0153-S	70-0153	Low-profile Locking Lateral Fibula Plate 13-Hole
70-0169-S	70-0169	Locking Lateral Fibula Plate 9-Hole
70-0171-S	70-0171	Locking Lateral Fibula Plate 11-Hole
70-0173-S	70-0173	Locking Lateral Fibula Plate 13-Hole
70-0227-S	70-0227	Low-profile Locking Medial Tibia Plate 7-Hole
70-0229-S	70-0229	Low-profile Locking Medial Tibia Plate 9-Hole
70-0245-S	70-0245	Low-profile Locking Anterior Tibia Plate 5-Hole
70-0247-S	70-0247	Low-profile Locking Anterior Tibia Plate 7-Hole
CA-4120-S	CA-4120	4.0 mm x 12.0 mm Cancellous Screw
CA-4140-S	CA-4140	4.0 mm x 14.0 mm Cancellous Screw
CA-4160-S	CA-4160	4.0 mm x 16.0 mm Cancellous Screw
CA-4180-S	CA-4180	4.0 mm x 18.0 mm Cancellous Screw
CA-4200-S	CA-4200	4.0 mm x 20.0 mm Cancellous Screw
CA-4220-S	CA-4220	4.0 mm x 22.0 mm Cancellous Screw
CA-4240-S	CA-4240	4.0 mm x 24.0 mm Cancellous Screw
CA-4260-S	CA-4260	4.0 mm x 26.0 mm Cancellous Screw
CA-4280-S	CA-4280	4.0 mm x 28.0 mm Cancellous Screw
CA-4300-S	CA-4300	4.0 mm x 30.0 mm Cancellous Screw
CA-4350-S	CA-4350	4.0 mm x 35.0 mm Cancellous Screw
CA-4400-S	CA-4400	4.0 mm x 40.0 mm Cancellous Screw
CA-4450-S	CA-4450	4.0 mm x 45.0 mm Cancellous Screw
CA-4500-S	CA-4500	4.0 mm x 50.0 mm Cancellous Screw
CA-4550-S	CA-4550	4.0 mm x 55.0 mm Cancellous Screw
CA-4600-S	CA-4600	4.0 mm x 60.0 mm Cancellous Screw
CO-2708-S	CO-2708	2.7 mm x 8.0 mm Cortical Screw
CO-2710-S	CO-2710	2.7 mm x 10.0 mm Cortical Screw
CO-2712-S	CO-2712	2.7 mm x 12.0 mm Cortical Screw
CO-2714-S	CO-2714	2.7 mm x 14.0 mm Cortical Screw

Foot and Ankle Portfolio [continued]

Ankle Plating System

Sterile Part Number	Non-Sterile Part Number	Description
CO-2716-S	CO-2716	2.7 mm x 16.0 mm Cortical Screw
CO-2718-S	CO-2718	2.7 mm x 18.0 mm Cortical Screw
CO-2720-S	CO-2720	2.7 mm x 20.0 mm Cortical Screw
CO-2722-S	CO-2722	2.7 mm x 22.0 mm Cortical Screw
CO-2724-S	CO-2724	2.7 mm x 24.0 mm Cortical Screw
CO-2726-S	CO-2726	2.7 mm x 26.0 mm Cortical Screw
CO-2728-S	CO-2728	2.7 mm x 28.0 mm Cortical Screw
CO-2730-S	CO-2730	2.7 mm x 30.0 mm Cortical Screw
CO-2732-S	CO-2732	2.7 mm x 32.0 mm Cortical Screw
CO-2734-S	CO-2734	2.7 mm x 34.0 mm Cortical Screw
CO-2736-S	CO-2736	2.7 mm x 36.0 mm Cortical Screw
CO-2738-S	CO-2738	2.7 mm x 38.0 mm Cortical Screw
CO-2740-S	CO-2740	2.7 mm x 40.0 mm Cortical Screw
CO-2745-S	CO-2745	2.7 mm x 45.0 mm Cortical Screw
CO-2750-S	CO-2750	2.7 mm x 50.0 mm Cortical Screw
CO-2755-S	CO-2755	2.7 mm x 55.0 mm Cortical Screw
CO-3080-S	CO-3080	3.5 mm x 8.0 mm Cortical Screw
CO-3100-S	CO-3100	3.5 mm x 10.0 mm Cortical Screw
CO-3120-S	CO-3120	3.5 mm x 12.0 mm Cortical Screw
CO-3140-S	CO-3140	3.5 mm x 14.0 mm Cortical Screw
CO-3160-S	CO-3160	3.5 mm x 16.0 mm Cortical Screw
CO-3180-S	CO-3180	3.5 mm x 18.0 mm Cortical Screw
CO-3200-S	CO-3200	3.5 mm x 20.0 mm Cortical Screw
CO-3220-S	CO-3220	3.5 mm x 22.0 mm Cortical Screw
CO-3240-S	CO-3240	3.5 mm x 24.0 mm Cortical Screw
CO-3260-S	CO-3260	3.5 mm x 26.0 mm Cortical Screw
CO-3280-S	CO-3280	3.5 mm x 28.0 mm Cortical Screw
CO-3300-S	CO-3300	3.5 mm x 30.0 mm Cortical Screw
CO-3320-S	CO-3320	3.5 mm x 32.0 mm Cortical Screw

Foot and Ankle Portfolio [continued]

Ankle Plating System

Sterile Part Number	Non-Sterile Part Number	Description
CO-3340-S	CO-3340	3.5 mm x 34.0 mm Cortical Screw
CO-3360-S	CO-3360	3.5 mm x 36.0 mm Cortical Screw
CO-3380-S	CO-3380	3.5 mm x 38.0 mm Cortical Screw
CO-3400-S	CO-3400	3.5 mm x 40.0 mm Cortical Screw
CO-3450-S	CO-3450	3.5 mm x 45.0 mm Cortical Screw
CO-3500-S	CO-3500	3.5 mm x 50.0 mm Cortical Screw
CO-3550-S	CO-3550	3.5 mm x 55.0 mm Cortical Screw
CO-3600-S	CO-3600	3.5 mm x 60.0 mm Cortical Screw
COL-2080-S	COL-2080	2.7 mm x 8.0 mm Locking Cortical Screw
COL-2100-S	COL-2100	2.7 mm x 10.0 mm Locking Cortical Screw
COL-2120-S	COL-2120	2.7 mm x 12.0 mm Locking Cortical Screw
COL-2140-S	COL-2140	2.7 mm x 14.0 mm Locking Cortical Screw
COL-2160-S	COL-2160	2.7 mm x 16.0 mm Locking Cortical Screw
COL-2180-S	COL-2180	2.7 mm x 18.0 mm Locking Cortical Screw
COL-2200-S	COL-2200	2.7 mm x 20.0 mm Locking Cortical Screw
COL-2220-S	COL-2220	2.7 mm x 22.0 mm Locking Cortical Screw
COL-2240-S	COL-2240	2.7 mm x 24.0 mm Locking Cortical Screw
COL-2260-S	COL-2260	2.7 mm x 26.0 mm Locking Cortical Screw
COL-2280-S	COL-2280	2.7 mm x 28.0 mm Locking Cortical Screw
COL-2300-S	COL-2300	2.7 mm x 30.0 mm Locking Cortical Screw
COL-2320-S	COL-2320	2.7 mm x 32.0 mm Locking Cortical Screw
COL-2340-S	COL-2340	2.7 mm x 34.0 mm Locking Cortical Screw
COL-2360-S	COL-2360	2.7 mm x 36.0 mm Locking Cortical Screw
COL-2380-S	COL-2380	2.7 mm x 38.0 mm Locking Cortical Screw
COL-2400-S	COL-2400	2.7 mm x 40.0 mm Locking Cortical Screw
COL-2450-S	COL-2450	2.7 mm x 45.0 mm Locking Cortical Screw
COL-2500-S	COL-2500	2.7 mm x 50.0 mm Locking Cortical Screw
COL-2550-S	COL-2550	2.7 mm x 55.0 mm Locking Cortical Screw
COL-3080-S	COL-3080	3.5 mm x 8.0 mm Locking Cortical Screw

Foot and Ankle Portfolio [continued]

Ankle Plating System

Sterile Part Number	Non-Sterile Part Number	Description
COL-3100-S	COL-3100	3.5 mm x 10.0 mm Locking Cortical Screw
COL-3120-S	COL-3120	3.5 mm x 12.0 mm Locking Cortical Screw
COL-3140-S	COL-3140	3.5 mm x 14.0 mm Locking Cortical Screw
COL-3160-S	COL-3160	3.5 mm x 16.0 mm Locking Cortical Screw
COL-3180-S	COL-3180	3.5 mm x 18.0 mm Locking Cortical Screw
COL-3200-S	COL-3200	3.5 mm x 20.0 mm Locking Cortical Screw
COL-3220-S	COL-3220	3.5 mm x 22.0 mm Locking Cortical Screw
COL-3240-S	COL-3240	3.5 mm x 24.0 mm Locking Cortical Screw
COL-3260-S	COL-3260	3.5 mm x 26.0 mm Locking Cortical Screw
COL-3280-S	COL-3280	3.5 mm x 28.0 mm Locking Cortical Screw
COL-3300-S	COL-3300	3.5 mm x 30.0 mm Locking Cortical Screw
COL-3320-S	COL-3320	3.5 mm x 32.0 mm Locking Cortical Screw
COL-3340-S	COL-3340	3.5 mm x 34.0 mm Locking Cortical Screw
COL-3360-S	COL-3360	3.5 mm x 36.0 mm Locking Cortical Screw
COL-3380-S	COL-3380	3.5 mm x 38.0 mm Locking Cortical Screw
COL-3400-S	COL-3400	3.5 mm x 40.0 mm Locking Cortical Screw
COL-3450-S	COL-3450	3.5 mm x 45.0 mm Locking Cortical Screw
COL-3500-S	COL-3500	3.5 mm x 50.0 mm Locking Cortical Screw
COL-3550-S	COL-3550	3.5 mm x 55.0 mm Locking Cortical Screw
COL-3600-S	COL-3600	3.5 mm x 60.0 mm Locking Cortical Screw

Ankle Plating System 3

Sterile Part Number	Non-Sterile Part Number	Description
N/A	30-0097	50.0 mm Tension Band Pin
N/A	30-0098	70.0 mm Tension Band Pin
N/A	30-0099	90.0 mm Tension Band Pin
N/A	30-0166	2 mm x 50 mm Tension Band Pin
N/A	30-0167	2 mm x 70 mm Tension Band Pin
N/A	30-0168	2 mm x 90 mm Tension Band Pin

Foot and Ankle Portfolio [continued]

Ankle Plating System 3

Sterile Part Number	Non-Sterile Part Number	Description
N/A	3005-40010	4.0 mm x 10 mm Cannulated Screw, Short Thread
N/A	3005-40012	4.0 mm x 12 mm Cannulated Screw, Short Thread
N/A	3005-40014	4.0 mm x 14 mm Cannulated Screw, Short Thread
N/A	3005-40016	4.0 mm x 16 mm Cannulated Screw, Short Thread
N/A	3005-40018	4.0 mm x 18 mm Cannulated Screw, Short Thread
N/A	3005-40020	4.0 mm x 20 mm Cannulated Screw, Short Thread
N/A	3005-40022	4.0 mm x 22 mm Cannulated Screw, Short Thread
N/A	3005-40024	4.0 mm x 24 mm Cannulated Screw, Short Thread
N/A	3005-40026	4.0 mm x 26 mm Cannulated Screw, Short Thread
N/A	3005-40028	4.0 mm x 28 mm Cannulated Screw, Short Thread
N/A	3005-40030	4.0 mm x 30 mm Cannulated Screw, Short Thread
N/A	3005-40032	4.0 mm x 32 mm Cannulated Screw, Short Thread
N/A	3005-40034	4.0 mm x 34 mm Cannulated Screw, Short Thread
N/A	3005-40036	4.0 mm x 36 mm Cannulated Screw, Short Thread
N/A	3005-40038	4.0 mm x 38 mm Cannulated Screw, Short Thread
N/A	3005-40040	4.0 mm x 40 mm Cannulated Screw, Short Thread
N/A	3005-40042	4.0 mm x 42 mm Cannulated Screw, Short Thread
N/A	3005-40044	4.0 mm x 44 mm Cannulated Screw, Short Thread
N/A	3005-40046	4.0 mm x 46 mm Cannulated Screw, Short Thread
N/A	3005-40048	4.0 mm x 48 mm Cannulated Screw, Short Thread
N/A	3005-40050	4.0 mm x 50 mm Cannulated Screw, Short Thread
N/A	3005-40052	4.0 mm x 52 mm Cannulated Screw, Short Thread
N/A	3005-40054	4.0 mm x 54 mm Cannulated Screw, Short Thread
N/A	3005-40056	4.0 mm x 56 mm Cannulated Screw, Short Thread
N/A	3005-40058	4.0 mm x 58 mm Cannulated Screw, Short Thread
N/A	3005-40060	4.0 mm x 60 mm Cannulated Screw, Short Thread
N/A	3005-40064	4.0 mm x 64 mm Cannulated Screw, Short Thread
N/A	3005-40068	4.0 mm x 68 mm Cannulated Screw, Short Thread
N/A	3005-40072	4.0 mm x 72 mm Cannulated Screw, Short Thread

Foot and Ankle Portfolio [continued]

Ankle Plating System 3

Sterile Part Number	Non-Sterile Part Number	Description
N/A	3006-40016	4.0 mm x 16 mm Cannulated Screw, Long Thread
N/A	3006-40018	4.0 mm x 18 mm Cannulated Screw, Long Thread
N/A	3006-40020	4.0 mm x 20 mm Cannulated Screw, Long Thread
N/A	3006-40022	4.0 mm x 22 mm Cannulated Screw, Long Thread
N/A	3006-40024	4.0 mm x 24 mm Cannulated Screw, Long Thread
N/A	3006-40026	4.0 mm x 26 mm Cannulated Screw, Long Thread
N/A	3006-40028	4.0 mm x 28 mm Cannulated Screw, Long Thread
N/A	3006-40030	4.0 mm x 30 mm Cannulated Screw, Long Thread
N/A	3006-40032	4.0 mm x 32 mm Cannulated Screw, Long Thread
N/A	3006-40034	4.0 mm x 34 mm Cannulated Screw, Long Thread
3006-40036-S	3006-40036	4.0 mm x 36 mm Cannulated Screw, Long Thread
N/A	3006-40038	4.0 mm x 38 mm Cannulated Screw, Long Thread
N/A	3006-40040	4.0 mm x 40 mm Cannulated Screw, Long Thread
3006-40042-S	3006-40042	4.0 mm x 42 mm Cannulated Screw, Long Thread
N/A	3006-40044	4.0 mm x 44 mm Cannulated Screw, Long Thread
N/A	3006-40046	4.0 mm x 46 mm Cannulated Screw, Long Thread
3006-40048-S	3006-40048	4.0 mm x 48 mm Cannulated Screw, Long Thread
N/A	3006-40050	4.0 mm x 50 mm Cannulated Screw, Long Thread
N/A	3006-40052	4.0 mm x 52 mm Cannulated Screw, Long Thread
N/A	3006-40054	4.0 mm x 54 mm Cannulated Screw, Long Thread
N/A	3006-40056	4.0 mm x 56 mm Cannulated Screw, Long Thread
N/A	3006-40058	4.0 mm x 58 mm Cannulated Screw, Long Thread
N/A	3006-40060	4.0 mm x 60 mm Cannulated Screw, Long Thread
N/A	3006-40064	4.0 mm x 64 mm Cannulated Screw, Long Thread
N/A	3006-40068	4.0 mm x 68 mm Cannulated Screw, Long Thread
N/A	3006-40072	4.0 mm x 72 mm Cannulated Screw, Long Thread
3014-23026-S	3014-23026	2.3 mm x 26 mm Locking Cortical Peg
3014-23036-S	3014-23036	2.3 mm x 36 mm Locking Cortical Peg
3014-23046-S	3014-23046	2.3 mm x 46 mm Locking Cortical Peg

Foot and Ankle Portfolio [continued]

Ankle Plating System 3

Sterile Part Number	Non-Sterile Part Number	Description
7003-07036-S	7003-07036	Cannulated Screw Washer 7.0 mm OD x 3.6 mm ID
7007-0104L-S	7007-0104L	Lateral Fibula Plate 4-Hole Left
7007-0104R-S	7007-0104R	Lateral Fibula Plate 4-Hole Right
7007-0105L-S	7007-0105L	Lateral Fibula Plate 5-Hole Left
7007-0105R-S	7007-0105R	Lateral Fibula Plate 5-Hole Right
7007-0106L-S	7007-0106L	Lateral Fibula Plate 6-Hole Left
7007-0106R-S	7007-0106R	Lateral Fibula Plate 6-Hole Right
7007-0107L-S	7007-0107L	Lateral Fibula Plate 7-Hole Left
7007-0107R-S	7007-0107R	Lateral Fibula Plate 7-Hole Right
7007-0109L-S	7007-0109L	Lateral Fibula Plate 9-Hole Left
7007-0109R-S	7007-0109R	Lateral Fibula Plate 9-Hole Right
7007-0111L-S	N/A	Lateral Fibula Plate 11-Hole Left
7007-0111R-S	N/A	Lateral Fibula Plate 11-Hole Right
7007-0113L-S	N/A	Lateral Fibula Plate 13-Hole Left
7007-0113R-S	N/A	Lateral Fibula Plate 13-Hole Right
7007-0203L-S	7007-0203L	Posterolateral Fibula Plate 3-Hole Left
7007-0203R-S	7007-0203R	Posterolateral Fibula Plate 3-Hole Right
7007-0204L-S	7007-0204L	Posterolateral Fibula Plate 4-Hole Left
7007-0204R-S	7007-0204R	Posterolateral Fibula Plate 4-Hole Right
7007-0205L-S	7007-0205L	Posterolateral Fibula Plate 5-Hole Left
7007-0205R-S	7007-0205R	Posterolateral Fibula Plate 5-Hole Right
7007-0206L-S	7007-0206L	Posterolateral Fibula Plate 6-Hole Left
7007-0206R-S	7007-0206R	Posterolateral Fibula Plate 6-Hole Right
7007-0207L-S	7007-0207L	Posterolateral Fibula Plate 7-Hole Left
7007-0207R-S	7007-0207R	Posterolateral Fibula Plate 7-Hole Right
7007-0303L-S	7007-0303L	Posterolateral Distal Tibia Plate 3-Hole Left
7007-0303R-S	7007-0303R	Posterolateral Distal Tibia Plate 3-Hole Right
7007-0304L-S	7007-0304L	Posterolateral Distal Tibia Plate 4-Hole Left
7007-0304R-S	7007-0304R	Posterolateral Distal Tibia Plate 4-Hole Right

Foot and Ankle Portfolio [continued]

Ankle Plating System 3

Sterile Part Number	Non-Sterile Part Number	Description
7007-0403L-S	7007-0403L	Posteromedial Distal Tibia Plate 3-Hole Left
7007-0403R-S	7007-0403R	Posteromedial Distal Tibia Plate 3-Hole Right
7007-0504-S	7007-0504	Medial Anti-Glide Plate 4-Hole
7007-0602-S	7007-0602	Hook Plate 2-Hole
7007-0603-S	7007-0603	Hook Plate 3-Hole
7007-0702-S	7007-0702	Locking Peg Hook Plate 2-Hole
7007-0703-S	7007-0703	Locking Peg Hook Plate 3-Hole
AI-0010-S	AI-0010	10 mm AcuTwist® Acutrak Compression Screw
AI-0012-S	AI-0012	12 mm AcuTwist Acutrak Compression Screw
AI-0014-S	AI-0014	14 mm AcuTwist Acutrak Compression Screw
AI-0016-S	AI-0016	16 mm AcuTwist Acutrak Compression Screw
AI-0018-S	AI-0018	18 mm AcuTwist Acutrak Compression Screw
AI-0020-S	AI-0020	20 mm AcuTwist Acutrak Compression Screw
AI-0022-S	AI-0022	22 mm AcuTwist Acutrak Compression Screw
AI-0024-S	AI-0024	24 mm AcuTwist Acutrak Compression Screw
AI-0026-S	AI-0026	26 mm AcuTwist Acutrak Compression Screw
AI-0028-S	AI-0028	28 mm AcuTwist Acutrak Compression Screw
AI-0030-S	AI-0030	30 mm AcuTwist Acutrak Compression Screw

Calcaneal Plating System

Sterile Part Number	Non-Sterile Part Number	Description
30-0232-S	30-0232	3.5 mm x 8 mm Locking Hexalobe Screw
30-0233-S	30-0233	3.5 mm x 10 mm Locking Hexalobe Screw
30-0234-S	30-0234	3.5 mm x 12 mm Locking Hexalobe Screw
30-0235-S	30-0235	3.5 mm x 14 mm Locking Hexalobe Screw
30-0236-S	30-0236	3.5 mm x 16 mm Locking Hexalobe Screw
30-0237-S	30-0237	3.5 mm x 18 mm Locking Hexalobe Screw
30-0238-S	30-0238	3.5 mm x 20 mm Locking Hexalobe Screw
30-0239-S	30-0239	3.5 mm x 22 mm Locking Hexalobe Screw

Foot and Ankle Portfolio [continued]

Calcaneal Plating System

Sterile Part Number	Non-Sterile Part Number	Description
30-0240-S	30-0240	3.5 mm x 24 mm Locking Hexalobe Screw
30-0241-S	30-0241	3.5 mm x 26 mm Locking Hexalobe Screw
30-0242-S	30-0242	3.5 mm x 28 mm Locking Hexalobe Screw
30-0243-S	30-0243	3.5 mm x 30 mm Locking Hexalobe Screw
30-0244-S	30-0244	3.5 mm x 32 mm Locking Hexalobe Screw
30-0245-S	30-0245	3.5 mm x 34 mm Locking Hexalobe Screw
30-0246-S	30-0246	3.5 mm x 36 mm Locking Hexalobe Screw
30-0247-S	30-0247	3.5 mm x 38 mm Locking Hexalobe Screw
30-0248-S	30-0248	3.5 mm x 40 mm Locking Hexalobe Screw
30-0249-S	30-0249	3.5 mm x 45 mm Locking Hexalobe Screw
30-0250-S	30-0250	3.5 mm x 50 mm Locking Hexalobe Screw
30-0251-S	30-0251	3.5 mm x 55 mm Locking Hexalobe Screw
30-0252-S	30-0252	3.5 mm x 60 mm Locking Hexalobe Screw
30-0255-S	30-0255	3.5 mm x 8 mm Nonlocking Hexalobe Screw
30-0256-S	30-0256	3.5 mm x 10 mm Nonlocking Hexalobe Screw
30-0257-S	30-0257	3.5 mm x 12 mm Nonlocking Hexalobe Screw
30-0258-S	30-0258	3.5 mm x 14 mm Nonlocking Hexalobe Screw
30-0259-S	30-0259	3.5 mm x 16 mm Nonlocking Hexalobe Screw
30-0260-S	30-0260	3.5 mm x 18 mm Nonlocking Hexalobe Screw
30-0261-S	30-0261	3.5 mm x 20 mm Nonlocking Hexalobe Screw
30-0262-S	30-0262	3.5 mm x 22 mm Nonlocking Hexalobe Screw
30-0263-S	30-0263	3.5 mm x 24 mm Nonlocking Hexalobe Screw
30-0264-S	30-0264	3.5 mm x 26 mm Nonlocking Hexalobe Screw
30-0265-S	30-0265	3.5 mm x 28 mm Nonlocking Hexalobe Screw
30-0266-S	30-0266	3.5 mm x 30 mm Nonlocking Hexalobe Screw
30-0267-S	30-0267	3.5 mm x 32 mm Nonlocking Hexalobe Screw
30-0268-S	30-0268	3.5 mm x 34 mm Nonlocking Hexalobe Screw
30-0269-S	30-0269	3.5 mm x 36 mm Nonlocking Hexalobe Screw
30-0270-S	30-0270	3.5 mm x 38 mm Nonlocking Hexalobe Screw

Foot and Ankle Portfolio [continued]

Calcaneal Plating System

Sterile Part Number	Non-Sterile Part Number	Description
30-0271-S	30-0271	3.5 mm x 40 mm Nonlocking Hexalobe Screw
30-0272-S	30-0272	3.5 mm x 45 mm Nonlocking Hexalobe Screw
30-0273-S	30-0273	3.5 mm x 50 mm Nonlocking Hexalobe Screw
30-0274-S	30-0274	3.5 mm x 55 mm Nonlocking Hexalobe Screw
30-0275-S	30-0275	3.5 mm x 60 mm Nonlocking Hexalobe Screw
30-0278-S	30-0278	3.0 mm x 8 mm Locking Hexalobe Screw
30-0279-S	30-0279	3.0 mm x 10 mm Locking Hexalobe Screw
30-0280-S	30-0280	3.0 mm x 12 mm Locking Hexalobe Screw
30-0281-S	30-0281	3.0 mm x 14 mm Locking Hexalobe Screw
30-0282-S	30-0282	3.0 mm x 16 mm Locking Hexalobe Screw
30-0283-S	30-0283	3.0 mm x 18 mm Locking Hexalobe Screw
30-0284-S	30-0284	3.0 mm x 20 mm Locking Hexalobe Screw
30-0285-S	30-0285	3.0 mm x 22 mm Locking Hexalobe Screw
30-0286-S	30-0286	3.0 mm x 24 mm Locking Hexalobe Screw
30-0287-S	30-0287	3.0 mm x 26 mm Locking Hexalobe Screw
30-0288-S	30-0288	3.0 mm x 28 mm Locking Hexalobe Screw
30-0289-S	30-0289	3.0 mm x 30 mm Locking Hexalobe Screw
30-0290-S	30-0290	3.0 mm x 32 mm Locking Hexalobe Screw
30-0291-S	30-0291	3.0 mm x 34 mm Locking Hexalobe Screw
30-0292-S	30-0292	3.0 mm x 36 mm Locking Hexalobe Screw
30-0293-S	30-0293	3.0 mm x 38 mm Locking Hexalobe Screw
30-0294-S	30-0294	3.0 mm x 40 mm Locking Hexalobe Screw
30-0295-S	30-0295	3.0 mm x 45 mm Locking Hexalobe Screw
30-0296-S	30-0296	3.0 mm x 50 mm Locking Hexalobe Screw
30-0297-S	30-0297	3.0 mm x 55 mm Locking Hexalobe Screw
30-0301-S	30-0301	3.0 mm x 8 mm Nonlocking Hexalobe Screw
30-0302-S	30-0302	3.0 mm x 10 mm Nonlocking Hexalobe Screw
30-0303-S	30-0303	3.0 mm x 12 mm Nonlocking Hexalobe Screw
30-0304-S	30-0304	3.0 mm x 14 mm Nonlocking Hexalobe Screw

Foot and Ankle Portfolio [continued]

Calcaneal Plating System

Sterile Part Number	Non-Sterile Part Number	Description
30-0305-S	30-0305	3.0 mm x 16 mm Nonlocking Hexalobe Screw
30-0306-S	30-0306	3.0 mm x 18 mm Nonlocking Hexalobe Screw
30-0307-S	30-0307	3.0 mm x 20 mm Nonlocking Hexalobe Screw
30-0308-S	30-0308	3.0 mm x 22 mm Nonlocking Hexalobe Screw
30-0309-S	30-0309	3.0 mm x 24 mm Nonlocking Hexalobe Screw
30-0310-S	30-0310	3.0 mm x 26 mm Nonlocking Hexalobe Screw
30-0311-S	30-0311	3.0 mm x 28 mm Nonlocking Hexalobe Screw
30-0312-S	30-0312	3.0 mm x 30 mm Nonlocking Hexalobe Screw
30-0313-S	30-0313	3.0 mm x 32 mm Nonlocking Hexalobe Screw
30-0314-S	30-0314	3.0 mm x 34 mm Nonlocking Hexalobe Screw
30-0315-S	30-0315	3.0 mm x 36 mm Nonlocking Hexalobe Screw
30-0316-S	30-0316	3.0 mm x 38 mm Nonlocking Hexalobe Screw
30-0317-S	30-0317	3.0 mm x 40 mm Nonlocking Hexalobe Screw
30-0318-S	30-0318	3.0 mm x 45 mm Nonlocking Hexalobe Screw
30-0319-S	30-0319	3.0 mm x 50 mm Nonlocking Hexalobe Screw
30-0320-S	30-0320	3.0 mm x 55 mm Nonlocking Hexalobe Screw
N/A	35-0002	4.0 mm Steinmann Pin, Smooth
N/A	35-0003	5.0 mm Steinmann Pin, Smooth
70-0022-S	70-0022	Locking Calcaneal Plate, Small Left
70-0023-S	70-0023	Locking Calcaneal Plate, Small Right
70-0024-S	70-0024	Locking Calcaneal Plate, Medium Left
70-0025-S	70-0025	Locking Calcaneal Plate, Medium Right
70-0026-S	70-0026	Locking Calcaneal Plate, Large Left
70-0027-S	70-0027	Locking Calcaneal Plate, Large Right
70-0386-S	70-0386	Anterior Process Calcaneal Plate, Medium, Left
70-0387-S	70-0387	Anterior Process Calcaneal Plate, Medium, Right
70-0388-S	70-0388	Anterior Process Calcaneal Plate, Large, Left
70-0389-S	70-0389	Anterior Process Calcaneal Plate, Large, Right
70-0394-S	70-0394	Posterior Tuberosity Calcaneal Plate, 5-Hole, Left

Foot and Ankle Portfolio [continued]

Calcaneal Plating System

Sterile Part Number	Non-Sterile Part Number	Description
70-0395-S	70-0395	Posterior Tuberosity Calcaneal Plate, 5-Hole, Right
70-0396-S	70-0396	Posterior Tuberosity Calcaneal Plate, 6-Hole, Left
70-0397-S	70-0397	Posterior Tuberosity Calcaneal Plate, 6-Hole, Right
70-0400-S	70-0400	Combo Calcaneal Plate, Medium, 8-Hole, Left
70-0401-S	70-0401	Combo Calcaneal Plate, Medium, 8-Hole, Right
70-0404-S	70-0404	Combo Calcaneal Plate, Large, 8-Hole, Left
70-0405-S	70-0405	Combo Calcaneal Plate, Large, 8-Hole, Right
70-0406-S	70-0406	Combo Calcaneal Plate, Large, 9-Hole, Left
70-0407-S	70-0407	Combo Calcaneal Plate, Large, 9-Hole, Right
CA-4120-S	CA-4120	4.0 mm x 12.0 mm Cancellous Screw
CA-4140-S	CA-4140	4.0 mm x 14.0 mm Cancellous Screw
CA-4160-S	CA-4160	4.0 mm x 16.0 mm Cancellous Screw
CA-4180-S	CA-4180	4.0 mm x 18.0 mm Cancellous Screw
CA-4200-S	CA-4200	4.0 mm x 20.0 mm Cancellous Screw
CA-4220-S	CA-4220	4.0 mm x 22.0 mm Cancellous Screw
CA-4240-S	CA-4240	4.0 mm x 24.0 mm Cancellous Screw
CA-4260-S	CA-4260	4.0 mm x 26.0 mm Cancellous Screw
CA-4280-S	CA-4280	4.0 mm x 28.0 mm Cancellous Screw
CA-4300-S	CA-4300	4.0 mm x 30.0 mm Cancellous Screw
CA-4350-S	CA-4350	4.0 mm x 35.0 mm Cancellous Screw
CA-4400-S	CA-4400	4.0 mm x 40.0 mm Cancellous Screw
CA-4450-S	CA-4450	4.0 mm x 45.0 mm Cancellous Screw
CA-4500-S	CA-4500	4.0 mm x 50.0 mm Cancellous Screw
CA-4550-S	CA-4550	4.0 mm x 55.0 mm Cancellous Screw
CA-4600-S	CA-4600	4.0 mm x 60.0 mm Cancellous Screw
CO-2708-S	CO-2708	2.7 mm x 8.0 mm Cortical Screw
CO-2710-S	CO-2710	2.7 mm x 10.0 mm Cortical Screw
CO-2712-S	CO-2712	2.7 mm x 12.0 mm Cortical Screw
CO-2714-S	CO-2714	2.7 mm x 14.0 mm Cortical Screw

Foot and Ankle Portfolio [continued]

Calcaneal Plating System

Sterile Part Number	Non-Sterile Part Number	Description
CO-2716-S	CO-2716	2.7 mm x 16.0 mm Cortical Screw
CO-2718-S	CO-2718	2.7 mm x 18.0 mm Cortical Screw
CO-2720-S	CO-2720	2.7 mm x 20.0 mm Cortical Screw
CO-2722-S	CO-2722	2.7 mm x 22.0 mm Cortical Screw
CO-2724-S	CO-2724	2.7 mm x 24.0 mm Cortical Screw
CO-2726-S	CO-2726	2.7 mm x 26.0 mm Cortical Screw
CO-2728-S	CO-2728	2.7 mm x 28.0 mm Cortical Screw
CO-2730-S	CO-2730	2.7 mm x 30.0 mm Cortical Screw
CO-2732-S	CO-2732	2.7 mm x 32.0 mm Cortical Screw
CO-2734-S	CO-2734	2.7 mm x 34.0 mm Cortical Screw
CO-2736-S	CO-2736	2.7 mm x 36.0 mm Cortical Screw
CO-2738-S	CO-2738	2.7 mm x 38.0 mm Cortical Screw
CO-2740-S	CO-2740	2.7 mm x 40.0 mm Cortical Screw
CO-2745-S	CO-2745	2.7 mm x 45.0 mm Cortical Screw
CO-2750-S	CO-2750	2.7 mm x 50.0 mm Cortical Screw
CO-2755-S	CO-2755	2.7 mm x 55.0 mm Cortical Screw
CO-3080-S	CO-3080	3.5 mm x 8.0 mm Cortical Screw
CO-3100-S	CO-3100	3.5 mm x 10.0 mm Cortical Screw
CO-3120-S	CO-3120	3.5 mm x 12.0 mm Cortical Screw
CO-3140-S	CO-3140	3.5 mm x 14.0 mm Cortical Screw
CO-3160-S	CO-3160	3.5 mm x 16.0 mm Cortical Screw
CO-3180-S	CO-3180	3.5 mm x 18.0 mm Cortical Screw
CO-3200-S	CO-3200	3.5 mm x 20.0 mm Cortical Screw
CO-3220-S	CO-3220	3.5 mm x 22.0 mm Cortical Screw
CO-3240-S	CO-3240	3.5 mm x 24.0 mm Cortical Screw
CO-3260-S	CO-3260	3.5 mm x 26.0 mm Cortical Screw
CO-3280-S	CO-3280	3.5 mm x 28.0 mm Cortical Screw
CO-3300-S	CO-3300	3.5 mm x 30.0 mm Cortical Screw
CO-3320-S	CO-3320	3.5 mm x 32.0 mm Cortical Screw

Foot and Ankle Portfolio [continued]

Calcaneal Plating System

Sterile Part Number	Non-Sterile Part Number	Description
CO-3340-S	CO-3340	3.5 mm x 34.0 mm Cortical Screw
CO-3360-S	CO-3360	3.5 mm x 36.0 mm Cortical Screw
CO-3380-S	CO-3380	3.5 mm x 38.0 mm Cortical Screw
CO-3400-S	CO-3400	3.5 mm x 40.0 mm Cortical Screw
CO-3450-S	CO-3450	3.5 mm x 45.0 mm Cortical Screw
CO-3500-S	CO-3500	3.5 mm x 50.0 mm Cortical Screw
CO-3550-S	CO-3550	3.5 mm x 55.0 mm Cortical Screw
CO-3600-S	CO-3600	3.5 mm x 60.0 mm Cortical Screw
COL-2080-S	COL-2080	2.7 mm x 8.0 mm Locking Cortical Screw
COL-2100-S	COL-2100	2.7 mm x 10.0 mm Locking Cortical Screw
COL-2120-S	COL-2120	2.7 mm x 12.0 mm Locking Cortical Screw
COL-2140-S	COL-2140	2.7 mm x 14.0 mm Locking Cortical Screw
COL-2160-S	COL-2160	2.7 mm x 16.0 mm Locking Cortical Screw
COL-2180-S	COL-2180	2.7 mm x 18.0 mm Locking Cortical Screw
COL-2200-S	COL-2200	2.7 mm x 20.0 mm Locking Cortical Screw
COL-2220-S	COL-2220	2.7 mm x 22.0 mm Locking Cortical Screw
COL-2240-S	COL-2240	2.7 mm x 24.0 mm Locking Cortical Screw
COL-2260-S	COL-2260	2.7 mm x 26.0 mm Locking Cortical Screw
COL-2280-S	COL-2280	2.7 mm x 28.0 mm Locking Cortical Screw
COL-2300-S	COL-2300	2.7 mm x 30.0 mm Locking Cortical Screw
COL-2320-S	COL-2320	2.7 mm x 32.0 mm Locking Cortical Screw
COL-2340-S	COL-2340	2.7 mm x 34.0 mm Locking Cortical Screw
COL-2360-S	COL-2360	2.7 mm x 36.0 mm Locking Cortical Screw
COL-2380-S	COL-2380	2.7 mm x 38.0 mm Locking Cortical Screw
COL-2400-S	COL-2400	2.7 mm x 40.0 mm Locking Cortical Screw
COL-2450-S	COL-2450	2.7 mm x 45.0 mm Locking Cortical Screw
COL-2500-S	COL-2500	2.7 mm x 50.0 mm Locking Cortical Screw
COL-2550-S	COL-2550	2.7 mm x 55.0 mm Locking Cortical Screw
COL-3080-S	COL-3080	3.5 mm x 8.0 mm Locking Cortical Screw

Foot and Ankle Portfolio [continued]

Calcaneal Plating System

Sterile Part Number	Non-Sterile Part Number	Description
COL-3100-S	COL-3100	3.5 mm x 10.0 mm Locking Cortical Screw
COL-3120-S	COL-3120	3.5 mm x 12.0 mm Locking Cortical Screw
COL-3140-S	COL-3140	3.5 mm x 14.0 mm Locking Cortical Screw
COL-3160-S	COL-3160	3.5 mm x 16.0 mm Locking Cortical Screw
COL-3180-S	COL-3180	3.5 mm x 18.0 mm Locking Cortical Screw
COL-3200-S	COL-3200	3.5 mm x 20.0 mm Locking Cortical Screw
COL-3220-S	COL-3220	3.5 mm x 22.0 mm Locking Cortical Screw
COL-3240-S	COL-3240	3.5 mm x 24.0 mm Locking Cortical Screw
COL-3260-S	COL-3260	3.5 mm x 26.0 mm Locking Cortical Screw
COL-3280-S	COL-3280	3.5 mm x 28.0 mm Locking Cortical Screw
COL-3300-S	COL-3300	3.5 mm x 30.0 mm Locking Cortical Screw
COL-3320-S	COL-3320	3.5 mm x 32.0 mm Locking Cortical Screw
COL-3340-S	COL-3340	3.5 mm x 34.0 mm Locking Cortical Screw
COL-3360-S	COL-3360	3.5 mm x 36.0 mm Locking Cortical Screw
COL-3380-S	COL-3380	3.5 mm x 38.0 mm Locking Cortical Screw
COL-3400-S	COL-3400	3.5 mm x 40.0 mm Locking Cortical Screw
COL-3450-S	COL-3450	3.5 mm x 45.0 mm Locking Cortical Screw
COL-3500-S	COL-3500	3.5 mm x 50.0 mm Locking Cortical Screw
COL-3550-S	COL-3550	3.5 mm x 55.0 mm Locking Cortical Screw
COL-3600-S	COL-3600	3.5 mm x 60.0 mm Locking Cortical Screw

Forefoot-Midfoot Plating System

Sterile Part Number	Non-Sterile Part Number	Description
N/A	30-0232	3.5 mm x 8 mm Locking Hexalobe Screw
N/A	30-0233	3.5 mm x 10 mm Locking Hexalobe Screw
N/A	30-0234	3.5 mm x 12 mm Locking Hexalobe Screw
N/A	30-0235	3.5 mm x 14 mm Locking Hexalobe Screw
N/A	30-0236	3.5 mm x 16 mm Locking Hexalobe Screw
N/A	30-0237	3.5 mm x 18 mm Locking Hexalobe Screw

Foot and Ankle Portfolio [continued]

Forefoot-Midfoot Plating System

Sterile Part Number	Non-Sterile Part Number	Description
N/A	30-0238	3.5 mm x 20 mm Locking Hexalobe Screw
N/A	30-0239	3.5 mm x 22 mm Locking Hexalobe Screw
N/A	30-0240	3.5 mm x 24 mm Locking Hexalobe Screw
N/A	30-0241	3.5 mm x 26 mm Locking Hexalobe Screw
N/A	30-0242	3.5 mm x 28 mm Locking Hexalobe Screw
N/A	30-0243	3.5 mm x 30 mm Locking Hexalobe Screw
N/A	30-0244	3.5 mm x 32 mm Locking Hexalobe Screw
N/A	30-0245	3.5 mm x 34 mm Locking Hexalobe Screw
N/A	30-0246	3.5 mm x 36 mm Locking Hexalobe Screw
N/A	30-0247	3.5 mm x 38 mm Locking Hexalobe Screw
N/A	30-0248	3.5 mm x 40 mm Locking Hexalobe Screw
N/A	30-0249	3.5 mm x 45 mm Locking Hexalobe Screw
N/A	30-0250	3.5 mm x 50 mm Locking Hexalobe Screw
N/A	30-0251	3.5 mm x 55 mm Locking Hexalobe Screw
N/A	30-0252	3.5 mm x 60 mm Locking Hexalobe Screw
30-0255-S	30-0255	3.5 mm x 8 mm Nonlocking Hexalobe Screw
30-0256-S	30-0256	3.5 mm x 10 mm Nonlocking Hexalobe Screw
30-0257-S	30-0257	3.5 mm x 12 mm Nonlocking Hexalobe Screw
30-0258-S	30-0258	3.5 mm x 14 mm Nonlocking Hexalobe Screw
30-0259-S	30-0259	3.5 mm x 16 mm Nonlocking Hexalobe Screw
30-0260-S	30-0260	3.5 mm x 18 mm Nonlocking Hexalobe Screw
30-0261-S	30-0261	3.5 mm x 20 mm Nonlocking Hexalobe Screw
30-0262-S	30-0262	3.5 mm x 22 mm Nonlocking Hexalobe Screw
30-0263-S	30-0263	3.5 mm x 24 mm Nonlocking Hexalobe Screw
30-0264-S	30-0264	3.5 mm x 26 mm Nonlocking Hexalobe Screw
30-0265-S	30-0265	3.5 mm x 28 mm Nonlocking Hexalobe Screw
30-0266-S	30-0266	3.5 mm x 30 mm Nonlocking Hexalobe Screw
30-0267-S	30-0267	3.5 mm x 32 mm Nonlocking Hexalobe Screw
30-0268-S	30-0268	3.5 mm x 34 mm Nonlocking Hexalobe Screw

Foot and Ankle Portfolio [continued]

Forefoot-Midfoot Plating System

Sterile Part Number	Non-Sterile Part Number	Description
30-0269-S	30-0269	3.5 mm x 36 mm Nonlocking Hexalobe Screw
30-0270-S	30-0270	3.5 mm x 38 mm Nonlocking Hexalobe Screw
30-0271-S	30-0271	3.5 mm x 40 mm Nonlocking Hexalobe Screw
30-0272-S	30-0272	3.5 mm x 45 mm Nonlocking Hexalobe Screw
30-0273-S	30-0273	3.5 mm x 50 mm Nonlocking Hexalobe Screw
30-0274-S	30-0274	3.5 mm x 55 mm Nonlocking Hexalobe Screw
30-0275-S	30-0275	3.5 mm x 60 mm Nonlocking Hexalobe Screw
30-0278-S	30-0278	3.0 mm x 8 mm Locking Hexalobe Screw
30-0279-S	30-0279	3.0 mm x 10 mm Locking Hexalobe Screw
30-0280-S	30-0280	3.0 mm x 12 mm Locking Hexalobe Screw
30-0281-S	30-0281	3.0 mm x 14 mm Locking Hexalobe Screw
30-0282-S	30-0282	3.0 mm x 16 mm Locking Hexalobe Screw
30-0283-S	30-0283	3.0 mm x 18 mm Locking Hexalobe Screw
30-0284-S	30-0284	3.0 mm x 20 mm Locking Hexalobe Screw
30-0285-S	30-0285	3.0 mm x 22 mm Locking Hexalobe Screw
30-0286-S	30-0286	3.0 mm x 24 mm Locking Hexalobe Screw
30-0287-S	30-0287	3.0 mm x 26 mm Locking Hexalobe Screw
30-0288-S	30-0288	3.0 mm x 28 mm Locking Hexalobe Screw
30-0289-S	30-0289	3.0 mm x 30 mm Locking Hexalobe Screw
30-0290-S	30-0290	3.0 mm x 32 mm Locking Hexalobe Screw
30-0291-S	30-0291	3.0 mm x 34 mm Locking Hexalobe Screw
30-0292-S	30-0292	3.0 mm x 36 mm Locking Hexalobe Screw
30-0293-S	30-0293	3.0 mm x 38 mm Locking Hexalobe Screw
30-0294-S	30-0294	3.0 mm x 40 mm Locking Hexalobe Screw
30-0295-S	30-0295	3.0 mm x 45 mm Locking Hexalobe Screw
30-0296-S	30-0296	3.0 mm x 50 mm Locking Hexalobe Screw
30-0297-S	30-0297	3.0 mm x 55 mm Locking Hexalobe Screw
30-0301-S	30-0301	3.0 mm x 8 mm Nonlocking Hexalobe Screw
30-0302-S	30-0302	3.0 mm x 10 mm Nonlocking Hexalobe Screw

Foot and Ankle Portfolio [continued]

Forefoot-Midfoot Plating System

Sterile Part Number	Non-Sterile Part Number	Description
30-0303-S	30-0303	3.0 mm x 12 mm Nonlocking Hexalobe Screw
30-0304-S	30-0304	3.0 mm x 14 mm Nonlocking Hexalobe Screw
30-0305-S	30-0305	3.0 mm x 16 mm Nonlocking Hexalobe Screw
30-0306-S	30-0306	3.0 mm x 18 mm Nonlocking Hexalobe Screw
30-0307-S	30-0307	3.0 mm x 20 mm Nonlocking Hexalobe Screw
30-0308-S	30-0308	3.0 mm x 22 mm Nonlocking Hexalobe Screw
30-0309-S	30-0309	3.0 mm x 24 mm Nonlocking Hexalobe Screw
30-0310-S	30-0310	3.0 mm x 26 mm Nonlocking Hexalobe Screw
30-0311-S	30-0311	3.0 mm x 28 mm Nonlocking Hexalobe Screw
30-0312-S	30-0312	3.0 mm x 30 mm Nonlocking Hexalobe Screw
30-0313-S	30-0313	3.0 mm x 32 mm Nonlocking Hexalobe Screw
30-0314-S	30-0314	3.0 mm x 34 mm Nonlocking Hexalobe Screw
30-0315-S	30-0315	3.0 mm x 36 mm Nonlocking Hexalobe Screw
30-0316-S	30-0316	3.0 mm x 38 mm Nonlocking Hexalobe Screw
30-0317-S	30-0317	3.0 mm x 40 mm Nonlocking Hexalobe Screw
30-0318-S	30-0318	3.0 mm x 45 mm Nonlocking Hexalobe Screw
30-0319-S	30-0319	3.0 mm x 50 mm Nonlocking Hexalobe Screw
30-0320-S	30-0320	3.0 mm x 55 mm Nonlocking Hexalobe Screw
70-0007-S	70-0007	5-Hole Locking 1st Tarsometatarsal Plate
70-0008-S	70-0008	4-Hole Locking 1st Tarsometatarsal Plate
70-0009-S	70-0009	4-Hole Locking 2nd & 3rd Tarsometatarsal Plate
70-0010-S	70-0010	Locking Proximal Metatarsophalangeal Wedge Plate, Left
70-0011-S	70-0011	Locking Proximal Metatarsophalangeal Wedge Plate, Right
70-0012-S	70-0012	Locking Dorsal Metatarsophalangeal Fusion Plate, Left
70-0013-S	70-0013	Locking Dorsal Metatarsophalangeal Fusion Plate, Right
70-0018-S	70-0018	Locking Metatarsophalangeal Revision Plate, Left
70-0019-S	70-0019	Locking Metatarsophalangeal Revision Plate, Right
70-0036-S	70-0036	Locking Metatarsophalangeal/Metatarsophalangeal Joint Combo Plate, Left

Foot and Ankle Portfolio [continued]

Forefoot-Midfoot Plating System

Sterile Part Number	Non-Sterile Part Number	Description
70-0037-S	70-0037	Locking Metatarsophalangeal/Metatarsophalangeal Joint Combo Plate, Right
70-0049-S	70-0049	8-Hole Locking 1st Tarsometatarsal Plate
70-0050-S	70-0050	7-Hole Locking 1st Tarsometatarsal Plate
70-0051-S	70-0051	7-Hole Locking 2nd & 3rd Tarsometatarsal Plate
70-0109-S	70-0109	Locking Metatarsophalangeal/Metatarsophalangeal Joint Combo Plate 4 Degree Left
70-0110-S	70-0110	Locking Metatarsophalangeal/Metatarsophalangeal Joint Combo Plate 4 Degree Right
70-0141-S	70-0141	Locking Dorsal Metatarsophalangeal Plate 4 Degree Left
70-0142-S	70-0142	Locking Dorsal Metatarsophalangeal Plate 4 Degree Right
70-0143-S	70-0143	Locking Metatarsophalangeal Revision Plate 4 Degree Left
70-0144-S	70-0144	Locking Metatarsophalangeal Revision Plate 4 Degree Right
70-0323-S	70-0323	Locking Proximal Metatarsal Wedge Plate, 0 Degree, Left
70-0324-S	70-0324	Locking Proximal Metatarsal Wedge Plate, 0 Degree, Right
70-0329-S	70-0329	Locking Metatarsophalangeal Joint Petite Plate 4 Degree, Left
70-0330-S	70-0330	Locking Metatarsophalangeal Joint Petite Plate 4 Degree, Right
CA-4120-S	CA-4120	4.0 mm x 12.0 mm Cancellous Screw
CA-4140-S	CA-4140	4.0 mm x 14.0 mm Cancellous Screw
CA-4160-S	CA-4160	4.0 mm x 16.0 mm Cancellous Screw
CA-4180-S	CA-4180	4.0 mm x 18.0 mm Cancellous Screw
CA-4200-S	CA-4200	4.0 mm x 20.0 mm Cancellous Screw
CA-4220-S	CA-4220	4.0 mm x 22.0 mm Cancellous Screw
CA-4240-S	CA-4240	4.0 mm x 24.0 mm Cancellous Screw
CA-4260-S	CA-4260	4.0 mm x 26.0 mm Cancellous Screw
CA-4280-S	CA-4280	4.0 mm x 28.0 mm Cancellous Screw
CA-4300-S	CA-4300	4.0 mm x 30.0 mm Cancellous Screw
CA-4350-S	CA-4350	4.0 mm x 35.0 mm Cancellous Screw
CA-4400-S	CA-4400	4.0 mm x 40.0 mm Cancellous Screw
CA-4450-S	CA-4450	4.0 mm x 45.0 mm Cancellous Screw
CA-4500-S	CA-4500	4.0 mm x 50.0 mm Cancellous Screw

Foot and Ankle Portfolio [continued]

Forefoot-Midfoot Plating System

Sterile Part Number	Non-Sterile Part Number	Description
CA-4550-S	CA-4550	4.0 mm x 55.0 mm Cancellous Screw
CA-4600-S	CA-4600	4.0 mm x 60.0 mm Cancellous Screw
CO-2708-S	CO-2708	2.7 mm x 8.0 mm Cortical Screw
CO-2710-S	CO-2710	2.7 mm x 10.0 mm Cortical Screw
CO-2712-S	CO-2712	2.7 mm x 12.0 mm Cortical Screw
CO-2714-S	CO-2714	2.7 mm x 14.0 mm Cortical Screw
CO-2716-S	CO-2716	2.7 mm x 16.0 mm Cortical Screw
CO-2718-S	CO-2718	2.7 mm x 18.0 mm Cortical Screw
CO-2720-S	CO-2720	2.7 mm x 20.0 mm Cortical Screw
CO-2722-S	CO-2722	2.7 mm x 22.0 mm Cortical Screw
CO-2724-S	CO-2724	2.7 mm x 24.0 mm Cortical Screw
CO-2726-S	CO-2726	2.7 mm x 26.0 mm Cortical Screw
CO-2728-S	CO-2728	2.7 mm x 28.0 mm Cortical Screw
CO-2730-S	CO-2730	2.7 mm x 30.0 mm Cortical Screw
CO-2732-S	CO-2732	2.7 mm x 32.0 mm Cortical Screw
CO-2734-S	CO-2734	2.7 mm x 34.0 mm Cortical Screw
CO-2736-S	CO-2736	2.7 mm x 36.0 mm Cortical Screw
CO-2738-S	CO-2738	2.7 mm x 38.0 mm Cortical Screw
CO-2740-S	CO-2740	2.7 mm x 40.0 mm Cortical Screw
CO-2745-S	CO-2745	2.7 mm x 45.0 mm Cortical Screw
CO-2750-S	CO-2750	2.7 mm x 50.0 mm Cortical Screw
CO-2755-S	CO-2755	2.7 mm x 55.0 mm Cortical Screw
CO-3080-S	CO-3080	3.5 mm x 8.0 mm Cortical Screw
CO-3100-S	CO-3100	3.5 mm x 10.0 mm Cortical Screw
CO-3120-S	CO-3120	3.5 mm x 12.0 mm Cortical Screw
CO-3140-S	CO-3140	3.5 mm x 14.0 mm Cortical Screw
CO-3160-S	CO-3160	3.5 mm x 16.0 mm Cortical Screw
CO-3180-S	CO-3180	3.5 mm x 18.0 mm Cortical Screw
CO-3200-S	CO-3200	3.5 mm x 20.0 mm Cortical Screw

Foot and Ankle Portfolio [continued]

Forefoot-Midfoot Plating System

Sterile Part Number	Non-Sterile Part Number	Description
CO-3220-S	CO-3220	3.5 mm x 22.0 mm Cortical Screw
CO-3240-S	CO-3240	3.5 mm x 24.0 mm Cortical Screw
CO-3260-S	CO-3260	3.5 mm x 26.0 mm Cortical Screw
CO-3280-S	CO-3280	3.5 mm x 28.0 mm Cortical Screw
CO-3300-S	CO-3300	3.5 mm x 30.0 mm Cortical Screw
CO-3320-S	CO-3320	3.5 mm x 32.0 mm Cortical Screw
CO-3340-S	CO-3340	3.5 mm x 34.0 mm Cortical Screw
CO-3360-S	CO-3360	3.5 mm x 36.0 mm Cortical Screw
CO-3380-S	CO-3380	3.5 mm x 38.0 mm Cortical Screw
CO-3400-S	CO-3400	3.5 mm x 40.0 mm Cortical Screw
CO-3450-S	CO-3450	3.5 mm x 45.0 mm Cortical Screw
CO-3500-S	CO-3500	3.5 mm x 50.0 mm Cortical Screw
CO-3550-S	CO-3550	3.5 mm x 55.0 mm Cortical Screw
CO-3600-S	CO-3600	3.5 mm x 60.0 mm Cortical Screw
COL-2080-S	COL-2080	2.7 mm x 8.0 mm Locking Cortical Screw
COL-2100-S	COL-2100	2.7 mm x 10.0 mm Locking Cortical Screw
COL-2120-S	COL-2120	2.7 mm x 12.0 mm Locking Cortical Screw
COL-2140-S	COL-2140	2.7 mm x 14.0 mm Locking Cortical Screw
COL-2160-S	COL-2160	2.7 mm x 16.0 mm Locking Cortical Screw
COL-2180-S	COL-2180	2.7 mm x 18.0 mm Locking Cortical Screw
COL-2200-S	COL-2200	2.7 mm x 20.0 mm Locking Cortical Screw
COL-2220-S	COL-2220	2.7 mm x 22.0 mm Locking Cortical Screw
COL-2240-S	COL-2240	2.7 mm x 24.0 mm Locking Cortical Screw
COL-2260-S	COL-2260	2.7 mm x 26.0 mm Locking Cortical Screw
COL-2280-S	COL-2280	2.7 mm x 28.0 mm Locking Cortical Screw
COL-2300-S	COL-2300	2.7 mm x 30.0 mm Locking Cortical Screw
COL-2320-S	COL-2320	2.7 mm x 32.0 mm Locking Cortical Screw
COL-2340-S	COL-2340	2.7 mm x 34.0 mm Locking Cortical Screw
COL-2360-S	COL-2360	2.7 mm x 36.0 mm Locking Cortical Screw

Foot and Ankle Portfolio [continued]

Forefoot-Midfoot Plating System

Sterile Part Number	Non-Sterile Part Number	Description
COL-2380-S	COL-2380	2.7 mm x 38.0 mm Locking Cortical Screw
COL-2400-S	COL-2400	2.7 mm x 40.0 mm Locking Cortical Screw
COL-2450-S	COL-2450	2.7 mm x 45.0 mm Locking Cortical Screw
COL-2500-S	COL-2500	2.7 mm x 50.0 mm Locking Cortical Screw
COL-2550-S	COL-2550	2.7 mm x 55.0 mm Locking Cortical Screw
COL-3080-S	COL-3080	3.5 mm x 8.0 mm Locking Cortical Screw
COL-3100-S	COL-3100	3.5 mm x 10.0 mm Locking Cortical Screw
COL-3120-S	COL-3120	3.5 mm x 12.0 mm Locking Cortical Screw
COL-3140-S	COL-3140	3.5 mm x 14.0 mm Locking Cortical Screw
COL-3160-S	COL-3160	3.5 mm x 16.0 mm Locking Cortical Screw
COL-3180-S	COL-3180	3.5 mm x 18.0 mm Locking Cortical Screw
COL-3200-S	COL-3200	3.5 mm x 20.0 mm Locking Cortical Screw
COL-3220-S	COL-3220	3.5 mm x 22.0 mm Locking Cortical Screw
COL-3240-S	COL-3240	3.5 mm x 24.0 mm Locking Cortical Screw
COL-3260-S	COL-3260	3.5 mm x 26.0 mm Locking Cortical Screw
COL-3280-S	COL-3280	3.5 mm x 28.0 mm Locking Cortical Screw
COL-3300-S	COL-3300	3.5 mm x 30.0 mm Locking Cortical Screw
COL-3320-S	COL-3320	3.5 mm x 32.0 mm Locking Cortical Screw
COL-3340-S	COL-3340	3.5 mm x 34.0 mm Locking Cortical Screw
COL-3360-S	COL-3360	3.5 mm x 36.0 mm Locking Cortical Screw
COL-3380-S	COL-3380	3.5 mm x 38.0 mm Locking Cortical Screw
COL-3400-S	COL-3400	3.5 mm x 40.0 mm Locking Cortical Screw
COL-3450-S	COL-3450	3.5 mm x 45.0 mm Locking Cortical Screw
COL-3500-S	COL-3500	3.5 mm x 50.0 mm Locking Cortical Screw
COL-3550-S	COL-3550	3.5 mm x 55.0 mm Locking Cortical Screw
COL-3600-S	COL-3600	3.5 mm x 60.0 mm Locking Cortical Screw

Foot and Ankle Portfolio [continued]

Small Fragment Base Set

Sterile Part Number	Non-Sterile Part Number	Description
N/A	30-0098	70.0 mm Tension Band Pin
N/A	30-0099	90.0 mm Tension Band Pin
30-0232-S	30-0232	3.5 mm x 8 mm Locking Hexalobe Screw
30-0233-S	30-0233	3.5 mm x 10 mm Locking Hexalobe Screw
30-0234-S	30-0234	3.5 mm x 12 mm Locking Hexalobe Screw
30-0235-S	30-0235	3.5 mm x 14 mm Locking Hexalobe Screw
30-0236-S	30-0236	3.5 mm x 16 mm Locking Hexalobe Screw
30-0237-S	30-0237	3.5 mm x 18 mm Locking Hexalobe Screw
30-0238-S	30-0238	3.5 mm x 20 mm Locking Hexalobe Screw
30-0239-S	30-0239	3.5 mm x 22 mm Locking Hexalobe Screw
30-0240-S	30-0240	3.5 mm x 24 mm Locking Hexalobe Screw
30-0241-S	30-0241	3.5 mm x 26 mm Locking Hexalobe Screw
30-0242-S	30-0242	3.5 mm x 28 mm Locking Hexalobe Screw
30-0243-S	30-0243	3.5 mm x 30 mm Locking Hexalobe Screw
30-0244-S	30-0244	3.5 mm x 32 mm Locking Hexalobe Screw
30-0245-S	30-0245	3.5 mm x 34 mm Locking Hexalobe Screw
30-0246-S	30-0246	3.5 mm x 36 mm Locking Hexalobe Screw
30-0247-S	30-0247	3.5 mm x 38 mm Locking Hexalobe Screw
30-0248-S	30-0248	3.5 mm x 40 mm Locking Hexalobe Screw
30-0250-S	30-0250	3.5 mm x 50 mm Locking Hexalobe Screw
30-0251-S	30-0251	3.5 mm x 55 mm Locking Hexalobe Screw
30-0252-S	30-0252	3.5 mm x 60 mm Locking Hexalobe Screw
30-0253-S	30-0253	3.5 mm x 65 mm Locking Hexalobe Screw
30-0255-S	30-0255	3.5 mm x 8 mm Nonlocking Hexalobe Screw
30-0256-S	30-0256	3.5 mm x 10 mm Nonlocking Hexalobe Screw
30-0257-S	30-0257	3.5 mm x 12 mm Nonlocking Hexalobe Screw
30-0258-S	30-0258	3.5 mm x 14 mm Nonlocking Hexalobe Screw
30-0259-S	30-0259	3.5 mm x 16 mm Nonlocking Hexalobe Screw
30-0260-S	30-0260	3.5 mm x 18 mm Nonlocking Hexalobe Screw

Foot and Ankle Portfolio [continued]

Small Fragment Base Set

Sterile Part Number	Non-Sterile Part Number	Description
30-0261-S	30-0261	3.5 mm x 20 mm Nonlocking Hexalobe Screw
30-0262-S	30-0262	3.5 mm x 22 mm Nonlocking Hexalobe Screw
30-0263-S	30-0263	3.5 mm x 24 mm Nonlocking Hexalobe Screw
30-0264-S	30-0264	3.5 mm x 26 mm Nonlocking Hexalobe Screw
30-0265-S	30-0265	3.5 mm x 28 mm Nonlocking Hexalobe Screw
30-0266-S	30-0266	3.5 mm x 30 mm Nonlocking Hexalobe Screw
30-0267-S	30-0267	3.5 mm x 32 mm Nonlocking Hexalobe Screw
30-0268-S	30-0268	3.5 mm x 34 mm Nonlocking Hexalobe Screw
30-0269-S	30-0269	3.5 mm x 36 mm Nonlocking Hexalobe Screw
30-0270-S	30-0270	3.5 mm x 38 mm Nonlocking Hexalobe Screw
30-0271-S	30-0271	3.5 mm x 40 mm Nonlocking Hexalobe Screw
30-0273-S	30-0273	3.5 mm x 50 mm Nonlocking Hexalobe Screw
30-0274-S	30-0274	3.5 mm x 55 mm Nonlocking Hexalobe Screw
30-0275-S	30-0275	3.5 mm x 60 mm Nonlocking Hexalobe Screw
30-0276-S	30-0276	3.5 mm x 65 mm Nonlocking Hexalobe Screw
30-0324-S	30-0324	2.7 mm x 8 mm Locking Hexalobe Screw
30-0325-S	30-0325	2.7 mm x 10 mm Locking Hexalobe Screw
30-0326-S	30-0326	2.7 mm x 12 mm Locking Hexalobe Screw
30-0327-S	30-0327	2.7 mm x 14 mm Locking Hexalobe Screw
30-0328-S	30-0328	2.7 mm x 16 mm Locking Hexalobe Screw
30-0329-S	30-0329	2.7 mm x 18 mm Locking Hexalobe Screw
30-0330-S	30-0330	2.7 mm x 20 mm Locking Hexalobe Screw
30-0331-S	30-0331	2.7 mm x 22 mm Locking Hexalobe Screw
30-0332-S	30-0332	2.7 mm x 24 mm Locking Hexalobe Screw
30-0333-S	30-0333	2.7 mm x 26 mm Locking Hexalobe Screw
30-0334-S	30-0334	2.7 mm x 28 mm Locking Hexalobe Screw
30-0335-S	30-0335	2.7 mm x 30 mm Locking Hexalobe Screw
30-0336-S	30-0336	2.7 mm x 32 mm Locking Hexalobe Screw
30-0337-S	30-0337	2.7 mm x 34 mm Locking Hexalobe Screw

Foot and Ankle Portfolio [continued]

Small Fragment Base Set

Sterile Part Number	Non-Sterile Part Number	Description
30-0338-S	30-0338	2.7 mm x 36 mm Locking Hexalobe Screw
30-0339-S	30-0339	2.7 mm x 38 mm Locking Hexalobe Screw
30-0340-S	30-0340	2.7 mm x 40 mm Locking Hexalobe Screw
30-0343-S	30-0343	2.7 mm x 8 mm Nonlocking Hexalobe Screw
30-0344-S	30-0344	2.7 mm x 10 mm Nonlocking Hexalobe Screw
30-0345-S	30-0345	2.7 mm x 12 mm Nonlocking Hexalobe Screw
30-0346-S	30-0346	2.7 mm x 14 mm Nonlocking Hexalobe Screw
30-0347-S	30-0347	2.7 mm x 16 mm Nonlocking Hexalobe Screw
30-0348-S	30-0348	2.7 mm x 18 mm Nonlocking Hexalobe Screw
30-0349-S	30-0349	2.7 mm x 20 mm Nonlocking Hexalobe Screw
30-0350-S	30-0350	2.7 mm x 22 mm Nonlocking Hexalobe Screw
30-0351-S	30-0351	2.7 mm x 24 mm Nonlocking Hexalobe Screw
30-0352-S	30-0352	2.7 mm x 26 mm Nonlocking Hexalobe Screw
30-0353-S	30-0353	2.7 mm x 28 mm Nonlocking Hexalobe Screw
30-0354-S	30-0354	2.7 mm x 30 mm Nonlocking Hexalobe Screw
30-0355-S	30-0355	2.7 mm x 32 mm Nonlocking Hexalobe Screw
30-0356-S	30-0356	2.7 mm x 34 mm Nonlocking Hexalobe Screw
30-0357-S	30-0357	2.7 mm x 36 mm Nonlocking Hexalobe Screw
30-0358-S	30-0358	2.7 mm x 38 mm Nonlocking Hexalobe Screw
30-0359-S	30-0359	2.7 mm x 40 mm Nonlocking Hexalobe Screw
30-0361-S	30-0361	2.7 mm x 50 mm Nonlocking Hexalobe Screw
30-0362-S	30-0362	2.7 mm x 55 mm Nonlocking Hexalobe Screw
30-0363-S	30-0363	2.7 mm x 60 mm Nonlocking Hexalobe Screw
3013-27010-S	3013-27010	2.7 mm x 10 mm Hexalobe Variable Angle Screw
3013-27012-S	3013-27012	2.7 mm x 12 mm Hexalobe Variable Angle Screw
3013-27014-S	3013-27014	2.7 mm x 14 mm Hexalobe Variable Angle Screw
3013-27016-S	3013-27016	2.7 mm x 16 mm Hexalobe Variable Angle Screw
3013-27018-S	3013-27018	2.7 mm x 18 mm Hexalobe Variable Angle Screw
3013-27020-S	3013-27020	2.7 mm x 20 mm Hexalobe Variable Angle Screw

Foot and Ankle Portfolio [continued]

Small Fragment Base Set

Sterile Part Number	Non-Sterile Part Number	Description
3013-27022-S	3013-27022	2.7 mm x 22 mm Hexalobe Variable Angle Screw
3013-27024-S	3013-27024	2.7 mm x 24 mm Hexalobe Variable Angle Screw
3013-27026-S	3013-27026	2.7 mm x 26 mm Hexalobe Variable Angle Screw
3013-27028-S	3013-27028	2.7 mm x 28 mm Hexalobe Variable Angle Screw
3013-27030-S	3013-27030	2.7 mm x 30 mm Hexalobe Variable Angle Screw
3013-27032-S	3013-27032	2.7 mm x 32 mm Hexalobe Variable Angle Screw
3013-27034-S	3013-27034	2.7 mm x 34 mm Hexalobe Variable Angle Screw
3013-27036-S	3013-27036	2.7 mm x 36 mm Hexalobe Variable Angle Screw
3013-27038-S	3013-27038	2.7 mm x 38 mm Hexalobe Variable Angle Screw
3013-27040-S	3013-27040	2.7 mm x 40 mm Hexalobe Variable Angle Screw
3013-27042-S	3013-27042	2.7 mm x 42 mm Hexalobe Variable Angle Screw
3013-27044-S	3013-27044	2.7 mm x 44 mm Hexalobe Variable Angle Screw
3013-27046-S	3013-27046	2.7 mm x 46 mm Hexalobe Variable Angle Screw
3013-27048-S	3013-27048	2.7 mm x 48 mm Hexalobe Variable Angle Screw
3013-27050-S	3013-27050	2.7 mm x 50 mm Hexalobe Variable Angle Screw
3013-27055-S	3013-27055	2.7 mm x 55 mm Hexalobe Variable Angle Screw
3013-27060-S	3013-27060	2.7 mm x 60 mm Hexalobe Variable Angle Screw
3013-35010-S	3013-35010	3.5 mm x 10 mm Hexalobe Variable Angle Screw
3013-35012-S	3013-35012	3.5 mm x 12 mm Hexalobe Variable Angle Screw
3013-35014-S	3013-35014	3.5 mm x 14 mm Hexalobe Variable Angle Screw
3013-35016-S	3013-35016	3.5 mm x 16 mm Hexalobe Variable Angle Screw
3013-35018-S	3013-35018	3.5 mm x 18 mm Hexalobe Variable Angle Screw
3013-35020-S	3013-35020	3.5 mm x 20 mm Hexalobe Variable Angle Screw
3013-35022-S	3013-35022	3.5 mm x 22 mm Hexalobe Variable Angle Screw
3013-35024-S	3013-35024	3.5 mm x 24 mm Hexalobe Variable Angle Screw
3013-35026-S	3013-35026	3.5 mm x 26 mm Hexalobe Variable Angle Screw
3013-35028-S	3013-35028	3.5 mm x 28 mm Hexalobe Variable Angle Screw
3013-35030-S	3013-35030	3.5 mm x 30 mm Hexalobe Variable Angle Screw
3013-35032-S	3013-35032	3.5 mm x 32 mm Hexalobe Variable Angle Screw

Foot and Ankle Portfolio [continued]

Small Fragment Base Set

Sterile Part Number	Non-Sterile Part Number	Description
3013-35034-S	3013-35034	3.5 mm x 34 mm Hexalobe Variable Angle Screw
3013-35036-S	3013-35036	3.5 mm x 36 mm Hexalobe Variable Angle Screw
3013-35038-S	3013-35038	3.5 mm x 38 mm Hexalobe Variable Angle Screw
3013-35040-S	3013-35040	3.5 mm x 40 mm Hexalobe Variable Angle Screw
3013-35042-S	3013-35042	3.5 mm x 42 mm Hexalobe Variable Angle Screw
3013-35044-S	3013-35044	3.5 mm x 44 mm Hexalobe Variable Angle Screw
3013-35046-S	3013-35046	3.5 mm x 46 mm Hexalobe Variable Angle Screw
3013-35048-S	3013-35048	3.5 mm x 48 mm Hexalobe Variable Angle Screw
3013-35050-S	3013-35050	3.5 mm x 50 mm Hexalobe Variable Angle Screw
3013-35055-S	3013-35055	3.5 mm x 55 mm Hexalobe Variable Angle Screw
3013-35060-S	3013-35060	3.5 mm x 60 mm Hexalobe Variable Angle Screw
3013-35065-S	3013-35065	3.5 mm x 65 mm Hexalobe Variable Angle Screw
3015-40010-S	3015-40010	4.0 mm x 10 mm Cancellous Hexalobe Screw
3015-40012-S	3015-40012	4.0 mm x 12 mm Cancellous Hexalobe Screw
3015-40014-S	3015-40014	4.0 mm x 14 mm Cancellous Hexalobe Screw
3015-40016-S	3015-40016	4.0 mm x 16 mm Cancellous Hexalobe Screw
3015-40018-S	3015-40018	4.0 mm x 18 mm Cancellous Hexalobe Screw
3015-40020-S	3015-40020	4.0 mm x 20 mm Cancellous Hexalobe Screw
3015-40022-S	3015-40022	4.0 mm x 22 mm Cancellous Hexalobe Screw
3015-40024-S	3015-40024	4.0 mm x 24 mm Cancellous Hexalobe Screw
3015-40026-S	3015-40026	4.0 mm x 26 mm Cancellous Hexalobe Screw
3015-40028-S	3015-40028	4.0 mm x 28 mm Cancellous Hexalobe Screw
3015-40030-S	3015-40030	4.0 mm x 30 mm Cancellous Hexalobe Screw
3015-40035-S	3015-40035	4.0 mm x 35 mm Cancellous Hexalobe Screw
3015-40040-S	3015-40040	4.0 mm x 40 mm Cancellous Hexalobe Screw
3015-40045-S	3015-40045	4.0 mm x 45 mm Cancellous Hexalobe Screw
3015-40050-S	3015-40050	4.0 mm x 50 mm Cancellous Hexalobe Screw
3015-40055-S	3015-40055	4.0 mm x 55 mm Cancellous Hexalobe Screw
3015-40060-S	3015-40060	4.0 mm x 60 mm Cancellous Hexalobe Screw

Foot and Ankle Portfolio [continued]

Small Fragment Base Set

Sterile Part Number	Non-Sterile Part Number	Description
3016-40012-S	3016-40012	4.0 mm x 12 mm Partially Threaded Cancellous Hexalobe Screw
3016-40014-S	3016-40014	4.0 mm x 14 mm Partially Threaded Cancellous Hexalobe Screw
3016-40016-S	3016-40016	4.0 mm x 16 mm Partially Threaded Cancellous Hexalobe Screw
3016-40018-S	3016-40018	4.0 mm x 18 mm Partially Threaded Cancellous Hexalobe Screw
3016-40020-S	3016-40020	4.0 mm x 20 mm Partially Threaded Cancellous Hexalobe Screw
3016-40022-S	3016-40022	4.0 mm x 22 mm Partially Threaded Cancellous Hexalobe Screw
3016-40024-S	3016-40024	4.0 mm x 24 mm Partially Threaded Cancellous Hexalobe Screw
3016-40026-S	3016-40026	4.0 mm x 26 mm Partially Threaded Cancellous Hexalobe Screw
3016-40028-S	3016-40028	4.0 mm x 28 mm Partially Threaded Cancellous Hexalobe Screw
3016-40030-S	3016-40030	4.0 mm x 30 mm Partially Threaded Cancellous Hexalobe Screw
3016-40035-S	3016-40035	4.0 mm x 35 mm Partially Threaded Cancellous Hexalobe Screw
3016-40040-S	3016-40040	4.0 mm x 40 mm Partially Threaded Cancellous Hexalobe Screw
3016-40045-S	3016-40045	4.0 mm x 45 mm Partially Threaded Cancellous Hexalobe Screw
3016-40050-S	3016-40050	4.0 mm x 50 mm Partially Threaded Cancellous Hexalobe Screw
3016-40055-S	3016-40055	4.0 mm x 55 mm Partially Threaded Cancellous Hexalobe Screw
3016-40060-S	3016-40060	4.0 mm x 60 mm Partially Threaded Cancellous Hexalobe Screw
30-2742-S	30-2742	2.7 mm x 42 mm Nonlocking Hexalobe Screw
30-2744-S	30-2744	2.7 mm x 44 mm Nonlocking Hexalobe Screw
30-2746-S	30-2746	2.7 mm x 46 mm Nonlocking Hexalobe Screw
30-2748-S	30-2748	2.7 mm x 48 mm Nonlocking Hexalobe Screw
30-2842-S	30-2842	2.7 mm x 42 mm Locking Hexalobe Screw
30-2844-S	30-2844	2.7 mm x 44 mm Locking Hexalobe Screw
30-2846-S	30-2846	2.7 mm x 46 mm Locking Hexalobe Screw

Foot and Ankle Portfolio [continued]

Small Fragment Base Set

Sterile Part Number	Non-Sterile Part Number	Description
30-2848-S	30-2848	2.7 mm x 48 mm Locking Hexalobe Screw
30-2850-S	30-2850	2.7 mm x 50 mm Locking Hexalobe Screw
30-2855-S	30-2855	2.7 mm x 55 mm Locking Hexalobe Screw
30-2860-S	30-2860	2.7 mm x 60 mm Locking Hexalobe Screw
30-3542-S	30-3542	3.5 mm x 42 mm Nonlocking Hexalobe Screw
30-3544-S	30-3544	3.5 mm x 44 mm Nonlocking Hexalobe Screw
30-3546-S	30-3546	3.5 mm x 46 mm Nonlocking Hexalobe Screw
30-3548-S	30-3548	3.5 mm x 48 mm Nonlocking Hexalobe Screw
30-3642-S	30-3642	3.5 mm x 42 mm Locking Hexalobe Screw
30-3644-S	30-3644	3.5 mm x 44 mm Locking Hexalobe Screw
30-3646-S	30-3646	3.5 mm x 46 mm Locking Hexalobe Screw
30-3648-S	30-3648	3.5 mm x 48 mm Locking Hexalobe Screw
7003-07036-S	7003-07036	Cannulated Screw Washer 7.0 mm OD x 3.6 mm ID
7008-0103-S	7008-0103	One-Third Tubular Plate 3-Hole
7008-0104-S	7008-0104	One-Third Tubular Plate 4-Hole
7008-0105-S	7008-0105	One-Third Tubular Plate 5-Hole
7008-0106-S	7008-0106	One-Third Tubular Plate 6-Hole
7008-0107-S	7008-0107	One-Third Tubular Plate 7-Hole
7008-0108-S	7008-0108	One-Third Tubular Plate 8-Hole
7008-0110-S	7008-0110	One-Third Tubular Plate 10-Hole
7008-0112-S	7008-0112	One-Third Tubular Plate 12-Hole
7010-0106N-S	7010-0106N	Fragment Plate 2.7 mm
7010-0107L-S	7010-0107L	L Fragment Plate 2.7 mm Left
7010-0107R-S	7010-0107R	L Fragment Plate 2.7 mm Right
7010-0108N-S	7010-0108N	T Fragment Plate 2.7 mm
AI-0020-S	AI-0020	20 mm AcuTwist® Acutrak Compression Screw
AI-0026-S	AI-0026	26 mm AcuTwist Acutrak Compression Screw
AI-0030-S	AI-0030	30 mm AcuTwist Acutrak Compression Screw

Foot and Ankle Portfolio [continued]

Conditions for safe scanning in the MRI

A patient with the plate, screw, pin, peg, and washer implants from the Ankle Plating, Ankle Plating 3, Calcaneal Plating, Forefoot-Midfoot Plating, and Small Fragment Base Set Systems may be safely scanned under the following conditions. Failure to follow these conditions may result in injury to the patient.

Nominal values(s) of Static Magnetic Field	1.5 T or 3 T
Maximum Spatial Field Gradient	19 T/m (1900 gauss/cm)
RF Excitation	Circularly Polarized (CP)
RF Transmit Coil Type	Whole body transmit coil, Head RF transmit-receive coil
Maximum Whole Body SAR	2.0 W/kg or 1.0 W/kg*
Limits on Scan Duration	2.0 W/kg whole body average SAR for 5 minutes of continuous RF (a sequence or back to back series/scan without breaks) / 1.0 W/kg whole body average SAR for 15 minutes of continuous RF (a sequence or back to back series/scan without breaks)
MR Image Artifact	The presence of this implant may produce an image artifact of 42 mm.
If more information about a specific parameter is not included, there are no conditions associated with that parameter.	

*Whole body SAR of 1.0 W/kg may limit the availability and diagnostic quality of some MR procedures

Hip and Pelvic Portfolio:

► Pelvic Plating System

Pelvic Plating System

Sterile Part Number	Non-Sterile Part Number	Description
30-0256-S	30-0256	3.5 mm x 10 mm Nonlocking Hexalobe Screw
30-0257-S	30-0257	3.5 mm x 12 mm Nonlocking Hexalobe Screw
30-0258-S	30-0258	3.5 mm x 14 mm Nonlocking Hexalobe Screw
30-0259-S	30-0259	3.5 mm x 16 mm Nonlocking Hexalobe Screw
30-0260-S	30-0260	3.5 mm x 18 mm Nonlocking Hexalobe Screw
30-0261-S	30-0261	3.5 mm x 20 mm Nonlocking Hexalobe Screw
30-0262-S	30-0262	3.5 mm x 22 mm Nonlocking Hexalobe Screw
30-0263-S	30-0263	3.5 mm x 24 mm Nonlocking Hexalobe Screw
30-0264-S	30-0264	3.5 mm x 26 mm Nonlocking Hexalobe Screw
30-0265-S	30-0265	3.5 mm x 28 mm Nonlocking Hexalobe Screw
30-0266-S	30-0266	3.5 mm x 30 mm Nonlocking Hexalobe Screw
30-0267-S	30-0267	3.5 mm x 32 mm Nonlocking Hexalobe Screw
30-0268-S	30-0268	3.5 mm x 34 mm Nonlocking Hexalobe Screw
30-0269-S	30-0269	3.5 mm x 36 mm Nonlocking Hexalobe Screw
30-0270-S	30-0270	3.5 mm x 38 mm Nonlocking Hexalobe Screw
30-0271-S	30-0271	3.5 mm x 40 mm Nonlocking Hexalobe Screw
30-0272-S	30-0272	3.5 mm x 45 mm Nonlocking Hexalobe Screw
30-0273-S	30-0273	3.5 mm x 50 mm Nonlocking Hexalobe Screw
30-0274-S	30-0274	3.5 mm x 55 mm Nonlocking Hexalobe Screw
30-0275-S	30-0275	3.5 mm x 60 mm Nonlocking Hexalobe Screw
30-0276-S	30-0276	3.5 mm x 65 mm Nonlocking Hexalobe Screw
30-0344-S	30-0344	2.7 mm x 10 mm Nonlocking Hexalobe Screw
30-0345-S	30-0345	2.7 mm x 12 mm Nonlocking Hexalobe Screw
30-0346-S	30-0346	2.7 mm x 14 mm Nonlocking Hexalobe Screw
30-0347-S	30-0347	2.7 mm x 16 mm Nonlocking Hexalobe Screw
30-0348-S	30-0348	2.7 mm x 18 mm Nonlocking Hexalobe Screw
30-0349-S	30-0349	2.7 mm x 20 mm Nonlocking Hexalobe Screw
30-0350-S	30-0350	2.7 mm x 22 mm Nonlocking Hexalobe Screw
30-0351-S	30-0351	2.7 mm x 24 mm Nonlocking Hexalobe Screw

Hip and Pelvic Portfolio [continued]

Pelvic Plating System

Sterile Part Number	Non-Sterile Part Number	Description
30-0352-S	30-0352	2.7 mm x 26 mm Nonlocking Hexalobe Screw
30-0353-S	30-0353	2.7 mm x 28 mm Nonlocking Hexalobe Screw
30-0354-S	30-0354	2.7 mm x 30 mm Nonlocking Hexalobe Screw
30-0355-S	30-0355	2.7 mm x 32 mm Nonlocking Hexalobe Screw
30-0877-S	30-0877	3.5 mm x 70 mm Nonlocking Hexalobe Screw
30-0878-S	30-0878	3.5 mm x 75 mm Nonlocking Hexalobe Screw
30-0879-S	30-0879	3.5 mm x 80 mm Nonlocking Hexalobe Screw
30-0880-S	30-0880	3.5 mm x 85 mm Nonlocking Hexalobe Screw
30-0881-S	30-0881	3.5 mm x 90 mm Nonlocking Hexalobe Screw
30-0882-S	30-0882	3.5 mm x 95 mm Nonlocking Hexalobe Screw
30-0883-S	30-0883	3.5 mm x 100 mm Nonlocking Hexalobe Screw
30-0884-S	30-0884	3.5 mm x 105 mm Nonlocking Hexalobe Screw
30-0885-S	30-0885	3.5 mm x 110 mm Nonlocking Hexalobe Screw
30-0886-S	30-0886	3.5 mm x 115 mm Nonlocking Hexalobe Screw
30-0887-S	30-0887	3.5 mm x 120 mm Nonlocking Hexalobe Screw
30-0888-S	30-0888	3.5 mm x 125 mm Nonlocking Hexalobe Screw
30-0889-S	30-0889	3.5 mm x 130 mm Nonlocking Hexalobe Screw
30-0890-S	30-0890	3.5 mm x 135 mm Nonlocking Hexalobe Screw
30-0891-S	30-0891	3.5 mm x 140 mm Nonlocking Hexalobe Screw
30-0892-S	30-0892	3.5 mm x 145 mm Nonlocking Hexalobe Screw
30-0893-S	30-0893	3.5 mm x 150 mm Nonlocking Hexalobe Screw
30-0894-S	30-0894	4.3 mm x 50 mm Hexalobe Column Screw
30-0895-S	30-0895	4.3 mm x 55 mm Hexalobe Column Screw
30-0896-S	30-0896	4.3 mm x 60 mm Hexalobe Column Screw
30-0897-S	30-0897	4.3 mm x 65 mm Hexalobe Column Screw
30-0898-S	30-0898	4.3 mm x 70 mm Hexalobe Column Screw
30-0899-S	30-0899	4.3 mm x 75 mm Hexalobe Column Screw
30-0900-S	30-0900	4.3 mm x 80 mm Hexalobe Column Screw
30-0901-S	30-0901	4.3 mm x 85 mm Hexalobe Column Screw

Hip and Pelvic Portfolio [continued]

Pelvic Plating System

Sterile Part Number	Non-Sterile Part Number	Description
30-0902-S	30-0902	4.3 mm x 90 mm Hexalobe Column Screw
30-0903-S	30-0903	4.3 mm x 95 mm Hexalobe Column Screw
30-0904-S	30-0904	4.3 mm x 100 mm Hexalobe Column Screw
30-0905-S	30-0905	4.3 mm x 105 mm Hexalobe Column Screw
30-0906-S	30-0906	4.3 mm x 110 mm Hexalobe Column Screw
30-0907-S	30-0907	4.3 mm x 115 mm Hexalobe Column Screw
30-0908-S	30-0908	4.3 mm x 120 mm Hexalobe Column Screw
30-0909-S	30-0909	4.3 mm x 125 mm Hexalobe Column Screw
30-0910-S	30-0910	4.3 mm x 130 mm Hexalobe Column Screw
30-0911-S	30-0911	4.3 mm x 135 mm Hexalobe Column Screw
30-0912-S	30-0912	4.3 mm x 140 mm Hexalobe Column Screw
30-0913-S	30-0913	4.3 mm x 145 mm Hexalobe Column Screw
30-0914-S	30-0914	4.3 mm x 150 mm Hexalobe Column Screw
70-0426-S	70-0426	Posterior Wall Acetabular Plate
70-0427-S	70-0427	Posterior Wall Acetabular Fragment Plate, Left
70-0428-S	70-0428	Posterior Wall Acetabular Fragment Plate, Right
70-0429-S	70-0429	Acetabular Spring Plate, 2-Hole
70-0430-S	70-0430	Acetabular Spring Plate, 3-Hole
70-0431-S	70-0431	Anterior Brim Plate, 12-Hole, Left
70-0432-S	70-0432	Anterior Brim Plate, 12-Hole, Right
70-0433-S	70-0433	Anterior Brim Plate, 14-Hole, Left
70-0434-S	70-0434	Anterior Brim Plate, 14-Hole, Right
70-0435-S	70-0435	Quadrilateral Surface Plate, Left
70-0436-S	70-0436	Quadrilateral Surface Plate, Right
70-0437-S	70-0437	Intrapelvic Plate, 5-Hole, Left
70-0438-S	70-0438	Intrapelvic Plate, 5-Hole, Right
70-0439-S	70-0439	Intrapelvic Plate, 9-Hole, Left
70-0440-S	70-0440	Intrapelvic Plate, 9-Hole, Right
70-0441-S	70-0441	3.5 mm Reconstruction Plate, 3-Hole

Hip and Pelvic Portfolio [continued]

Pelvic Plating System

Sterile Part Number	Non-Sterile Part Number	Description
70-0442-S	70-0442	3.5 mm Reconstruction Plate, 4-Hole
70-0443-S	70-0443	3.5 mm Reconstruction Plate, 6-Hole
70-0444-S	70-0444	3.5 mm Reconstruction Plate, 8-Hole
70-0445-S	70-0445	3.5 mm Reconstruction Plate, 10-Hole
70-0446-S	70-0446	3.5 mm Reconstruction Plate, 12-Hole
70-0447-S	70-0447	3.5 mm Reconstruction Plate, 14-Hole
70-0448-S	70-0448	3.5 mm Reconstruction Plate, 16-Hole
70-0449-S	70-0449	3.5 mm Interlocking Recon Plate, 11-Hole
70-0450-S	70-0450	Pubic Symphysis Plate, 4-Hole
70-0451-S	70-0451	Pubic Symphysis Plate, 6-Hole
70-0452-S	70-0452	Superior Sacroiliac Plate, 4-Hole
70-0458-S	70-0458	Curved Posterior Wall Acetabular Fragment Plate, Left
70-0459-S	70-0459	Curved Posterior Wall Acetabular Fragment Plate, Right

Conditions for safe scanning in the MRI

A patient with the plate and screw implants from the Pelvic Plating System may be safely scanned under the following conditions. Failure to follow these conditions may result in injury to the patient.

Nominal values(s) of Static Magnetic Field	1.5 T or 3 T
Maximum Spatial Field Gradient	19 T/m (1900 gauss/cm)
RF Excitation	Circularly Polarized (CP)
RF Transmit Coil Type	Whole body transmit coil, Head RF transmit-receive coil.
Maximum Whole Body SAR [W/kg]	2.0 W/kg
Limits on Scan Duration	2.0 W/kg whole body average SAR for 13 minutes of continuous RF (a sequence or back to back series/scan without breaks)
MR Image Artifact	The presence of this implant may produce an image artifact of 28 mm.
If more information about a specific parameter is not included, there are no conditions associated with that parameter.	


External Fixation Systems

- ▶ Stableloc External Fixation System
- ▶ Small Bone Distractor
- ▶ Small Bone Fixator

The Small Bone External Fixation System (Small Bone Distractor and Small Bone Fixator) and the Stableloc External Fixation System (shown below) are MR Unsafe.

MR Unsafe

A patient with the screw and other implants from the Stableloc External Fixation System and Small Bone External Fixation System (Small Bone Distractor and Small Bone Fixator) should not receive an MR exam. Failure to follow these conditions may result in injury to the patient.

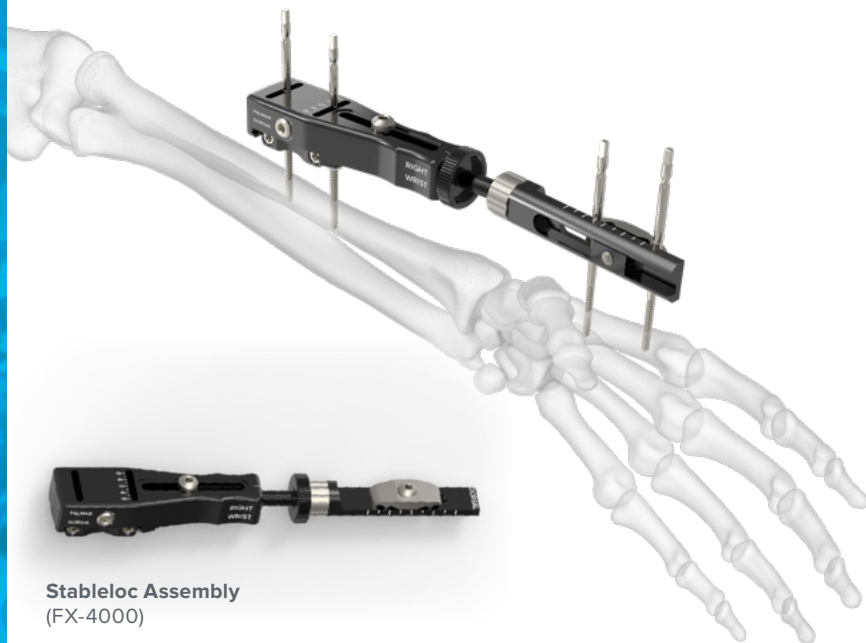


Stableloc Assembly
(FX-4000)



**3.2 mm x 30 mm
Self-Drilling Pin**
(FX-4004-S)



**2.7 mm x 25 mm
Self-Drilling Pin**
(FX-4054-S)

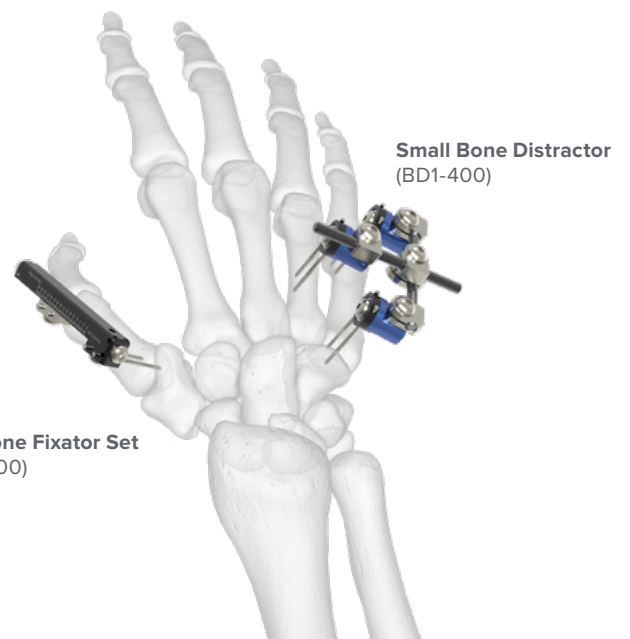

Small Bone Distractor
(BD1-400)

Small Bone Fixator Set
(SBF-0000)

Scapholunate Repair System

Scapholunate Repair System

Sterile Part Number	Non-Sterile Part Number	Description
N/A	55-0029	Transrotational Screw, 4.5 mm x 22 mm
N/A	55-0030	Transrotational Screw, 4.5 mm x 23 mm
N/A	55-0031	Transrotational Screw, 4.5 mm x 24 mm
N/A	55-0032	Transrotational Screw, 4.5 mm x 25 mm
N/A	55-0033	Transrotational Screw, 4.5 mm x 26 mm

Conditions for safe scanning in the MRI

A patient with the screw implants from the Scapholunate Repair System may be safely scanned under the following conditions. Failure to follow these conditions may result in injury to the patient.

Nominal values(s) of Static Magnetic Field	1.5 T or 3 T
Maximum Spatial Field Gradient	19 T/m
RF Excitation	Circularly Polarized (CP)
RF Transmit Coil Type	Whole body transmit coil
Maximum Whole Body SAR at 1.5 T and 3.0 T [W/kg]	2.0 W/kg
Limits on Scan Duration 1.5 T	N/A
Limits on Scan Duration 3 T	N/A
MR Image Artifact	The presence of this implant may produce an image artifact length of 46 mm \pm 5 mm including the screw and a maximal artifact width of 10 mm \pm 3mm surrounding the screw.
If more information about a specific parameter is not included, there are no conditions associated with that parameter.	

Tendon and Ligament Portfolio

- ▶ Acu-Sinch® Knotless Ankle Syndesmosis Repair System
- ▶ PEEK Suture Anchor System

Acu-Sinch Knotless Ankle Syndesmosis Repair System

Sterile Part Number	Non-Sterile Part Number	Description
46-0023-S	N/A	Acu-Sinch Knotless System w/Inserter 3.5 mm Kit
46-0024-S	N/A	Acu-Sinch Knotless System w/Inserter 3.5 mm
5502-0205-S	N/A	Acu-Sinch Washer 3.5 mm Kit

Conditions for safe scanning in the MRI

A patient with the button and washer implants from the Ankle Syndesmosis Repair System may be safely scanned under the following conditions. Failure to follow these conditions may result in injury to the patient.

Nominal values(s) of Static Magnetic Field	1.5 T or 3 T
Maximum Spatial Field Gradient	20 T/m (2000 gauss/cm)
RF Excitation	Circularly Polarized (CP)
RF Transmit Coil Type	Whole body transmit coil, Head RF transmit-receive coil
Maximum Whole Body SAR [W/kg]	4.0 W/kg
Limits on Scan Duration	4.0 W/kg whole body average SAR for 15 minutes of continuous RF (a sequence or back to back series/scan without breaks)
MR Image Artifact	The presence of this implant may produce an image artifact of 18 mm.
If more information about a specific parameter is not included, there are no conditions associated with that parameter.	

Tendon and Ligament Portfolio [continued]

PEEK Suture Anchor System

Sterile Part Number	Non-Sterile Part Number	Description
46-0008-S	N/A	4.5 mm PEEK Suture Anchor Kit
46-0009-S	N/A	5.5 mm PEEK Suture Anchor Kit
46-0012-S	N/A	4.5 mm PEEK Suture Anchor w/o Drill, Tap
46-0013-S	N/A	5.5 mm PEEK Suture Anchor w/o Drill, Tap

MRI Compatibility Claim:

PEEK is a thermoplastic polymer that is nonconductive, nonmetallic, nonmagnetic with no known hazards resulting from exposure of this device to any magnetic resonance environment, hence the device is considered MR Safe.


Chest Wall Portfolio

► RibLoc® U Plus Chest Wall Plating System

RibLoc U Plus wall plating system **when used for indications on the sternum** is MR unsafe. Otherwise, MR conditional.

MR Unsafe

A patient with the RibLoc U Plus wall plating system **implanted on the sternum** should not receive an MR exam. Failure to follow these conditions may result in injury to the patient.

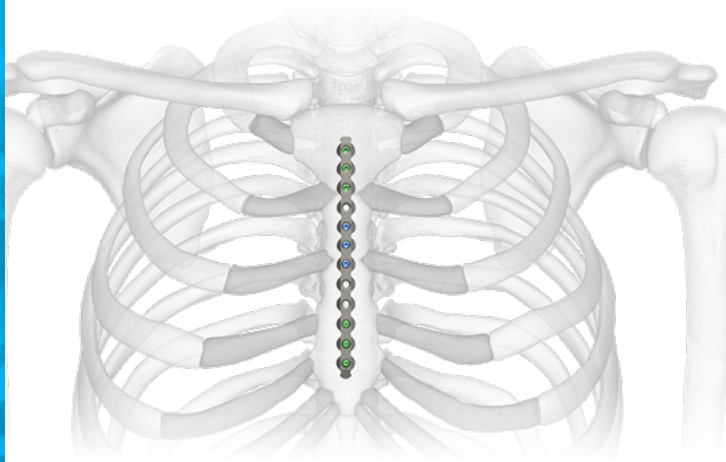

RibLoc U Plus Chest Wall Plating System

Sterile Part Number	Non-Sterile Part Number	Description
N/A	RBL1221	6 mm x 2.7 mm Locking Screw
N/A	RBL1222	8 mm x 2.7 mm Locking Screw
N/A	RBL1223	10 mm x 2.7 mm Locking Screw
N/A	RBL1224	12 mm x 2.7 mm Locking Screw
N/A	RBL1225	14 mm x 2.7 mm Locking Screw
N/A	RBL1301	50 mm Rib Plate
N/A	RBL1302	75 mm Rib Plate
N/A	RBL1303	115 mm Rib Plate
N/A	RBL1304	155 mm Rib Plate
N/A	RBL1305	215 mm Rib Plate
N/A	RBL1401	126 mm Straight Plate

Chest Wall Portfolio [continued]

Conditions for safe scanning in the MRI

A patient with the screw and plate implants from the RibLoc U Plus Chest Wall Plating System may be safely scanned under the following conditions. Failure to follow these conditions may result in injury to the patient.

Nominal values(s) of Static Magnetic Field	1.5 T or 3 T
Maximum Spatial Field Gradient	40 T/m (4000 gauss/cm)
RF Excitation	Circularly Polarized (CP)
RF Transmit Coil Type	Whole body transmit coil, Head RF transmit-receive coil
Maximum Whole Body SAR [W/kg]	2.0 W/kg
Limits on Scan Duration	2.0 W/kg whole body average SAR for 15 minutes of continuous RF (a sequence or back to back series/scan without breaks)
MR Image Artifact	The presence of this implant may produce an image artifact of 10 mm.

Note: MRI Conditional scan parameters are applicable to *rib plating only*. Sternum plating is MRI Unsafe due to temperature rise.

If more information about a specific parameter is not included, there are no conditions associated with that parameter.



www.acumed.net

Acumed USA Campus
5885 NE Cornelius Pass Road
Hillsboro, OR 97124
+1.888.627.9957

OsteoMed USA Campus
3885 Arapaho Road
Addison, TX 75001
+1.800.456.7779

Acumed Iberica Campus
C. Proci3n, 1
Edificio Oficor
28023 Madrid, Spain
+34.913.51.63.57

GEN70-56-M | Effective: 2024/03 | © 2024 Acumed® LLC

Reference

1. Woods T, Delfino J, Rajan S. Assessment of magnetically induced displacement force and torque on metal alloys used in medical devices. *J Test Eval* 49. Published online July 3, 2019. <https://doi.org/10.1520/JTE20190096>.

These materials contain information about products that may or may not be available in any particular country or may be available under different trademarks in different countries. The products may be approved or cleared by governmental regulatory organizations for sale or use with different indications or restrictions in different countries. Products may not be approved for use in all countries. Nothing contained in these materials should be construed as a promotion or solicitation for any product or for the use of any product in a particular way that is not authorized under the laws and regulations of the country where the reader is located. Nothing in these materials should be construed as a representation or warranty as to the efficacy or quality of any product, nor the appropriateness of any product to treat any specific condition. Physicians may direct questions about the availability and use of the products described in these materials to their authorized Acumed distributor. Specific questions patients may have about the use of the products described in these materials or the appropriateness for their own conditions should be directed to their own physician.

Acumed®, Polarus®, Acu-Loc®, Acutrak®, Acutrak 2®, AcuTwist®, Biotrak®, Acu-Sinch®, FlexBraid®, Hub Cap®, Tap-Loc®, Frag-Loc®, and RibLoc® are registered trademarks of Acumed LLC